



## **MITRACLIP™ G4 SYSTEM**

**MitraClip™ G4 Clip Delivery System CDS0706**

**MitraClip™ G4 Steerable Guide Catheter SGC0701**

## **MITRACLIP™ and TRICLIP™ ACCESSORIES**

**Stabilizer SZR01ST**

**Lift LFT01ST**

**Support Plate PLT01ST**

### ***Instructions for Use***

**WARNING:** Read all instructions carefully. Failure to follow these instructions, warnings and precautions may lead to device damage or patient injury. Use of the MitraClip™ G4 System should be restricted to those physicians trained to perform invasive endovascular and transseptal procedures and to those physicians trained in the proper use of the system.

### *Table of Contents*

<b>1.0</b>	<b>INDICATION FOR USE</b>
<b>2.0</b>	<b>CONTRAINDICATIONS</b>
<b>3.0</b>	<b>CLIP SIZE CONSIDERATIONS</b>
<b>4.0</b>	<b>WARNINGS</b>
<b>5.0</b>	<b>PRECAUTIONS</b>
<b>6.0</b>	<b>SPECIAL PATIENT POPULATIONS</b>
<b>7.0</b>	<b>POTENTIAL COMPLICATIONS AND ADVERSE EVENTS</b>
<b>8.0</b>	<b>PATIENT COUNSELING AND PROVIDING INFORMATION TO THE PATIENT</b>
<b>9.0</b>	<b>HOW SUPPLIED</b>
9.1	Contents
9.2	Sterile
9.3	Non-Sterile
<b>10.0</b>	<b>STORAGE</b>
<b>11.0</b>	<b>MITRACLIP™ G4 SYSTEM DIMENSIONS</b>
<b>12.0</b>	<b>GLOSSARY OF ACRONYMS</b>
<b>13.0</b>	<b>DEVICE DESCRIPTION</b>
13.1	MRI Safety Information
13.2	MitraClip™ and TriClip™ Accessories Overview
<b>14.0</b>	<b>REQUIRED ACCESSORIES</b>
<b>15.0</b>	<b>ADDITIONAL REQUIRED EQUIPMENT NOT INCLUDED</b>

---

## **16.0 OVERVIEW OF CLINICAL STUDIES**

- 16.1 EVEREST I Trial (Feasibility)
- 16.2 EVEREST II Randomized Clinical Trial (RCT)
- 16.3 EVEREST II High Risk Registry and EVEREST II REALISM Continued Access Study – High Risk (EVEREST II HRR and REALISM HR)
- 16.4 COAPT (Cardiovascular Outcomes Assessment of the MitraClip™ Percutaneous Therapy for Heart Failure Patients with Functional Mitral Regurgitation)

## **17.0 CLINICAL RESULTS IN PROHIBITIVE RISK PRIMARY MR PATIENTS**

## **18.0 SUMMARY OF THE COAPT RESULTS**

- 18.1 Study Design
- 18.2 Accountability of PMA Cohort
- 18.3 Study Population Demographics and Baseline Parameters
- 18.4 Safety and Effectiveness Results

## **19.0 POST-APPROVAL STUDY – MITRACLIP ANALYSIS COHORT (MAC)**

- 19.1 Study Objectives
- 19.2 Study Design
- 19.3 Study Population
- 19.4 Data Source
- 19.5 Key Study Endpoints
- 19.6 Total Number of Enrolled Study Sites and Subjects, Follow-Up Rate
- 19.7 Study Visits and Length of Follow-Up

## **20.0 REAL-WORLD USE SURVEILLANCE STUDIES – EXPAND AND EXPAND G4**

- 20.1 EXPAND
- 20.2 EXPAND G4

## **21.0 MITRACLIP™ G4 PROCEDURE STEP-BY-STEP INSTRUCTIONS**

- 21.1 Definition of Terms

## **22.0 PATIENT PREPARATION**

## **23.0 MITRACLIP™ G4 SYSTEM PREPARATION BEFORE USE**

- 23.1 Steerable Guide Catheter Preparation
- 23.2 Stabilizer Preparation
- 23.3 Clip Delivery System Preparation

## **24.0 ACCESS TO THE MITRAL VALVE**

## **25.0 STEERABLE GUIDE CATHETER INSERTION**

## **26.0 CLIP DELIVERY SYSTEM INSERTION**

## **27.0 INITIAL MITRACLIP™ G4 SYSTEM POSITIONING IN THE LEFT ATRIUM**

## **28.0 FINAL MITRACLIP™ G4 SYSTEM POSITIONING**

## **29.0 GRASPING THE LEAFLETS AND VERIFYING THE GRASP**

## **30.0 CLOSING THE CLIP AND EVALUATING CLIP POSITION**

## **31.0 MITRACLIP™ G4 IMPLANT PRE-DEPLOYMENT CLIP ASSESSMENT**

## **32.0 CLIP DEPLOYMENT**

- 32.1 Deployment Step 1: Lock Line Removal

- 
- 32.2 Deployment Step 2: Delivery Catheter Shaft Detachment
  - 33.0 ADDITIONAL MITRACLIP™ G4 IMPLANT PLACEMENT**
  - 34.0 MITRACLIP™ G4 SYSTEM REMOVAL**
    - 34.1 MitraClip™ G4 System Removal After Clip Deployment
    - 34.2 MitraClip™ G4 System Removal With Clip Attached
  - GRAPHICAL SYMBOLS FOR MEDICAL DEVICE LABELING**

---

## 1.0 INDICATION FOR USE

The MitraClip™ G4 System is indicated for the percutaneous reduction of significant symptomatic mitral regurgitation (MR  $\geq$  3+) due to primary abnormality of the mitral apparatus [degenerative MR] in patients who have been determined to be at prohibitive risk for mitral valve surgery by a heart team, which includes a cardiac surgeon experienced in mitral valve surgery and a cardiologist experienced in mitral valve disease, and in whom existing comorbidities would not preclude the expected benefit from reduction of the mitral regurgitation.

The MitraClip™ G4 System, when used with maximally tolerated guideline-directed medical therapy (GDMT), is indicated for the treatment of symptomatic, moderate-to-severe or severe secondary (or functional) mitral regurgitation (MR; MR  $\geq$  Grade III per American Society of Echocardiography criteria) in patients with a left ventricular ejection fraction (LVEF)  $\geq$  20% and  $\leq$  50%, and a left ventricular end systolic dimension (LVESD)  $\leq$  70 mm whose symptoms and MR severity persist despite maximally tolerated GDMT as determined by a multidisciplinary heart team experienced in the evaluation and treatment of heart failure and mitral valve disease.

## 2.0 CONTRAINDICATIONS

The MitraClip™ G4 System is contraindicated in patients with the following conditions:

- Patients who cannot tolerate, including allergy or hypersensitivity to, procedural anticoagulation or post procedural anti-platelet regimen
- Patients with known hypersensitivity to clip components (nickel / titanium, cobalt, chromium, polyester), or with contrast sensitivity
- Active endocarditis of the mitral valve
- Rheumatic mitral valve disease
- Evidence of intracardiac, inferior vena cava (IVC) or femoral venous thrombus

## 3.0 CLIP SIZE CONSIDERATIONS

When selecting a Clip length, consider the leaflet length needed to engage all the frictional elements for leaflet insertion.

- MitraClip™ G4 NT and MitraClip™ G4 NTW Clip arm length is maximum 9 mm with 6mm of leaflet insertion needed for complete frictional element engagement.
- MitraClip™ G4 XT and MitraClip™ G4 XTW Clip arm length is maximum 12 mm with 9 mm of leaflet insertion needed for complete frictional element engagement.
- Note: Mitral Valve Area and gradient concerns also need to be considered when selecting a longer Clip size.

When selecting a Clip width, consider selecting a wider Clip for broader jets and greater MR severity.

Note: Mitral Valve Area and gradient concerns also need to be considered when selecting a wider Clip size.

## 4.0 WARNINGS

- **DO NOT use the MitraClip™ Implant outside of the labeled indication.**
- The MitraClip™ G4 Implant should be implanted with sterile techniques using fluoroscopy and echocardiography (e.g. transesophageal [TEE] and transthoracic [TTE]) in a facility with on-site cardiac surgery and immediate access to a cardiac operating room.

- 
- Read all instructions carefully. Use universal precautions for biohazards and sharps while handling the MitraClip™ G4 System to avoid user injury. Failure to follow these instructions, warnings and precautions may lead to device damage, user injury, or patient injury including:
    - MitraClip™ G4 Implant erosion, migration or malposition
    - Failure to deliver MitraClip™ G4 Implant to the intended site
    - Difficulty or failure to retrieve MitraClip™ G4 system components
  - Use caution when treating patients with hemodynamic instability requiring inotropic support or mechanical heart assistance due to the increased risk of mortality in this patient population. The safety and effectiveness of the MitraClip™ Therapy in these patients has not been evaluated.
  - Patients with a rotated heart due to prior cardiac surgery in whom the System is used may have a potential risk of experiencing adverse events such as atrial perforation, cardiac tamponade, tissue damage, and embolism which may be avoided with preoperative evaluation and proper device usage.
  - For the Steerable Guide Catheter and Delivery Catheter only:
    - The Guide Catheter: the distal 65 cm of the Steerable Guide Catheter with the exception of the distal soft tip, is coated with a hydrophilic coating.
    - The Delivery Catheter: coated with a hydrophilic coating for a length of approximately 131 cm.
    - Failure to prepare the device as stated in these instructions and failure to handle the device with care could lead to additional intervention or serious adverse event.
  - The Clip Delivery System is provided sterile and designed for single use only. Cleaning, re-sterilization and / or re-use may result in infections, malfunction of the device and other serious injury or death.
  - Note the product “Use by” date specified on the package.
  - Inspect all product prior to use. Do not use if the package is open or damaged, or if product is damaged.

## 5.0 PRECAUTIONS

- Prohibitive Risk Primary (or degenerative) Mitral Regurgitation
  - Prohibitive risk is determined by the clinical judgment of a heart team, including a cardiac surgeon experienced in mitral valve surgery and a cardiologist experienced in mitral valve disease, due to the presence of one or more of the following documented surgical risk factors:
    - 30-day STS predicted operative mortality risk score of
      - ≥8% for patients deemed likely to undergo mitral valve replacement or
      - ≥6% for patients deemed likely to undergo mitral valve repair
    - Porcelain aorta or extensively calcified ascending aorta.
    - Frailty (assessed by in-person cardiac surgeon consultation)
    - Hostile chest
    - Severe liver disease / cirrhosis (MELD Score > 12)
    - Severe pulmonary hypertension (systolic pulmonary artery pressure > 2/3 systemic pressure)

- 
- Unusual extenuating circumstance, such as right ventricular dysfunction with severe tricuspid regurgitation, chemotherapy for malignancy, major bleeding diathesis, immobility, AIDS, severe dementia, high risk of aspiration, internal mammary artery (IMA) at high risk of injury, etc.
  - Evaluable data regarding safety or effectiveness is not available for prohibitive risk Primary MR patients with an LVEF < 20% or an LVESD > 60 mm. The MitraClip™ Implant should be used only when criteria for clip suitability for DMR have been met.
  - The heart team should include a cardiac surgeon experienced in mitral valve surgery and a cardiologist experienced in mitral valve disease and may also include appropriate physicians to assess the adequacy of heart failure treatment and valvular anatomy.
  - Secondary Mitral Regurgitation
    - Evaluable data regarding safety or effectiveness is not available for secondary MR patients with an LVEF < 20% or an LVESD > 70 mm.
    - The multidisciplinary heart team should be experienced in the evaluation and treatment of heart failure and mitral valve disease and determine that symptoms and MR severity persist despite maximally tolerated GDMT.

## 6.0 SPECIAL PATIENT POPULATIONS

### *Pregnancy*

The MitraClip™ G4 System has not been tested in pregnant women. Effects on the developing fetus have not been studied. The risks and reproductive effects are unknown at this time.

### *Gender*

No safety or effectiveness related gender differences were observed in clinical studies.

### *Ethnicity*

Insufficient subject numbers prevent ethnicity-related analyses on the clinical safety and effectiveness.

### *Pediatrics*

Safety and effectiveness of the MitraClip™ G4 System has not been established in pediatric patients.

## 7.0 POTENTIAL COMPLICATIONS AND ADVERSE EVENTS

The following ANTICIPATED EVENTS have been identified as possible complications of the MitraClip™ G4 procedure.

<ul style="list-style-type: none"><li>• Allergic reactions or hypersensitivity to latex, contrast agent, anaesthesia, device materials (nickel / titanium, cobalt, chromium, polyester), and drug reactions to anticoagulation, or antiplatelet drugs</li><li>• Vascular access complications which may require transfusion or vessel repair including:<ul style="list-style-type: none"><li>○ wound dehiscence,</li><li>○ catheter site reactions,</li><li>○ Bleeding (including ecchymosis, oozing, hematoma, hemorrhage, retroperitoneal hemorrhage),</li><li>○ Arteriovenous fistula, pseudoaneurysm, aneurysm, dissection, perforation / rupture, vascular occlusion</li><li>○ Emboli (air thrombotic material, implant, device component)</li><li>○ Peripheral Nerve Injury</li></ul></li><li>• Lymphatic complications</li><li>• Pericardial complications which may require additional intervention, including:<ul style="list-style-type: none"><li>○ Pericardial effusion</li><li>○ Cardiac tamponade</li><li>○ Pericarditis</li></ul></li><li>• Cardiac complications which may require additional interventions or emergency cardiac surgery, including:<ul style="list-style-type: none"><li>○ Cardiac perforation</li><li>○ Atrial septal defect</li></ul></li><li>• Mitral valve complications, which may complicate or prevent later surgical repair, including:<ul style="list-style-type: none"><li>○ Chordal entanglement / rupture</li><li>○ Single Leaflet Device Attachment (SLDA)</li><li>○ Thrombosis</li><li>○ Dislodgement of previously implanted devices</li><li>○ Tissue damage</li><li>○ Mitral valve stenosis</li><li>○ Persistent or residual mitral regurgitation</li><li>○ Endocarditis</li></ul></li></ul>	<ul style="list-style-type: none"><li>• Cardiac arrhythmias (including conduction disorders, atrial arrhythmias, ventricular arrhythmias)</li><li>• Cardiac ischemic conditions (including myocardial infarction, myocardial ischemia, and unstable / stable angina)</li><li>• Venous thromboembolism (including deep vein thrombosis, pulmonary embolism, post procedure pulmonary embolism)</li><li>• Stroke / Cerebrovascular accident (CVA) and Transient Ischemic Attack (TIA)</li><li>• System organ failure:<ul style="list-style-type: none"><li>○ Cardio-respiratory arrest</li><li>○ Worsening heart failure</li><li>○ Pulmonary congestion</li><li>○ Respiratory dysfunction / failure / atelectasis</li><li>○ Renal insufficiency or failure</li><li>○ Shock (including cardiogenic and anaphylactic)</li></ul></li><li>• Blood cell disorders (including coagulopathy, hemolysis, and Heparin Induced Thrombocytopenia (HIT))</li><li>• Hypotension / hypertension</li><li>• Infection including:<ul style="list-style-type: none"><li>○ Urinary Tract Infection (UTI)</li><li>○ Pneumonia</li><li>○ Septicemia</li></ul></li><li>• Nausea / vomiting</li><li>• Chest pain</li><li>• Dyspnea</li><li>• Edema</li><li>• Fever or hyperthermia</li><li>• Pain</li><li>• Death</li><li>• Fluoroscopy, Transesophageal echocardiogram (TEE) and Transthoracic echocardiogram (TTE) -related complications:<ul style="list-style-type: none"><li>• Skin injury or tissue changes due to exposure to ionizing radiation</li><li>• Esophageal irritation</li><li>• Esophageal perforation</li><li>• Gastrointestinal bleeding</li></ul></li></ul>
---	---

---

## 8.0 PATIENT COUNSELING AND PROVIDING INFORMATION TO THE PATIENT

Patients undergoing any procedures known to potentially be associated with bacteremia after implantation of the MitraClip™ G4 Implant should be prescribed prophylactic antibiotic therapy prior to such procedures.

Short-term anticoagulation therapy may be necessary after mitral valve repair with the MitraClip™ G4 Implant. Prescribe anticoagulation and other medical therapy per institutional guidelines.

After placement of a MitraClip™ G4 Implant, the Clip Patient Implant Card should be filled out and the patient should be instructed to carry it at all times.

All patients should be advised to limit strenuous physical activity for at least the first month post-procedure or longer if warranted.

Physicians should consider the following in counseling patients about the MitraClip™ G4 Implant:

- Discuss the risks associated with MitraClip™ G4 Implant placement.
- Discuss why surgery is not an option for the patient.
- Discuss the risk / benefit considerations for the patient.

## 9.0 HOW SUPPLIED

### 9.1 Contents

One (1) Clip Delivery System with the MitraClip™ G4 Implant, one (1) Implant Patient Card

### 9.2 Sterile

The Clip Delivery System and Steerable Guide Catheter are sterilized with ethylene oxide gas and provided in a thermoformed tray with lid, in a sealed pouch.

Parts of the devices that are in either direct or indirect contact with circulating blood are non-pyrogenic.

Note the product “Use By” date specified on the package. DO NOT use if the “Use by” date has passed. These devices are intended for single-use only. Do not reuse. Do not resterilize. This single use device cannot be reused on another patient, as it is not designed to perform as intended after the first usage. Changes in mechanical, physical, and / or chemical characteristics introduced under conditions of repeated use, cleaning, and / or resterilization may compromise the integrity of the design and / or materials, leading to contamination due to narrow gaps and / or spaces and diminished safety and / or performance of the device. Absence of original labeling may lead to misuse and eliminate traceability. Absence of original packaging may lead to device damage, loss of sterility, and risk of injury to the patient and / or user. Inspect all product prior to use. Do not use if the package is open or damaged, or if product is damaged. Do not reinsert the CDS or Clip after single use in the patient.

The white Guide tip shape retainer and transparent protective tubing are provided sterile and pre-installed on the distal tip of the Steerable Guide Catheter. The Fasteners and the Silicone Pad used with the Stabilizer are provided sterile with the Steerable Guide Catheter. The Dilator, Fasteners and the Silicone Pad are intended for single use only. Do not reuse. Do not resterilize, as this can compromise device performance and increase the risk of cross contamination due to inappropriate reprocessing.



### 9.3 Non-Sterile

**CAUTION:** The Stabilizer, Support Plate and Lift are provided non-sterile. Follow the cleaning and sterilization instructions provided with the Stabilizer, Support Plate and Lift.

### 10.0 STORAGE


Handle with care. Store in original packaging. Keep dry. Keep away from sunlight.

### 11.0 MITRACLIP™ G4 SYSTEM DIMENSIONS


**Table 1: MitraClip™ G4 System Dimensions**

Component	Dimension			
Delivery Catheter				
Extended Length (from Sleeve curved at 90 degrees)	> 65 mm			
Steerable Sleeve				
Working Length	109.5 cm			
Catheter Distal Shaft Outer Diameter	5.3 mm (16 Fr)			
MitraClip™ G4 Implant	Clip Sizes			
	G4 NT	G4 NTW	G4 XT	G4 XTW
Grasping Width at 120 degrees (Figure 1A, Figure 2A)	17 mm minimum		22 mm minimum	
Clip Width at 180 degrees (Figure 1B and 2B)	20 mm nominal		25 mm nominal	
Arm Width (Figures 1D and 1E, Figures 2D and 2E)	4 mm maximum	6 mm maximum	4 mm maximum	6 mm maximum
Arm Length (Coaptation Length) (Figure 1C, Figure 2C)	9 mm maximum		12 mm maximum	
Steerable Guide Catheter				
Working Length	80.0 cm			
Catheter Proximal Shaft Outer Diameter	8.4 mm (25 Fr)			
Catheter Distal Shaft / Tip Outer Diameter (Septal Crossing)	7.7 mm (23 Fr)			
Dilator				
Working Length	122.0 cm			
Shaft Inner Diameter	1.0 mm (3 Fr)			
Shaft Outer Diameter	5.4 mm (16 Fr)			
Distal Tip Outer Diameter	1.5 mm (4 Fr)			

**Figure 1: MitraClip™ G4 NT and NTW Implant Dimensions**

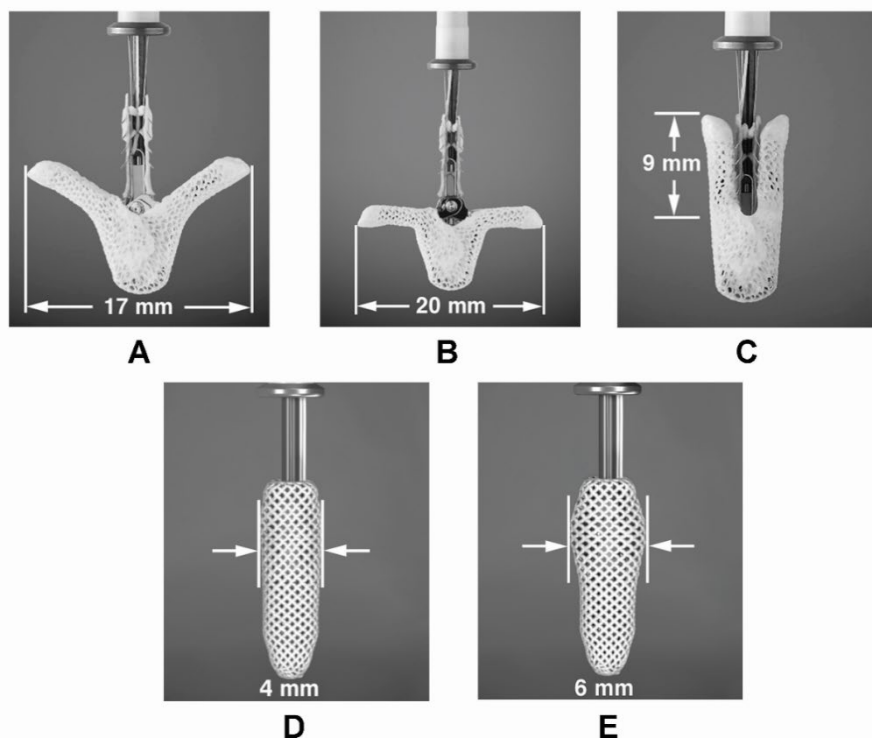

**Figure 2: MitraClip™ G4 XT and XTW Implant Dimensions**

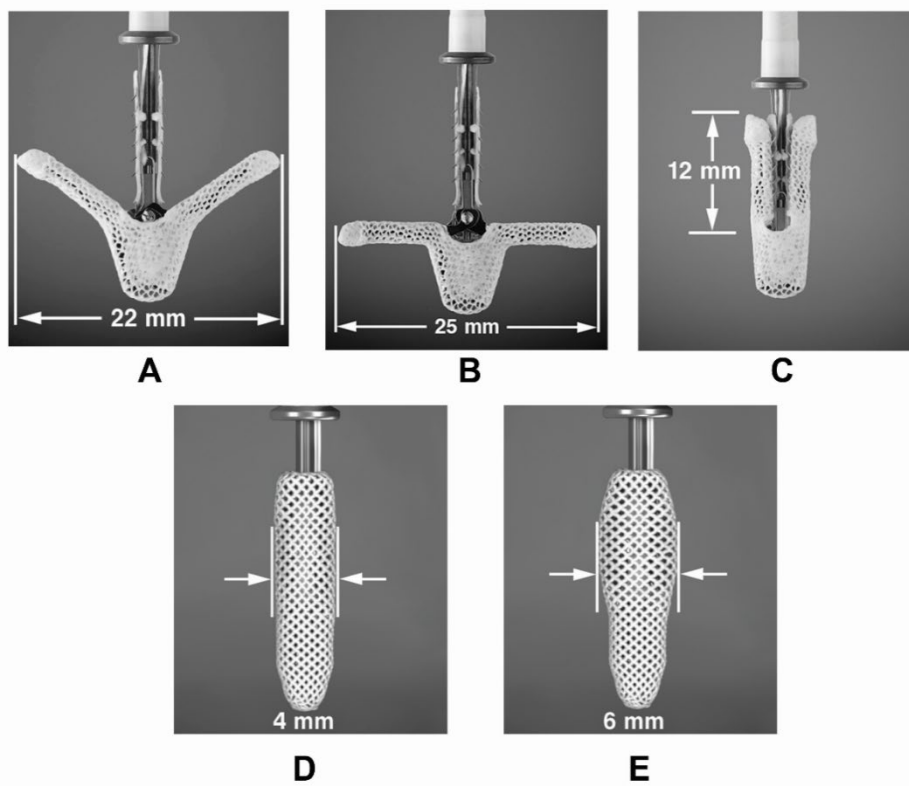


## 12.0 GLOSSARY OF ACRONYMS

Guide:	Steerable Guide Catheter	RA:	Right Atrium
CDS:	Clip Delivery System	LA:	Left Atrium
Sleeve:	Steerable Sleeve	LV:	Left Ventricle
DC:	Delivery Catheter	RO:	Radiopaque
Clip:	MitraClip™ G4 Implant	MR:	Mitral Regurgitation

## 13.0 DEVICE DESCRIPTION

The MitraClip™ G4 System consists of two parts: 1) the Clip Delivery System and 2) the Steerable Guide Catheter.

The Clip Delivery System consists of three major components: 1) the Delivery Catheter 2) the Steerable Sleeve and 3) the MitraClip™ G4 Implant. The Clip Delivery System is introduced into the body through a Steerable Guide Catheter which includes a dilator. The Clip Delivery System and Steerable Guide Catheter constitute the MitraClip™ G4 System.

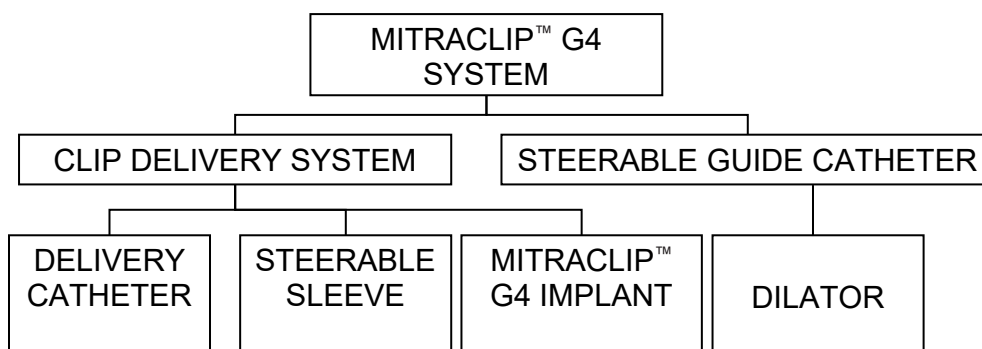

The Clip Delivery System (Figures 4 and 5) is used to advance and manipulate the implantable MitraClip™ G4 Implant for proper positioning and placement on the mitral valve leaflets. The Clip Delivery System is designed to deploy the Implant in a way that requires multiple steps to ensure safe delivery of the device.

The outer surfaces of the Delivery Catheter and the Steerable Guide Catheter have a hydrophilic coating. The coatings are intended to decrease friction during insertion into the vasculature (Steerable Guide Catheter) and decrease friction between the Steerable Sleeve and the Delivery Catheter.

The MitraClip™ G4 Implant (Figure 6) is a percutaneously implanted mechanical Clip. There are four Implant sizes in the MitraClip™ G4 product family (G4 NT, G4 XT, G4 NTW, G4 XTW). The MitraClip™ G4 Implant grasps and coapts the mitral valve leaflets resulting in fixed approximation of the mitral leaflets throughout the cardiac cycle. The MitraClip™ G4 Implant is placed without the need for arresting the heart or cardiopulmonary bypass. The implantable MitraClip™ G4 Implant is manufactured with metal alloys and polyester fabric (Clip cover) that are commonly used in cardiovascular implants.

---

The MitraClip™ G4 Implant Arms can be adjusted to any position from fully opened, fully inverted and fully closed. These positions are designed to allow the MitraClip™ G4 Implant to grasp and approximate the leaflets of the mitral valve using controls on the Delivery Catheter Handle. The MitraClip™ G4 Implant can be locked, unlocked, and repeatedly opened and closed. The Grippers can be raised or lowered repeatedly, either simultaneously or independently.

The MitraClip™ G4 Implant can be removed using standard surgical techniques and can be disposed of according to institutional guidelines.

The Steerable Guide Catheter (Figure 3a) is used to introduce the Clip Delivery System into the left side of the heart through the interatrial septum. The Steerable Guide Catheter is also used to position and orient the Clip Delivery System to the appropriate location above the mitral valve. The Steerable Guide Catheter can also be attached to a fluid filled pressure monitoring system (not provided) to measure left atrial pressures. The Dilator (Figure 3b) is used for the introduction of the Steerable Guide Catheter into the femoral vein and left atrium.

### 13.1 MRI Safety Information

Non-clinical testing has demonstrated that the MitraClip™ G4 Implants are MR Conditional. A patient with this device can be safely scanned in an MR system meeting the following conditions:



- Static magnetic field of 1.5-Tesla (1.5 T) or 3-Tesla (3.0 T)
- Maximum spatial field gradient of 4,000 Gauss/cm (40 T/m)
- Maximum MR system reported, whole body averaged specific absorption rate (SAR) of 2 W/kg (Normal Operating Mode).

Under the scan conditions defined above, MitraClip™ G4 Implants are expected to produce a maximum temperature rise of less than or equal to 3.1°C after 15 minutes of continuous scanning.

In non-clinical testing, the image artifact caused by a pair of MitraClip™ Implants extends approximately 40 mm beyond the MitraClip™ G4 Implants when imaged with a spin echo or gradient echo pulse sequence in a 3 T magnetic resonance imaging system. The presence of additional implants in a patient's valve may increase the image artifact size when imaged in an MRI system.


### 13.2 MitraClip™ and TriClip™ Accessories Overview

Several accessories are used in conjunction with the MitraClip™ G4 System including: 1) a Stabilizer, 2) a Lift, 3) a Support Plate, 4) a Silicone Pad, and 5) Fasteners. The Stabilizer is provided separately as a non-sterile reusable device and must be cleaned and sterilized prior to each use. The Stabilizer is used on the sterile field to support and position the Steerable Guide Catheter and Clip Delivery System during the procedure. The Lift and Support Plate are provided separately as non-sterile reusable devices and must be cleaned prior to each use. The Lift and Support Plate are used outside the sterile field to provide a stable platform for the Stabilizer and MitraClip™ G4 System during the procedure. Follow the cleaning and sterilization instructions provided with the Stabilizer, Support Plate and Lift. The Silicone Pad and Fasteners are single use accessories and are provided sterile with the Steerable Guide Catheter packaging. The Silicone Pad is used on the sterile field under the Stabilizer to prevent incidental movement of the Stabilizer during the procedure. The Fasteners are used on the sterile field to secure the Steerable Guide Catheter and Clip Delivery System to the Stabilizer.


## Legend of Figure Labels

<p><b>Figure 3a:</b> <b>Steerable Guide Catheter</b></p> <p>1 Hemostasis Valve 2 Alignment Marker 3 Flush Port 4 +/- Knob 5 Proximal Shaft 6 Distal Shaft 7 Radiopaque Tip Ring</p> <p><b>Figure 3b:Dilator</b></p> <p>8 Rotating Hemostatic Valve 9 Flush Port 10 Echogenic Spiral Groove</p> <p><b>Figure 4: CDS Handles</b></p> <p>11 Actuator Knob 12 Release Pin 13 Arm Positioner 14 Lock Lever Cap</p>	<p>15 Lock Lever 16 Gripper Levers 17 Gripper Lever Latch 18 Delivery Catheter Top Flush Port (Bottom Flush Port Not Shown) 19 Delivery Catheter Handle 20 Delivery Catheter Fastener 21 Sleeve Flush Port 22 A/P Knob 23 M/L Knob 24 Steerable Sleeve Handle</p> <p><b>Figure 5: CDS Distal End</b></p> <p>25 Longitudinal Alignment Marker 26 Key 27 Steerable Sleeve Shaft 28 Radiopaque Alignment Markers 28a Proximal 28b Distal</p>	<p>29 Sleeve Radiopaque Tip Ring 30 Delivery Catheter Shaft 31 Delivery Catheter Radiopaque Ring 32 MitraClip™ G4 Implant</p> <p><b>Figure 6: MitraClip™ G4 Implant Device Positions</b></p> <p>A Clip fully closed (low profile) B Clip at 180 degrees C Clip at 120 degrees D Clip at 60 degrees E Clip at 20 degrees F Clip fully inverted</p>
---	---	---


**Figure 3a: Steerable Guide Catheter**

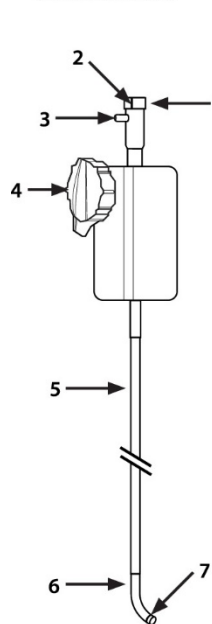


**Figure 3b: Dilator**

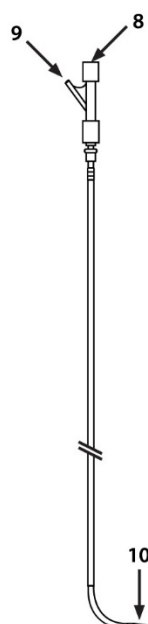

**Figure 4: CDS Handles**

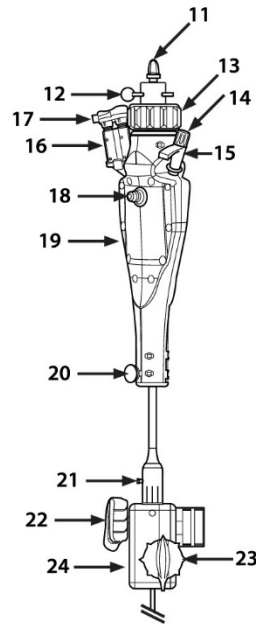

**Figure 5: CDS Distal End**

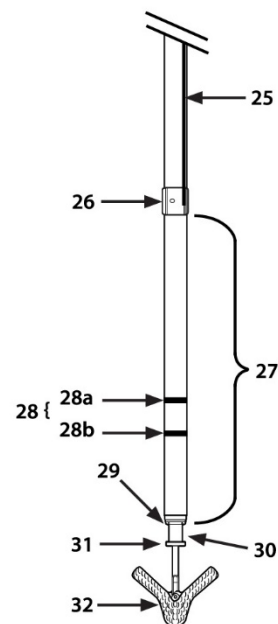

---


**Figure 6: MitraClip™ G4 Implant Positions**




**A: Clip fully closed**



**B: Clip at 180 degrees**

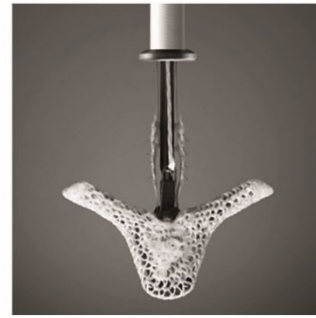


**C: Clip at 120 degrees**

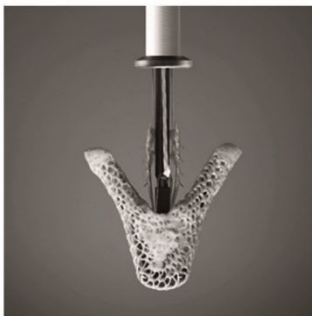

**D: Clip at 60 degrees**



**E: Clip at 20 degrees**

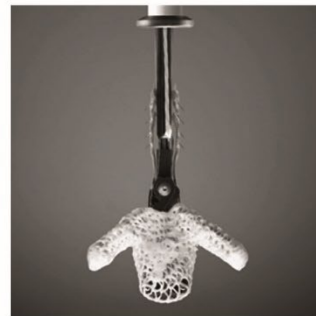

**F: Clip fully inverted**

The Steerable Guide Catheter and Clip Delivery System (Steerable Sleeve, Delivery Catheter and Clip) are steered and actuated by the use of control knobs, levers and fasteners located on the handles.

**Table 2: MitraClip™ System Handle Controls**

Device	Control	Function
Steerable Guide Catheter	+/- Knob	Tip deflection
Steerable Sleeve	M/L Knob	Tip deflection
	A/P Knob	Tip deflection
Delivery Catheter	Lock Lever	Locks-Unlocks Clip
	Gripper Levers	Raises-Lowers Grippers
	Arm Positioner	Opens-Closes Clip Arms
	DC Fastener	Locks-Unlocks DC translation / torque
	Lock Lever Cap	Controls removal of Lock Line
	Release Pin	Prevents Actuator Knob turning
	Actuator Knob	Rotates for Clip deployment

## 14.0 REQUIRED ACCESSORIES

SZR01: One (1) Stabilizer  
 LFT01: One (1) Lift  
 PLT01: One (1) Support Plate  
 One (1) Silicone Pad, three (3) Fasteners (All are included sterile with the Steerable Guide Catheter)

## 15.0 ADDITIONAL REQUIRED EQUIPMENT NOT INCLUDED

Transseptal sheath and guidewire  
 Transseptal needle  
 Exchange length supportive guidewire  
 High pressure three way stopcocks (5)  
 Arterial high pressure extension tubing (3)  
 50-60 cc syringes with luer fitting (2)  
 1000 ml pressure bags (2)  
 Sterile IV tubing with thumbwheel occluders (2)  
 Heparinized sterile saline solution (2) 1 liter bags  
 Rolling IV Pole  
 Sterile Basin

## 16.0 OVERVIEW OF CLINICAL STUDIES

Table 3 presents an overview of the MitraClip clinical program in the United States including study design, enrollment criteria, endpoints, and sample size.

**Table 3: Overview of MitraClip US Clinical Trials**

Type	Study	Key Inclusion Criteria	Key Exclusion Criteria	Endpoint	sites	patients
Feasibility	EVEREST I enrollment 2003-2006	<ul style="list-style-type: none"> <li>MR≥3+</li> <li>Symptomatic or asymptomatic with*: LVEF 30-50% and / or LVESD 50-55mm or LVEF 50-60% and LVESD &lt; 45 mm or LVEF&gt;60 and LVESD 45-55 mm</li> <li>Candidate for mitral valve surgery including cardiopulmonary bypass</li> </ul>	<ul style="list-style-type: none"> <li>LVEF&lt;30%, and / or LVESD &gt;55mm</li> <li>Mitral valve orifice area &lt;4.0 cm<sup>2</sup></li> <li>Leaflet anatomy which may preclude MitraClip device implantation, proper MitraClip device positioning on the leaflets or sufficient reduction in MR</li> </ul>	<ul style="list-style-type: none"> <li>Primary: Major Adverse Event rate through 30 days</li> </ul>	11	55
Randomized Control Trial	EVEREST II RCT enrollment 2005-2008	<ul style="list-style-type: none"> <li>MR≥3+</li> <li>Symptomatic with LVEF &gt; 25% and LVESD ≤ 55 mm or asymptomatic with*: LVEF 25% to 60% LVESD ≥ 40 mm New onset of atrial fibrillation PASP&gt;50mmHg at rest of &gt;60 mmHg with exercise</li> </ul>	<ul style="list-style-type: none"> <li>LVEF≤25%, and / or LVESD &gt;55mm</li> <li>Mitral valve orifice area &lt;4.0 cm<sup>2</sup></li> <li>Leaflet anatomy which may preclude MitraClip device implantation, proper MitraClip device positioning on the leaflets or sufficient reduction in MR</li> </ul>	<ul style="list-style-type: none"> <li>Primary Safety: Major Adverse Event rate through 30 days or discharge, whichever is greater</li> <li>Primary Effectiveness: Freedom from death, MV surgery (for Device group) or re-operation (for Control group), and MR &gt; 2+ at 12 months</li> <li>Secondary Effectiveness: <ul style="list-style-type: none"> <li>Measures of LV Function</li> <li>SF-36 quality of life</li> <li>NYHA Functional Class</li> </ul> </li> </ul>	37	60 roll-in  178 <sup>b</sup> Device  80 <sup>b</sup> Surgical Control
Single-Arm Registry	EVEREST II High Risk Registry Study  enrollment 2007-2008	<ul style="list-style-type: none"> <li>MR≥3+</li> <li>Predicted procedural mortality risk calculated using the STS surgical risk calculator of ≥ 12% or in the judgment of a cardiac surgeon the patient is considered a high risk surgical candidate due to the presence of one of the following indications: <ol style="list-style-type: none"> <li>Porcelain aorta, mobile ascending aortic atheroma</li> <li>Post-radiation mediastinum</li> <li>Previous mediastinitis</li> <li>Functional MR with EF&lt;40</li> <li>Over 75 years old with EF&lt;40</li> <li>Re-operation with patent grafts</li> <li>Two or more prior chest surgeries</li> <li>Hepatic cirrhosis</li> <li>Three or more of the following STS high risk factors <ol style="list-style-type: none"> <li>1 Creatinine &gt; 2.5 mg/dL</li> <li>2 Prior chest surgery</li> <li>3 Age over 75</li> <li>4 EF&lt;35</li> </ol> </li> </ol> </li> </ul>	<ul style="list-style-type: none"> <li>LVEF&lt;20% and / or LVESD&gt;60mm</li> <li>Mitral valve orifice area &lt;4.0 cm<sup>2</sup></li> <li>Leaflet anatomy which may preclude MitraClip device implantation, proper MitraClip device positioning on the leaflets or sufficient reduction in MR</li> </ul>	<ul style="list-style-type: none"> <li>Primary Safety: Procedural mortality at 30 days</li> <li>Major Secondary: <ul style="list-style-type: none"> <li>Measures of LV Function</li> <li>SF-36 quality of life</li> <li>NYHA Functional Class</li> <li>CHF Hospitalizations</li> </ul> </li> <li>Secondary Safety: <ul style="list-style-type: none"> <li>Major Adverse Event rate at 30 days and 12 months</li> </ul> </li> </ul>	25	78
Continued Access Registry	REALISM High Risk enrollment 2009-2013	<ul style="list-style-type: none"> <li>Same as High Risk Registry with the exception of the requirement for predicted procedural mortality risk ≥ 12%</li> </ul>	<ul style="list-style-type: none"> <li>Same as High Risk Registry</li> </ul>	<ul style="list-style-type: none"> <li>Same as High Risk Registry</li> </ul>	39	581 <sup>c</sup>
	REALISM Non-High Risk  enrollment 2009-2011	<ul style="list-style-type: none"> <li>Same as RCT</li> </ul>	<ul style="list-style-type: none"> <li>Same as RCT</li> </ul>	<ul style="list-style-type: none"> <li>Same as RCT</li> <li>6 Minute Walk Test (6MWT) Distance<sup>d</sup></li> </ul>	39	272



Type	Study	Key Inclusion Criteria	Key Exclusion Criteria	Endpoint	sites	patients
Pivotal RCT	COAPT Enrollment 2012-2017	<ul style="list-style-type: none"> <li>Ischemic or non-ischemic cardiomyopathy with LVEF 20%-50% and LVESD ≤70 mm</li> <li>Moderate-to-severe (3+) or severe (4+) secondary MR confirmed by an independent echo core laboratory prior to enrollment</li> <li>NYHA Functional Class II, III or ambulatory IV despite stable maximally-tolerated GDMT regimen and CRT (if appropriate)</li> <li>At least one hospitalization for heart failure in the 12 months prior to subject registration and/or a corrected BNP ≥300 pg/ml or corrected NT-proBNP ≥1500 pg/ml measured within 90 days prior to subject registration</li> <li>Surgery per local heart team assessment</li> <li>(LVESD) is ≤70 mm assessed by site based on a transthoracic echocardiographic (TTE)</li> <li>Treating interventionalist believes secondary MR can be successfully treated with MitraClip</li> </ul>	<ul style="list-style-type: none"> <li>ACC/AHA stage D heart failure, hemodynamic instability or cardiogenic shock</li> <li>Untreated clinically significant coronary artery disease requiring revascularization</li> <li>COPD requiring continuous home oxygen therapy or chronic oral steroid use</li> <li>Severe pulmonary hypertension or moderate or severe right ventricular dysfunction</li> <li>Aortic or tricuspid valve disease requiring surgery or transcatheter intervention</li> <li>Mitral valve orifice area &lt;4.0 cm<sup>2</sup></li> <li>Life expectancy &lt; 12 months due to non-cardiac reasons</li> </ul>	<ul style="list-style-type: none"> <li>Primary Safety: Composite of device-related complications at 12 months</li> <li>Primary Effectiveness: Recurrent HF hospitalizations through 24 months</li> </ul>	84	614 (51 roll-n)

<sup>a</sup> Inclusion criteria based on the current indication for mitral valve surgery for mitral regurgitation in the ACC/AHA guidelines for management of valvular dysfunction.

<sup>b</sup> Of the 184 patients randomized to Device, 178 received Device. Of the 95 patients randomized to Control, 80 underwent mitral valve surgery.

<sup>c</sup> As of July 12, 2013.

<sup>d</sup> In protocol version dated November 17, 2008, only patients with NYHA Functional Class III or IV in the Non-High Risk arm were considered for a 6-minute walk test. In the amended protocol version dated September 14, 2010, all patients enrolled in REALISM are required to perform the 6-minute walk test.

---

### **16.1 EVEREST I Trial (Feasibility)**

The EVEREST I trial was a prospective, multi-center, registry trial designed to evaluate the preliminary safety and effectiveness of the MitraClip device in the treatment of moderate-to-severe (3+) or severe (4+) chronic MR using up to 2 MitraClip devices per patient. The EVEREST I trial demonstrated the preliminary safety and feasibility of the MitraClip device as a percutaneous method for the reduction of MR severity. EVEREST I enrolled 55 patients at 12 US sites. Enrolled patients were required to complete clinical follow up at 30 days, 6, 18 and 24 months, and 3, 4, and 5 years. The primary safety endpoint of EVEREST I was MAE rate through 30 days (acute safety). Multiple additional secondary endpoints were pre-specified for safety and effectiveness for reporting with descriptive statistics. The study is now closed.

### **16.2 EVEREST II Randomized Clinical Trial (RCT)**

The EVEREST II RCT was a landmark trial, being the first randomized trial to compare a percutaneous intervention for the reduction of MR to standard of care mitral valve surgery. The EVEREST II RCT was a prospective, blinded, randomized, controlled, multi-center study of 279 patients (184 MitraClip, 95 surgical control) comparing the safety and effectiveness of the MitraClip to the standard of care mitral valve surgery. The intended population was patients with significant symptomatic mitral regurgitation (MR  $\geq$  3+) of either secondary MR or primary MR etiology that were non-high risk candidates indicated for and who could undergo mitral valve surgery. Study design elements including key inclusion/exclusion criteria and endpoints are provided in Table 3. Patients were evaluated at baseline, discharge, 30 days, 6, 12, 18 and 24 months, and annually thereafter through 5 years. Results of this study showed that the safety advantages of the percutaneous procedure were offset by the diminution of MR reduction with MitraClip compared to surgery, and therefore good surgical candidates should continue to receive surgical intervention.

### **16.3 EVEREST II High Risk Registry and EVEREST II REALISM Continued Access Study - High Risk (EVEREST II HRR and REALISM HR)**

The EVEREST II High Risk Registry (HRR) was a prospective, multi-center, registry designed to be adjunctive to the RCT and to evaluate the safety and effectiveness of the MitraClip device in the treatment of high surgical risk ( $\geq$  12%) patients with moderate-to-severe (3+) or severe (4+) chronic MR using up to 2 MitraClip devices per patient. The EVEREST II HRR enrolled 78 patients at 35 North American sites. Enrolled patients were required to complete clinical follow up at 30 days, 6, 12, 18 and 24 months, and 3, 4, and 5 years. The primary safety endpoint of the EVEREST II HRR was procedural mortality at 30 days or prior to discharge, whichever is longer. REALISM HR was a single-arm, self-controlled adjunctive study enrolling the same patient population as the EVEREST II HRR and designed to continue to collect safety and effectiveness data and allow patients continued access to the MitraClip during review of the PMA application.

## 16.4 COAPT (Cardiovascular Outcomes Assessment of the MitraClip™ Percutaneous Therapy for Heart Failure Patients with Functional Mitral Regurgitation)

The COAPT Trial was a prospective, randomized (1:1; MitraClip + GDMT vs. GDMT alone), open-label, multicenter investigational study intended to demonstrate: (1) MitraClip was safe in subjects with secondary MR, and (2) MitraClip could reduce recurrent HF hospitalization as compared to the GDMT Control group. The randomization was further stratified by study site and by cardiomyopathy etiology (e.g., ischemic or non-ischemic). The planned sample size of the trial was 760, including 150 roll-in subjects.

## 17.0 CLINICAL RESULTS IN PROHIBITIVE RISK PRIMARY MR PATIENTS

Data on 127 patients with significant symptomatic mitral regurgitation due to primary abnormality of the mitral apparatus (primary mitral regurgitation) determined to be at prohibitive risk for mitral valve surgery (Prohibitive Risk primary MR Cohort or PR primary MR Cohort) that were collected from the EVEREST II HRR and REALISM HR studies are provided in detail below. The analysis cohort of 127 subjects was developed post-hoc; this severely limits the statistical interpretability of reported data. These data were determined to adequately establish the safety, effectiveness, and positive benefit-risk profile of the MitraClip for the indicated population (PR primary MR) and are the basis for PMA approval. The totality of evidence demonstrates reasonable assurance of safety and effectiveness of MitraClip to reduce MR and provide patient benefit in this discreet and specific patient population.

Prohibitive Risk primary MR patients treated with the MitraClip were elderly with a high rate of serious comorbidities (Table 4).

**Table 4: Prohibitive Risk Primary MR MitraClip Cohort – Key Baseline Characteristics**

Baseline Characteristic <sup>a</sup>	Prohibitive Risk Primary MR MitraClip Patients % (n/N) (N = 127)
Age (years), Mean±SD (N)	82.4±8.7 (127)
Patients over 75 years of age	83.5% (106/127)
Female Gender	44.9% (57/127)
Body Mass Index (kg/m <sup>2</sup> ), Mean±SD (N)	25.0±5.7 (127)
Coronary Artery Disease	72.8% (91/125)
Prior Myocardial Infarction	24.4% (31/127)
Atrial Fibrillation History	70.5% (86/122)
Prior Stroke	10.2% (13/127)
Diabetes	29.9% (38/127)
Moderate to Severe Renal Disease	28.3% (36/127)
Cardiomyopathy	23.6% (30/127)
Chronic Obstructive Pulmonary Disease (w/ or w/o Home O <sub>2</sub> )	31.5% (40/127)
Hypertension	88.2% (112/127)

Baseline Characteristic <sup>a</sup>	Prohibitive Risk Primary MR MitraClip Patients % (n/N) (N = 127)
Previous Cardiovascular Surgery	48.0% (61/127)
Previous Percutaneous Coronary Intervention	33.3% (42/126)
NYHA Functional Class III/IV Heart Failure	86.6% (110/127)
LV Ejection Fraction (%), Mean±SD (N)	60.6±9.5 (112)
LV Internal Diameter, systole (cm), Mean±SD (N)	3.4±0.8 (113)
STS Mortality Risk (determined at enrollment for replacement) <sup>b</sup> , Mean±SD (N)	13.6±7.9 (127)

<sup>a</sup> Sample sizes or denominators smaller than the N reported for the group reflect missing data.

<sup>b</sup> STS replacement score calculated using the version of the calculator at the time of enrollment.

The mean Procedure Time, defined as the start time of the transseptal procedure to the time the Steerable Guide Catheter is removed, was approximately 2.5 hours (Table 5). Device time, defined as the time of insertion of the Steerable Guide Catheter to the time the MitraClip Delivery Catheter is retracted into the Steerable Guide Catheter, averaged 125 minutes. The mean fluoroscopy duration was 46 minutes. Fluoroscopy time was a relatively short proportion of the overall Procedure Time (29%). There were no intra-procedural deaths.

**Table 5: Prohibitive Risk Primary MR MitraClip Cohort - Procedural Results**

Procedural Result <sup>a</sup>	Mean±SD (N) Median (Min, Max) <sup>a</sup>
Procedure Time <sup>b</sup> (min)	157±81 (124) 134 (39, 524)
Device Time <sup>c</sup> (min)	125±75 (124) 110 (9, 511)
Fluoroscopy Duration (min)	46±26 (126) 39 (3, 167)

<sup>a</sup> Sample sizes or denominators smaller than 127 reflect missing data.

<sup>b</sup> Procedure time is measured from the time the transseptal procedure starts until the time the Steerable Guide Catheter is removed.

<sup>c</sup> Device time is measured from the time the Steerable Guide Catheter is placed in the intra-atrial septum until the time the MitraClip Delivery System is retracted into the Steerable Guide Catheter.

The MitraClip device was implanted successfully in a majority (95.3%) of patients (Table 6).

**Table 6: Prohibitive Risk MitraClip Primary MR Cohort – Number of MitraClip Devices Implanted**

# Devices Implanted	% (n/N)
0	4.7% (6/127)
1	44.1% (56/127)
2	51.2% (65/127)

Procedural mortality rate was 6.3%, which was less than both the mean and median predicted STS mortality risk using either the repair or replacement calculator.

**Table 7: Prohibitive Risk Primary MR MitraClip Cohort - Procedural Mortality**

<b>Observed Procedural Mortality, % (n/N)</b>	6.3% (8/127)
95% CI <sup>a,c</sup>	(2.8%, 12.0%)
<b>STS v2.73 Replacement Risk Score</b>	
Mean (95% CI <sup>b,c</sup> )	13.2% (11.9%, 14.5%)
Median (95% CI <sup>b,c</sup> )	12.4% (11.3%, 13.7%)
<b>STS v2.73 Repair Risk Score</b>	
Mean (95% CI <sup>b,c</sup> )	9.5% (8.5%, 10.6%)
Median (95% CI <sup>b,c</sup> )	8.5% (7.6%, 9.3%)

<sup>a</sup> Based on Clopper-Pearson method.

<sup>b</sup> CI for mean is calculated based on two-sample t-distribution and CI for median is based on non-parametric methods.

<sup>c</sup> Confidence intervals are provided to illustrate the variability of the corresponding summary statistic. They should not be used to draw statistical inference.

At 12 months, MAEs occurred at a rate of 35.4% among Prohibitive Risk primary MR MitraClip patients, with deaths (23.6%) and transfusions (19.7%) comprising the majority of events. The rate of stroke was 2.4% and rate of non-elective cardiovascular surgery was 0.8% at 12 months.

**Table 8: Prohibitive Risk Primary MR MitraClip Cohort -  
CEC Adjudicated Major Adverse Events at 30 Days and 12 Months**

Description of Event	Prohibitive Risk Primary MR MitraClip Patients (N = 127)	
	30 Days % (n/N)	12 Months % (n/N)
Death	6.3% (8/127)	23.6% (30/127)
Myocardial infarction	0.8% (1/127)	0.8% (1/127)
Re-operation for failed surgical repair or replacement	0	0
Non-elective cardiovascular surgery for adverse events	0.8% (1/127)	0.8% (1/127)
Stroke	2.4% (3/127)	2.4% (3/127)
Renal Failure	1.6% (2/127)	3.9% (5/127)
Deep wound infection	0	0.0% (0/127)
Ventilation > 48 hours	3.1% (4/127)	4.7% (6/127)
GI complication requiring surgery	0.8% (1/127)	2.4% (3/127)
New onset of permanent AF	0	0.0% (0/127)
Septicemia	0	4.7% (6/127)
Transfusion ≥ 2 units	12.6% (16/127)	19.7% (25/127)
<b>Total<sup>a</sup></b>	<b>18.9% (24/127)</b>	<b>35.4% (45/127)</b>
<b>Total<sup>a</sup> (Excluding Transfusions ≥ 2 units)</b>	<b>9.4% (12/127)</b>	<b>26.0% (33/127)</b>

<sup>a</sup> Total number of patients may not equal the sum of patients in each row since one patient may experience multiple events.

Other secondary safety endpoints occurred at a relatively low rate, consistent with access to the mitral valve achieved via the femoral vein and inferior vena cava. Major vascular complications occurred in 5.5% of patients at 30 days and in 7.1% of patients at 12 months. Major bleeding complications, defined as procedure-related bleeding requiring transfusions of at least 2 units or surgery, occurred at a rate of 12.6% at 30 days. The majority of bleeding events required transfusions rather than surgery. Bleeding events that occurred after 30 days were unrelated to the MitraClip procedure. Clinically significant atrial septal defect requiring treatment occurred at a rate of 2.4% at 12 months. A low rate (2.4%) of mitral stenosis was observed at 12 months, with a total of 3 patients reported to have experienced mitral stenosis defined as Echocardiography Core Laboratory assessed mitral valve area less than 1.5 cm<sup>2</sup> through 12 months. The site did not report mitral stenosis for these patients and none of these patients underwent mitral valve surgery for stenosis.

**Table 9: Prohibitive Risk Primary MR Cohort - Other Secondary Safety Events at 30 Days and 12 Months**

Description of Event	30 Days % (n/N)	12 Months % (n/N)
Major Vascular Complications	5.5% (7/127)	7.1% (9/127)
Major Bleeding Complications	12.6% (16/127)	15.7% (20/127)
Non-Cerebral Thromboembolism	1.6% (2/127)	1.6% (2/127)
New Onset of Persistent Atrial Fibrillation	3.9% (5/127)	3.9% (5/127)
Heart Block / Other Arrhythmia requiring Permanent Pacemaker	0.0% (0/127)	1.6% (2/127)
Endocarditis	0.0% (0/127)	0.0% (0/127)
Thrombosis	0.0% (0/127)	0.0% (0/127)
Hemolysis	0.0% (0/127)	0.0% (0/127)
Atrial Septal Defect	1.6% (2/127)	2.4% (3/127)
Mitral Valve Stenosis	0.0% (0/127)	2.4% (3/127)

The Duke University Medical Center database, which consists of patient-level data with echocardiographic, medical history and follow-up data on a large number of patients with MR  $\geq$  3+ provides a descriptive comparator for mortality. This database allowed for characterization of survival in patients deemed high risk for surgery and managed non-surgically at the Duke University Medical Center despite clear Class I indications for surgery according to the ACC/AHA Guidelines for the Management of Patients with Valvular Heart Disease. Nine hundred and fifty-three (953) patients in the Duke database with 3+ or 4+ MR were identified as too high risk for surgery based on the same high risk criteria as those in the EVEREST II HRR and REALISM studies (i.e., STS mortality risk  $\geq$  12% or protocol-specified surgical risk factors) and managed non-surgically. This made up the Duke High Risk Cohort, of which 65 patients were identified as primary MR. Table 10 shows both groups were comprised of elderly patients, with a majority of patients over the age of 75 years. The Duke High Risk primary MR Cohort reported a lower LVEF at baseline and a higher proportion of female patients than the Prohibitive Risk primary MR Cohort. The Prohibitive Risk primary MR Cohort reported a higher proportion of patients with COPD and NYHA III/IV symptoms at baseline. Both groups had high rates of previous MI, atrial fibrillation, and previous cardiovascular surgery.


Figure 7 and Table 11 display Kaplan-Meier curves comparing survival in the Prohibitive Risk primary MR patients to the Duke High Risk primary MR patients. Based on these Kaplan-Meier curves, mortality in the Prohibitive Risk primary MR Cohort was 6.4% at 30 days and 24.8% at 12 months compared to 10.9% at 30 days and 30.6% at 12 months in the Duke High Risk primary MR patients. While these results are descriptive and limited by differences described above, they suggest that there is no elevated risk of mortality in Prohibitive Risk primary MR patients who undergo the MitraClip procedure over non-surgical management.

**Table 10: Baseline and Demographic Characteristics –  
Prohibitive Risk Primary MR and Duke High Risk Primary MR Cohorts**

Baseline Characteristic	Prohibitive Risk Primary MR MitraClip Cohort % (n/N) (N = 127)	Duke High Risk Primary MR Medical Therapy Cohort % (n/N) (N = 65)
Age (years), Mean±SD (N)	82.4±8.7 (127)	76.8±11.3 (65)
Patients over 75 years of age	83.5% (106/127)	67.7% (44/65)
Male Gender	55.1% (70/127)	36.9% (24/65)
Body Mass Index (kg/m <sup>2</sup> ), Mean±SD (N)	25.0±5.7 (127)	25.4±5.0(65)
Prior Myocardial Infarction	24.4% (31/127)	33.8% (22/65)
Atrial Fibrillation History	70.5% (86/122)	58.5% (38/65)
Prior Stroke	10.2% (13/127)	18.5% (12/65)
COPD with Home Oxygen	13.4% (17/127)	6.2% (4/65)
Hypertension	88.2% (112/127)	75.4% (49/65)
Diabetes	29.9% (38/127)	36.9% (24/65)
Moderate to Severe Renal Disease	28.3% (36/127)	20.0% (13/65)
Previous Cardiovascular Surgery	48.0% (61/127)	56.9% (37/65)
Previous Percutaneous Coronary Intervention	33.3% (42/126)	58.5% (38/65)
NYHA Functional Class III/IV	86.6% (110/127)	43.8% (28/65)
STS Predicted Mortality Risk	13.2±7.3 (127)	13.3±9.0
LV Ejection Fraction (%), Mean±SD (N)	60.6±9.5 (112)	44.9±11.7 (65)
LV Internal Diameter, systole (cm), Mean±SD (N)	3.4±0.8 (113)	3.4±0.9 (65)



**Figure 7: Kaplan-Meier Freedom from Mortality –  
Prohibitive Risk Primary MR MitraClip and Duke High Risk Primary MR Medical Therapy  
Patients**

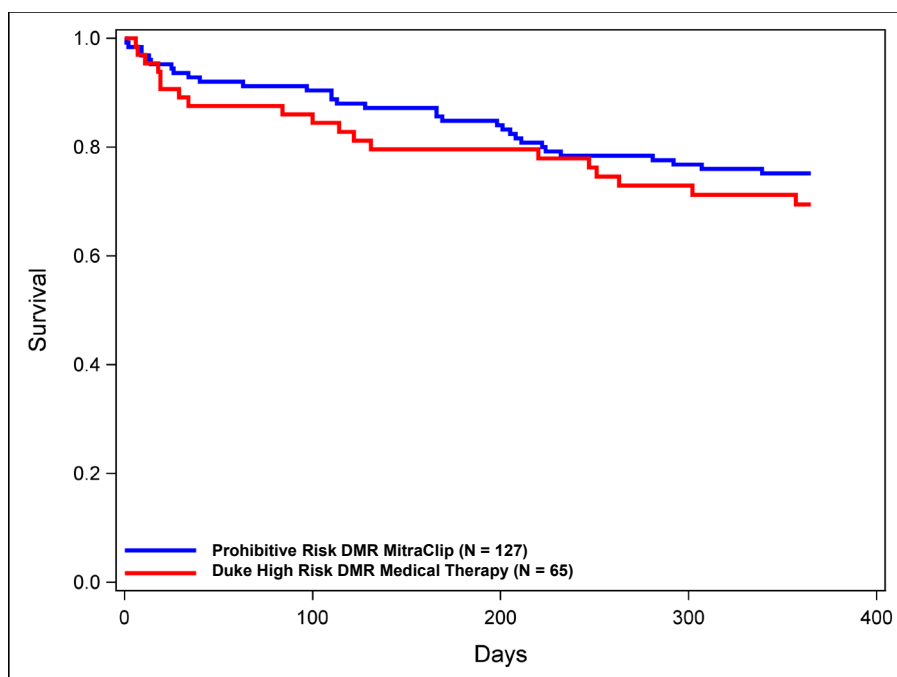

**Table 11: Number at Risk, Kaplan-Meier Estimates and 95% CIs**

Time Post Index Procedure	Baseline	30 Days	6 Months	12 Months
<b>Prohibitive Risk Primary MR MitraClip Patients (N = 127)</b>				
# At Risk	127	117	106	85
# Events	0	8	19	31
% Event Free	100%	93.6%	84.8%	75.2%
95% CI <sup>a</sup>	-	[87.6%, 96.8%]	[77.2%, 90.0%]	[66.1%, 82.1%]
<b>Duke High Risk Primary MR Medical Therapy Patients (N = 65)</b>				
# At Risk	65	57	49	39
# Events	0	7	13	19
% Event Free	100%	89.1%	79.6%	69.4%
95% CI <sup>a</sup>	-	[78.5%, 94.7%]	[67.4%, 87.6%]	[56.3%, 79.3%]

<sup>a</sup> Confidence intervals are provided to illustrate the variability of the corresponding summary statistic. They should not be used to draw statistical inference.

MR severity at baseline, discharge and 12 months are presented in Table 12 for patients with data available at each follow-up (Completers Analysis). Immediate improvement in MR severity was noted at discharge with 82.1% and 53.7% of surviving patients reporting MR severity  $\leq 2+$  and  $\leq 1+$ , respectively. This improvement was sustained at 12 months, with the majority (83.3%) of surviving patients reporting MR severity  $\leq 2+$  and 36.9% reporting MR severity  $\leq 1+$ . At 12 months, freedom from death and MR  $> 2+$  was 61.4% and freedom from death and MR  $> 1+$  was 27.2% patients.

**Table 12: Prohibitive Risk Primary MR MitraClip Cohort - MR Severity at Baseline and Follow-up Completers Analysis**

MR Severity	Baseline % (n/N)	Discharge <sup>a</sup> % (n/N)	12 Months % (n/N)
0 : None	0	1.6% (2/123)	0
1+: Mild	0	52.0% (64/123)	36.9% (31/84)
2+: Moderate	9.7% (12/124)	28.5% (35/123)	46.4% (39/84)
3+: Moderate-to-severe	58.9% (73/124)	13.0% (16/123)	13.1% (11/84)
4+: Severe	31.5% (39/124)	4.9% (6/123)	3.6% (3/84)
Missing	3	3	13
Death	0	1	30
MR $\leq 2+$ in surviving patients	9.7% (12/124)	82.1% (101/123)	83.3% (70/84)
MR $\leq 1+$ in surviving patients	0.0% (0/124)	53.7% (66/123)	36.9% (31/84)
Freedom from Death and MR $> 2+$	9.7% (12/124)	81.5% (101/124)	61.4% (70/114)
Freedom from Death and MR $> 1+$	0.0% (0/124)	53.2% (66/124)	27.2% (31/114)

<sup>a</sup> 30-day MR severity was used if discharge MR was unavailable.

Reduced preload as a result of the reduction in MR severity achieved with the MitraClip device resulted in reverse left ventricular remodeling (Table 13), characterized largely by a clinically important decrease in diastolic volume (-16.6 ml) and dimension (-0.2 cm).

**TABLE 13: PROHIBITIVE RISK PRIMARY MR MITRACLIP COHORT –  
LV Measurements at Baseline and 12 Months  
Patients with Paired Data<sup>a</sup>**

LV Measurement	N	Baseline	12-month	Difference (12-month - Baseline)	%Change (12-month - Baseline)
<b>LVEDV, ml</b>					
Mean±SD	69	125.1±40.1	108.5±37.9	-16.6±22.9	-11.5±17.9
Median		119.7	104.7	-12.3	-10.2
95% CI <sup>b,c</sup>				(-22.1, -11.1)	(-15.9, -7.2)
<b>LVIDd, cm</b>					
Mean±SD	80	5.0±0.6	4.8±0.6	-0.2±0.4	-3.7±8.2
Median		5.1	4.9	-0.2	-4.0
95% CI <sup>b,c</sup>				(-0.3, -0.1)	(-5.6, -1.9)
<b>LVESV, ml</b>					
Mean±SD	69	49.1±24.5	46.1±21.4	-3.0±13.7	-1.3±27.0
Median		45.7	41.0	-1.5	-2.7
95% CI <sup>b,c</sup>				(-6.3, 0.3)	(-7.7, 5.2)
<b>LVIDs, cm</b>					
Mean±SD	75	3.4±0.7	3.3±0.7	-0.1±0.5	-0.2±16.4
Median		3.2	3.3	-0.1	-2.3
95% CI <sup>b,c</sup>				(-0.2, 0.1)	(-4.0, 3.6)

<sup>a</sup> Only patients who had a measurement at both Baseline and 12 months are included.

<sup>b</sup> 95% CI is based on a t-distribution.

<sup>c</sup> Confidence intervals are provided to illustrate the variability of the corresponding summary statistic. They should not be used to draw statistical inference.

Improvement in LV function resulted in improvements in heart failure symptoms. NYHA Functional Class at baseline and follow-up are presented in Table 14 for patients with data available at each follow-up (Completers Analysis). Immediate improvement in NYHA Class was noted at 30 days with 82.3% of surviving patients reporting NYHA Class I or II symptoms. This improvement was sustained at 12 months, with the majority (86.9%) of surviving patients reporting NYHA Class I or II symptoms. At 12 months, freedom from death and NYHA Class III or IV symptoms was 64.0%. This improvement in NYHA Class symptoms is clinically important given that the majority of these patients (86.6%) were enrolled with NYHA Class III or IV symptoms.

**Table 14: Prohibitive Risk Primary MR MitraClip Cohort - NYHA Functional Class at Baseline and Follow-up Completers Analysis**

NYHA Functional Class	Baseline % (n/N)	30 Days % (n/N)	12 Months % (n/N)
I	2.4% (3/127)	33.6% (38/113)	40.5% (34/84)
II	11.0% (14/127)	48.7% (55/113)	46.4% (39/84)
III	63.8% (81/127)	15.9% (18/113)	10.7% (9/84)
IV	22.8% (29/127)	1.8% (2/113)	2.4% (2/84)
Missing	0	5	13
Death	0	9	30
NYHA I/II in surviving patients	13.4% (17/127)	82.3% (93/113)	86.9% (73/84)
Freedom from Death and NYHA Class III/IV	13.4% (17/127)	76.2% (93/122)	64.0% (73/114)

Table 15 shows the change in NYHA Class at 12 months from baseline. The table shows that 73 of 83 (88%) surviving patients improved by at least 1 class, and 30 of 83 (36.1%) surviving patients improved by at least 2 classes. Inclusion of deaths in the denominator results in 64.6% of patients alive and improved by at least 1 class and 26.5% alive and improved by at least 2 classes.

**Table 15: Prohibitive Risk Primary MR MitraClip Cohort – Change in NYHA Class at 12 Months from Baseline**

NYHA Class Change	Number of Patients
3 Class Improvement	4
2 Class Improvement	26
1 Class Improvement	43
No Change	9
1 Class Worsening	2
Death	30
Missing	13

Table 16 shows a mean change of +6.0 points in the Physical Component Summary (PCS) score and +5.6 points in the Mental Component Summary (MCS) score from baseline to 12 months after the MitraClip procedure. These changes are well above the 2-3 point minimally important difference (MID) threshold reported in the literature.

**Table 16: Prohibitive Risk Primary MR MitraClip Cohort – SF-36 Quality of Life at Baseline and 12 Months Completers Analysis<sup>a</sup>**

Component	N	Baseline	12-month	Difference (12-month - Baseline)
<b>Physical Component Summary Score</b>				
Mean±SD	73	33.4±8.6	39.4±10.5	6.0±8.6
Median		32.4	40.7	5.6
95% CI <sup>b,c</sup>				(4.0, 8.0)
<b>Mental Component Summary Score</b>				
Mean±SD	73	46.6±13.4	52.2±10.2	5.6±14.0
Median		49.8	54.0	3.2
95% CI <sup>b,c</sup>				(2.3, 8.9)

<sup>a</sup> Only patients who had a measurement at both Baseline and 12 months are included.

<sup>b</sup> 95% CI is based on a t-distribution.

<sup>c</sup> Confidence intervals are provided to illustrate the variability of the corresponding summary statistic. They should not be used to draw statistical inference.

The proportion of responders for both the PCS and MCS scores are shown in Table 17 based on distribution-based methods recommended by the SF-36 authors (Significant Change Criteria, SCC) and the Standard Error of Measurement (SEM) method suggested by the FDA in its 2009 PRO Guidance. The proportion of responders was 63-68% for PCS and 49-53% for MCS.

**Table 17: Prohibitive Risk Primary MR MitraClip Cohort – SF-36 QOL Responder Rate**

Component	Minimally Important Difference	Completers Analysis
<b>Physical Component Summary Score</b>	SCC <sup>a</sup> (3.1)	63.0% (46/73)
	SEM <sup>b</sup> (2.2)	68.5% (50/73)
<b>Mental Component Summary Score</b>	SCC <sup>a</sup> (3.8)	49.3% (36/73)
	SEM <sup>b</sup> (2.7)	53.4% (39/73)

<sup>a</sup> SCC (Significant Change Criteria): Significant change assuming baseline-follow-up correlation of .4 and using a 80% CI.

<sup>b</sup> SEM (Standard Error of Measurement): One SEM equals 68% CI.

A clinically important decrease in the rate of hospitalization for heart failure was observed following discharge from the MitraClip procedure (0.67 to 0.18 per patient-year, a 73% reduction, Table 18) between the pre-enrollment and the post-discharge 12-month periods.

**Table 18: Prohibitive Risk Primary MR MitraClip Cohort - Heart Failure Hospitalizations**

	12 months Pre-enrollment	Post-discharge through 12 months
# Patients for Analysis	127	120
# Patients with Events	48	13
# Events	85	17
Follow-up (Patient-Years)	127	97
Rate <sup>a</sup>	0.67	0.18
(95% Two-sided CI <sup>a,b</sup> )	(0.54, 0.83)	(0.11, 0.28)
# days hospitalized (Mean±SD)	6.0±4.5	5.9±3.8

<sup>a</sup> CI is obtained from a Poisson regression model.

<sup>b</sup> Confidence intervals are provided to illustrate the variability of the corresponding summary statistic. They should not be used to draw statistical inference.

Effectiveness results demonstrate that 82.1% (101/123) of completers experienced MR reduction from 3+ or 4+ to 2+ or less at discharge following the MitraClip procedure (Table 12). Reduction of MR at 12 months was sustained to ≤ 2+ in 83.3% (70/84), and to ≤ 1+ in 36.9% (31/84) of patients for whom echocardiographic data was available. Reduction in MR severity was associated with reverse left ventricular remodeling characterized largely by clinically important decreases in diastolic volume and dimension.

Patients also experienced clinically important improvement in NYHA Functional Class at 12 months; more than 80% of patients experienced NYHA Class III or Class IV symptoms at baseline, which reduced to less than 15% at 12 months. Despite the elderly and highly co-morbid nature of the population, quality of life as measured by the SF-36 quality of life physical and mental component scores showed clinically important improvement. Sensitivity analyses showed that these effectiveness results are robust to missing data. Finally, heart failure hospitalizations showed clinically important reduction in the 12 months post-MitraClip procedure from the 12 months pre-MitraClip procedure, including in a sensitivity analysis where death is included in the analysis as a heart failure hospitalization.

**Table 19: Effectiveness in Prohibitive Risk Primary MR MitraClip Cohort**

Effectiveness Measure <sup>§</sup>	Prohibitive Risk DMR MitraClip Cohort (N=127)
Improvement in LVEDV at 1 year	-17±23
Improvement in LVESV at 1 year	-3±14
Improvement in SF-36 PCS at 1 year	6.0±8.6
Improvement in SF-36 MCS at 1 year	5.6±14.0
NYHA Class III or IV: Baseline → 1 year	85% → 13%

<sup>§</sup> LVEDV, LVESV, SF-36 PCS and MCS results are in patients with paired data, and NYHA Class results are in Completers.


Reduction in MR severity was assessed in patients who have 2-year follow-up available. Table 20 shows that MR reduction in surviving patients to  $\leq 2+$  and  $\leq 1+$  is 82.5% (33/40) and 35.0% (14/40) respectively, at 2 years. Therefore, there is no evidence of deterioration of MR severity from 12 months to 2 years in surviving patients.

**Table 20: Prohibitive Risk Primary MR MitraClip Cohort - Durability of MR Reduction**

MR Severity	Baseline % (n/N)	12 Months % (n/N)	2 Years % (n/N)
0 : None	0	0	0
1+: Mild	0	36.9% (31/84)	35.0% (14/40)
2+: Moderate	9.7% (12/124)	46.4% (39/84)	47.5% (19/40)
3+: Moderate-to-severe	58.9% (73/124)	13.1% (11/84)	15.0% (6/40)
4+: Severe	31.5% (39/124)	3.6% (3/84)	2.5% (1/40)
MR $\leq 2+$ in surviving patients	9.7% (12/124)	83.3% (70/84)	82.5% (33/40)
MR $\leq 1+$ in surviving patients	0.0% (0/124)	36.9% (31/84)	35.0% (14/40)


In order to evaluate the relationship between MR severity and measures of effectiveness, statistical models were fit to the effectiveness data. MR severity was importantly associated with LVEDV in the Prohibitive Risk primary MR MitraClip patients (Figure 8). Reduction of MR severity to  $\leq 2+$  at 12 months resulted in clinically important decreases in LVEDV. No clinically important difference in LVEDV reduction is observed between MR 1+ and 2+. Reduction of MR to 2+ or less is associated with a decrease in left ventricular size that is not observed with ongoing MR of 3+ or greater.

**Figure 8: Prohibitive Risk Primary MR MitraClip Cohort - Dose-Response Relationship between MR Severity and Change in LVEDV 12 Months over Baseline**

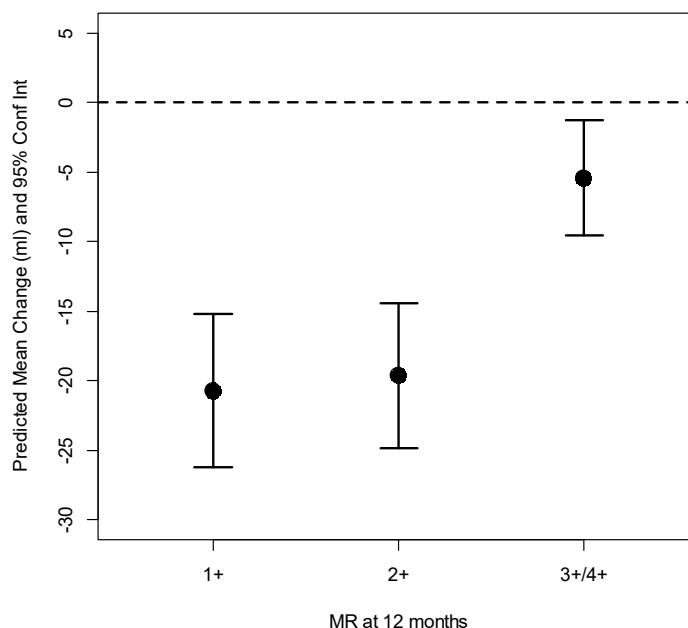

MR severity was importantly associated with PCS and MCS scores in Prohibitive Risk primary MR MitraClip patients. Reduction of MR severity to  $\leq 2+$  at 12 months resulted in clinically important improvement in PCS and MCS scores. When MR severity remained 3+/4+, the changes in PCS and MCS scores were small and not clinically important (Figure 9). Reduction of MR to 2+ or less is thus associated with an improvement in quality of life that is not observed with ongoing MR of 3+ or greater.

**Figure 9: Prohibitive Risk Primary MR MitraClip Cohort - Dose-Response Relationship between MR Severity and Change in SF-36 12 Months over Baseline**

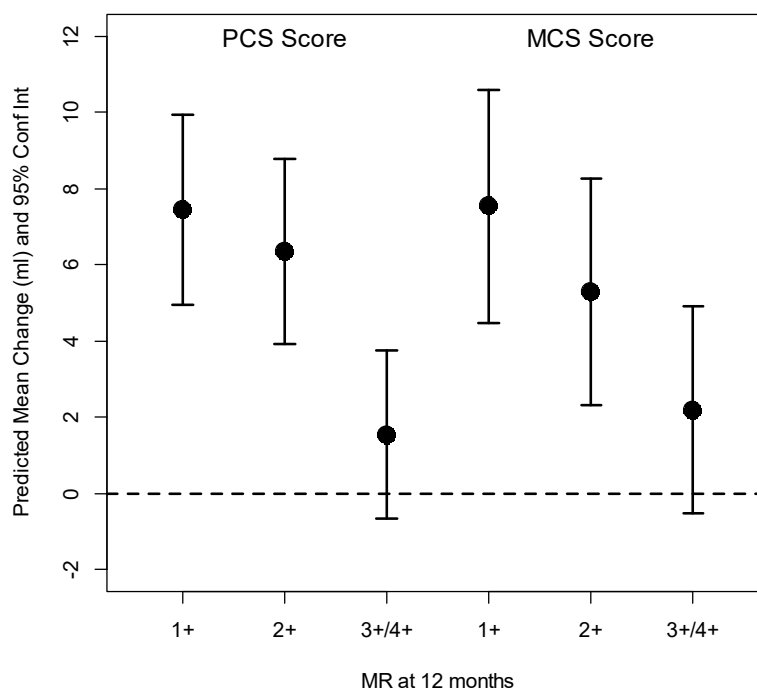

The observed number and corresponding estimated proportions of NYHA Classes at 12 months by two discharge MR groups are summarized in Table 21. The results demonstrate that reduction of MR to 2+ or less at discharge is associated with improved NYHA Functional Class that is not observed with MR of 3+ or greater at discharge.


**Table 21: Prohibitive Risk Primary MR MitraClip Cohort – Summary of Binary NYHA Functional Class Data By Discharge MR Severity**

Discharge MR	NYHA Functional Class at 12 Months	
	I/II	III/IV/Death
≤ 2+	66/93 (0.710)	27/93 (0.290)
3+/4+	7/19 (0.368)	12/19 (0.632)



Kaplan-Meier survival curves are plotted by discharge MR severity (Figure 10). There was no clinically important difference between the “≤1+” discharge MR group and the “2+” discharge MR group; however, there was a clinically important difference between the “≤1+” discharge MR group and the “3+/4+” discharge MR group and between the “2+” discharge MR group and the “3+/4+” discharge MR group. Reduction of MR to 2+ or less is associated with decreased mortality compared to ongoing MR of 3+ or greater.

**Figure 10: Prohibitive Risk Primary MR MitraClip Cohort - Kaplan-Meier Survival Curves by Discharge MR Severity (≤1+, 2+, 3+/4+)**

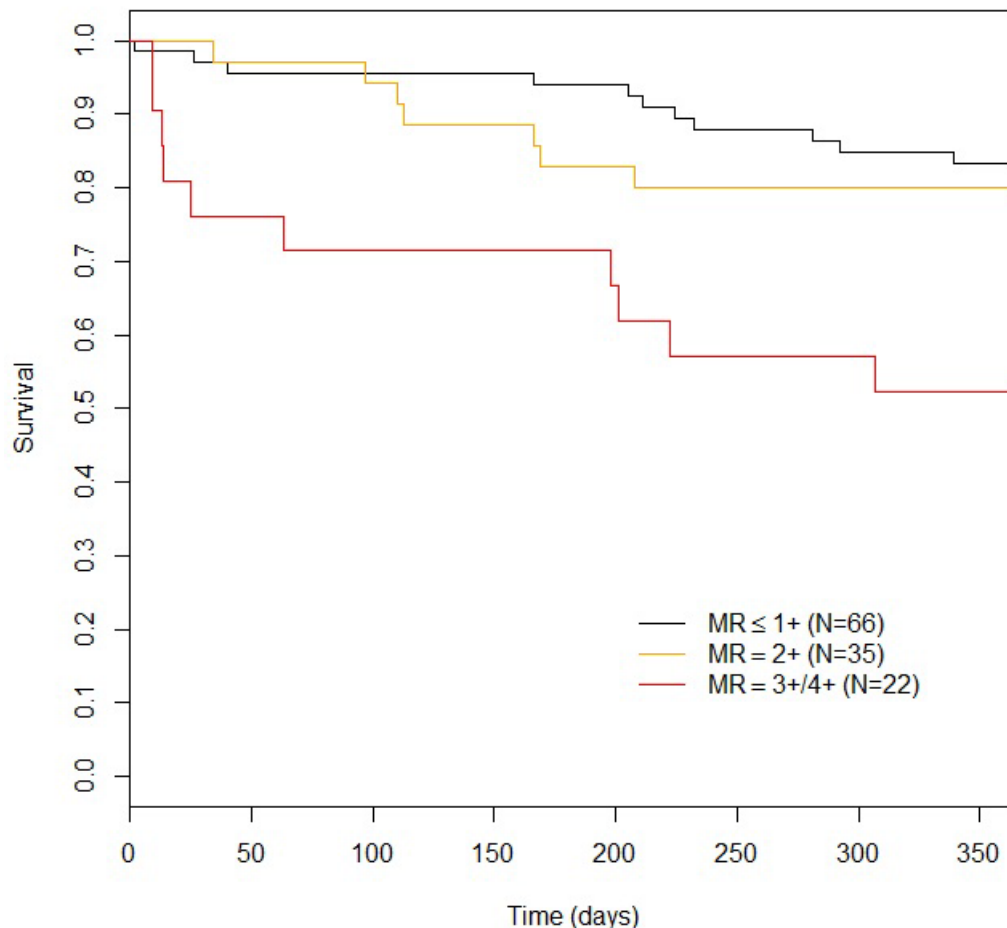

## 18.0 SUMMARY OF THE COAPT RESULTS

### 18.1 Study Design

Patients were enrolled between December 27, 2012, and June 23, 2017. The database for this Panel Track Supplement reflected data collected through September 21, 2022, and included 614 randomized patients. There were 78 investigational sites.

The COAPT Trial was a prospective, randomized (1:1; MitraClip + GDMT vs. GDMT alone), open-label, multicenter investigational study intended to demonstrate: (1) MitraClip was safe in subjects with secondary MR, and (2) MitraClip could reduce recurrent HF hospitalization as compared to the GDMT Control group.

---

The randomization was further stratified by study site and by cardiomyopathy etiology (e.g., ischemic or non-ischemic). The planned sample size of the trial was 760, including 150 roll-in subjects.

The COAPT Trial was conducted under the oversight of several independent committees, including: (1) a Steering Committee, which provided scientific and medical input on trial design, data collection, data analyses, and interpretation of results; (2) an independent Eligibility Committee, which confirmed that each subject was on optimal therapy including GDMT prior to being considered for the trial and that the subject was not appropriate for mitral valve surgery, even if randomized to the Control group; (3) a Central Echocardiography Core Laboratory (ECL), which was responsible for reviewing subject's screening echocardiography images to determine if the subject met the MR severity eligibility criterion prior to the subject being considered eligible for the trial, and for assessing MR severity and left ventricular measurements, along with other measures, at baseline and follow-ups; (4) a Clinical Events Committee (CEC), which adjudicated all adverse events per pre-established definitions (blinding was maintained whenever feasible); (5) a Data Monitoring Committee (DMC), which monitored the safety of subjects throughout the trial; and (6) a Contract Research Organization, which participated in source data verification.

#### 18.1.1 Clinical Inclusion and Exclusion Criteria

Enrollment in the COAPT Trial was limited to patients who met the following inclusion criteria:

- Symptomatic functional MR ( $\geq 3+$ ) due to cardiomyopathy of either ischemic or non-ischemic etiology determined by assessment of a qualifying transthoracic echocardiogram (TTE) obtained within 90 days and transesophageal echocardiogram (TEE) obtained within 180 days prior to subject registration, with MR severity based principally on the TTE study, confirmed by the ECL. The ECL may request a transesophageal echocardiogram (TEE) to confirm MR etiology.
- In the judgment of the HF specialist investigator at the site, the subject has been adequately treated per applicable standards, including for coronary artery disease, left ventricular dysfunction, mitral regurgitation, and HF (e.g., with CRT, revascularization, and/or GDMT. The Eligibility Committee must concur that the subject has been adequately treated.
- New York Heart Association (NYHA) Functional Class II, III or ambulatory IV.
- The Local Site Heart Team (cardiothoracic surgeon and HF specialist investigators) and the Central Eligibility Committee concur that surgery will not be offered as a treatment option and that medical therapy was the intended therapy for the subject, even if the subject was randomized to the Control group.
- Left Ventricular Ejection Fraction (LVEF) was  $\geq 20\%$  and  $\leq 50\%$  within 90 days prior to subject registration, assessed by the site using any one of the following methods: echocardiography, contrast left ventriculography, gated blood pool scan or cardiac magnetic resonance imaging (MRI).
- Left Ventricular End Systolic Dimension (LVESD) was  $\leq 70$  mm assessed by site based on a TTE obtained within 90 days prior to subject registration.
- The primary regurgitant jet was non-commissural, and in the opinion of the MitraClip implanting investigator can successfully be treated by the MitraClip. If a secondary jet exists, it must be considered clinically insignificant.

- 
- Creatine Kinase-MB (CK-MB) obtained within prior 14 days < local laboratory ULN (Upper Limit of Normal).
  - Transseptal catheterization and femoral vein access was determined to be feasible by the MitraClip implanting investigator.
  - Age 18 years or older.
  - The subject or the subject's legal representative understands and agrees that should he/she be assigned to the Control group, he/she will be treated with medical therapy and conservative management without surgery and without the MitraClip, either domestically or abroad. If the subject would actively contemplate surgery and/or MitraClip if randomized to Control, he/she should not be registered in this trial.
  - The subject or the subject's legal representative has been informed of the nature of the trial and agrees to its provisions, including the possibility of randomization to the Control group and returning for all required post-procedure follow-up visits, and has provided written informed consent.

Patients were not permitted to enroll in the COAPT Trial if they met any of the following clinical or anatomical exclusion criteria:

- Chronic Obstructive Pulmonary Disease (COPD) requiring continuous home oxygen therapy or chronic outpatient oral steroid use.
- Untreated clinically significant coronary artery disease requiring revascularization.
- Coronary artery bypass grafting (CABG) within 30 days prior to subject registration.
- Percutaneous coronary intervention within 30 days prior to subject registration.
- Transcatheter aortic valve replacement (TAVR) within 30 days prior to subject registration.
- Tricuspid valve disease requiring surgery.
- Aortic valve disease requiring surgery or transcatheter intervention.
- Cerebrovascular accident within 30 days prior to subject registration.
- Severe symptomatic carotid stenosis (> 70% by ultrasound).
- Carotid surgery or stenting within 30 days prior to subject registration.
- American College of Cardiology (ACC)/American Heart Association (AHA) Stage D heart failure.
- Presence of any of the following:
  - Estimated pulmonary artery systolic pressure (PASP) > 70 mmHg assessed by site based on echocardiography or right heart catheterization, unless active vasodilator therapy in the catheterization laboratory was able to reduce the pulmonary vascular resistance (PVR) to < 3 Wood Units or between 3 and 4.5 Wood Units with v wave less than twice the mean of the pulmonary capillary wedge pressure
  - Hypertrophic cardiomyopathy, restrictive cardiomyopathy, constrictive pericarditis, or any other structural heart disease causing heart failure other than dilated cardiomyopathy of either ischemic or non-ischemic etiology
  - Infiltrative cardiomyopathies (e.g., amyloidosis, hemochromatosis, sarcoidosis)

- 
- Hemodynamic instability requiring inotropic support or mechanical heart assistance
  - Physical evidence of right-sided congestive heart failure with echocardiographic evidence of moderate or severe right ventricular dysfunction, as assessed by site.
  - Implant of any CRT or CRT with cardioverter-defibrillator (CRT-D) within the last 30 days prior to subject registration.
  - Mitral valve orifice area  $< 4.0 \text{ cm}^2$  assessed by site based on a TTE within 90 days prior to subject registration.
  - Leaflet anatomy which may preclude MitraClip implantation, proper MitraClip positioning on the leaflets or sufficient reduction in MR by the MitraClip. This evaluation was based on TEE evaluation of the mitral valve within 180 days prior to subject registration and includes:
    - Insufficient mobile leaflet available for grasping with the MitraClip device
    - Evidence of calcification in the grasping area
    - Presence of a significant cleft in the grasping area
    - Lack of both primary and secondary chordal support in the grasping area
    - Leaflet mobility length  $< 1 \text{ cm}$
  - Hemodynamic instability defined as systolic pressure  $< 90 \text{ mmHg}$  with or without afterload reduction, cardiogenic shock or the need for inotropic support or intra-aortic balloon pump or other hemodynamic support device.
  - Need for emergent or urgent surgery for any reason or any planned cardiac surgery within the next 12 months.
  - Life expectancy  $< 12$  months due to non-cardiac conditions.
  - Modified Rankin Scale (MRS)  $\geq 4$  disability.
  - Status 1 heart transplant or prior orthotopic heart transplantation.
  - Prior mitral valve leaflet surgery or any currently implanted prosthetic mitral valve, or any prior transcatheter mitral valve procedure.
  - Echocardiographic evidence of intracardiac mass, thrombus, or vegetation.
  - Active endocarditis or active rheumatic heart disease or leaflets degenerated from rheumatic disease (i.e., noncompliant, perforated).
  - Active infections requiring current antibiotic therapy.
  - Subjects in whom TEE was contraindicated or high risk.
  - Known hypersensitivity or contraindication to procedural medications which cannot be adequately managed medically.
  - Pregnant or planning pregnancy within next 12 months.
  - Currently participating in an investigational drug or another device study that has not reached its primary endpoint. Note: Trials requiring extended follow-up for products that were investigational, but have since become commercially available, are not considered investigational trials.
  - Subject belongs to a vulnerable population per investigator's judgment or subject has any kind of disorder that compromises his/her ability to give written informed consent and/or to comply with study procedures.

### 18.1.2 Follow-up Schedule

All patients were scheduled for follow-up examinations at 1 week (phone contact), 30 days, 6 months, 12 months, 18 months, 24 months, and annually thereafter through 5 years. Preoperative and post-operative assessments included physical assessment and patient interview, laboratory measurements, imaging tests, and health status/quality of life (QoL) questionnaire. Adverse events and complications were recorded at all visits.

### 18.1.3 Clinical Endpoints

#### Primary Safety Endpoint:

The primary safety endpoint was a composite of SLDA, device embolizations, endocarditis requiring surgery, ECL confirmed mitral stenosis requiring surgery, left ventricular assist device (LVAD) implant, heart transplant, and any device related complications requiring non-elective cardiovascular surgery at 12 months.

The proportion of subjects free from the primary safety endpoint events was tested against a pre-specified performance goal (PG) of 88% for the Safety Analysis population, as defined in Section X.B.

#### Primary Effectiveness Endpoint:

The primary effectiveness endpoint was recurrent HF hospitalizations through 24 months, with the following null and alternative hypotheses:

$$H_0: RRR \leq 0$$

$$H_A: RRR > 0$$

where  $RRR$  is the relative risk reduction in the rate of recurrent HF hospitalization due to treatment with the MitraClip device as compared to the Control group. The primary effectiveness endpoint was analyzed when the last subject completed 12 months of follow-up. Hypothesis testing was performed using the Joint Frailty Model to adjust for the competing risk of death.<sup>1-3</sup>

#### Secondary Endpoints:

An ordered list of powered secondary endpoints, as shown in Table 22, was included in a hierarchical testing scheme, which were carried out after both the primary safety and effectiveness endpoints were met.

**Table 22: Ordered List of Secondary Endpoints for Hierarchical Testing**

Order	Secondary Endpoint	Alternative Hypothesis
#1	Proportion of MR severity $\leq 2+$ at 12 months	$H_A: P_D - P_C \neq 0$
#2	All-cause mortality at 12 months	$H_A: HR < 1.5$
#3	Hierarchical composite of all-cause mortality and recurrent HF hospitalization (analyzed when the last subject completes 12 months of follow-up)	$H_A$ : Either rate of death or rate of recurrent HF hospitalization is lower in the Device group compared to the Control group.
#4	Change in quality of life (QoL) at 12 months from baseline, as measured by the Kansas City Cardiomyopathy Questionnaire (KCCQ)	$H_A: \mu_D - \mu_C \neq 0$
#5	Change in distance walked on the 6 Minute Walk Test (6MWT distance or 6MWD) at 12 months over baseline	$H_A: \mu_D - \mu_C \neq 0$

Order	Secondary Endpoint	Alternative Hypothesis
#6	Recurrent hospitalizations - all-cause (analyzed when the last subject completes 12 months of follow-up)	$H_A: RRR \neq 0$
#7	Proportion of New York Heart Association (NYHA) Functional Class I/II at 12 months	$H_A: P_D - P_C \neq 0$
#8	Change in Left Ventricular End Diastolic Volume (LVEDV) at 12 months over baseline	$H_A: \mu_D - \mu_C \neq 0$
#9	All-cause mortality at 24 months	$H_A: HR \neq 1$
#10	Freedom from all-cause mortality, stroke, myocardial infarction, or non-elective cardiovascular surgery for device related complications in the MitraClip group at 30 days	$H_A: P_D(30) > 0.80$
<i>P</i> : proportion; $\mu$ : mean. <i>HR</i> : hazard ratio; <i>RRR</i> : relative risk reduction. Subscript D: Device; Subscript C: Control.		

#### Long-Term Endpoints up to 60 Months:

Key descriptive secondary endpoints through 60 months collected in the trial included the following:

- Kaplan-Meier freedom from the components of the primary safety composite at 12 months, and yearly through 60 months (Device group only)
- Kaplan-Meier freedom from the primary safety composite at 24 months and yearly through 60 months (Device group only)
- Kaplan-Meier freedom from all-cause mortality at 12 months, 24 months and yearly through 60 months
- Kaplan-Meier freedom from the first HF related hospitalization
- Proportion of MR Severity Grade < 2+ at 30 days, 6 months, 12 months, 24 months, and then yearly through 60 months
- NYHA Functional Class at baseline, 30 days, 6 months, 12 months, 24 months and then yearly through 60 months

## 18.2 Accountability of PMA Cohort

At the time of database lock, a total of 614 subjects were randomized in this trial, including 302 Device subjects and 312 Control subjects.

There were four different analysis populations defined in the protocol: Intention-to-Treat (ITT) population, Per Protocol (PP) population, As Treated (AT) population, and Safety Analysis (SA) population, as summarized in Table 23 and Figure 11. The primary analysis for safety was the Safety Analysis, and that for effectiveness was the ITT analysis.

**Table 23: Analysis Populations**

Analysis Population	Definition	Number of Patients	
		Device	Control
Intention-to-Treat (ITT)	All randomized subjects	302	312
As Treated (AT)	Randomized subjects who received the treatment as randomized	294	320
Per Protocol (PP)	Subjects who met major inclusion and none of the major exclusion criteria and received the treatment as randomized	270	289


---

Safety Analysis (SA)	All ITT subjects in the Device group with an attempted implant procedure*	293	
*Attempted implant procedure is defined as administration of anesthesia for the MitraClip procedure.			

### 18.3 **Study Population Demographics and Baseline Parameters**

The demographics and baseline characteristics of the study population are typical for an HF study performed in the U.S., as shown in Table 24. The two study groups were well-balanced, with no significant difference in patient demographics and baseline characteristics.

**Figure 11: Disposition of COAPT Randomized Subjects**

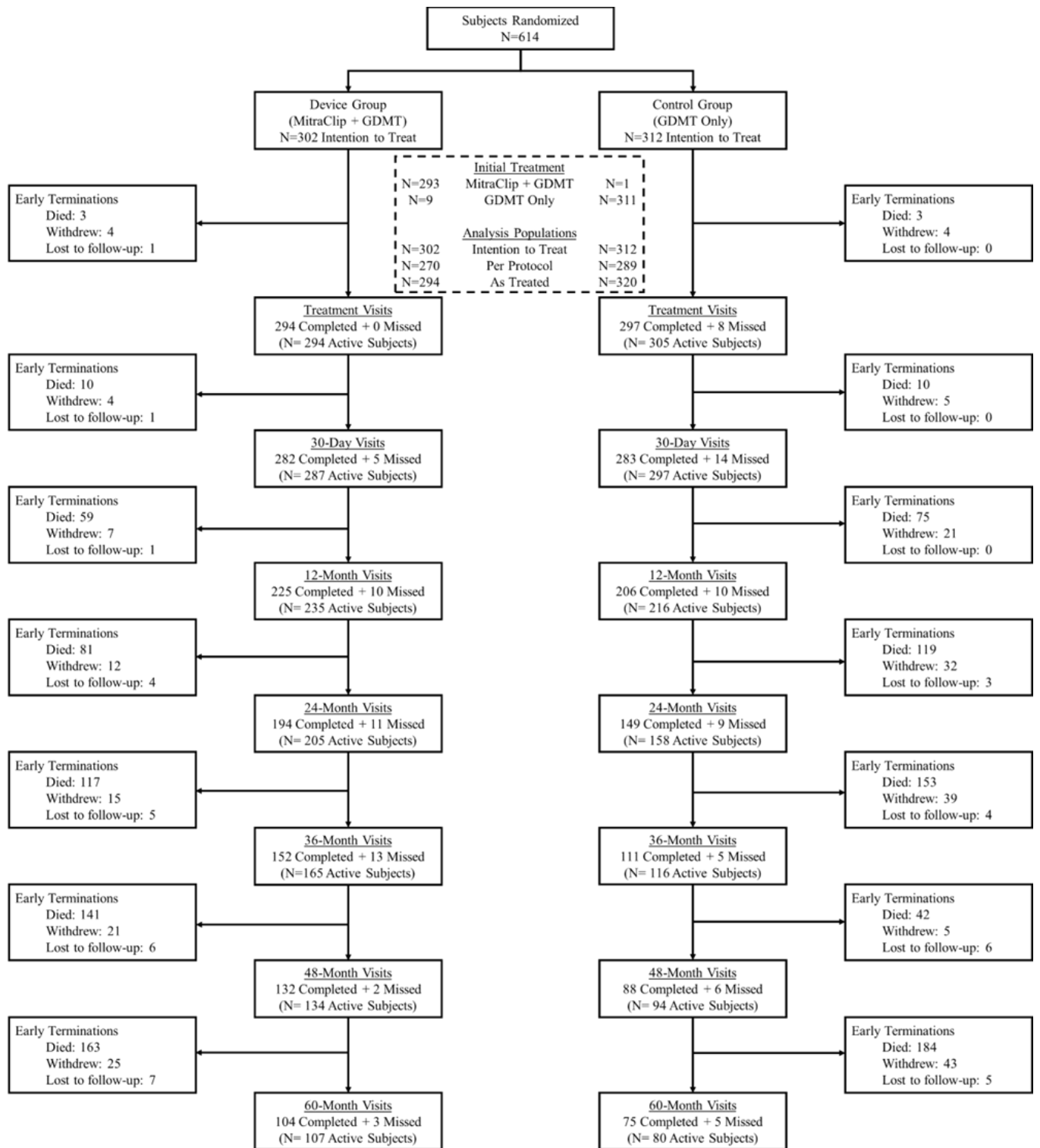



**Table 24: Patient Demographics and Baseline Characteristics (ITT Population)**

Demographics and Baseline Characteristics	Summary Statistics*		P-Value <sup>†</sup>
	Device (N=302)	Control (N=312)	
Age at Registration (year)	71.7 ± 11.8 (302)	72.8 ± 10.5 (312)	0.22
Male	66.6% (201/302)	61.5% (192/312)	0.2
Race/Ethnicity			
White or Caucasian	74.5% (225/302)	74.4% (232/312)	0.97
Non-white	25.5% (77/302)	25.6% (80/312)	
Height (cm)	170.8 ± 10.4 (301)	169.9 ± 10.8 (306)	0.25
Weight (kg)	78.8 ± 17.2 (301)	78.4 ± 20.1 (307)	0.80
Body Mass Index (kg/m <sup>2</sup> )	27.0 ± 5.8 (300)	27.1 ± 5.9 (305)	0.99
Serum Creatinine (mg/dL)	1.8 ± 1.2 (300)	1.8 ± 1.4 (306)	0.84
Creatinine Clearance (mL/min)	50.9 ± 28.5 (299)	47.8 ± 25.0 (302)	0.16
Creatinine Clearance ≤ 60 mL/min	71.6% (214/299)	75.2% (227/302)	0.32
BNP (pg/mL)	1014.8 ± 1086.0 (208)	1017.1 ± 1212.8 (209)	0.98
NT-proBNP (pg/mL)	5174.3 ± 6566.6 (74)	5943.9 ± 8437.6 (85)	0.52
Elevated BNP or NT-proBNP prior to Enrollment	93.4% (267/286)	93.1% (282/303)	0.89
Extremely High Risk for MV Surgery	68.6% (205/299)	69.9% (218/312)	0.73
KCCQ Overall Summary Score	53.2 ± 22.8 (302)	51.6 ± 23.3 (309)	0.39
Six Minute Walk Test Distance (meters)	249.6 ± 123.8 (296)	234.5 ± 123.5 (305)	0.13
NYHA Functional Class			
Class I	0.3% (1/302)	0.0% (0/311)	0.49
Class II	42.7% (129/302)	35.4% (110/311)	0.06
Class III	51.0% (154/302)	54.0% (168/311)	0.45
Class IV	6.0% (18/302)	10.6% (33/311)	0.04
SF-36 Quality of Life Physical Component Score	33.0 ± 9.1 (299)	32.6 ± 10.0 (308)	0.63
SF-36 Quality of Life Mental Component Score	46.7 ± 12.7 (299)	45.3 ± 13.0 (308)	0.19
Cardiovascular Event History			
Ischemic Cardiomyopathy	60.9% (184/302)	60.6% (189/312)	0.93
Non-Ischemic Cardiomyopathy	39.1% (118/302)	39.4% (123/312)	
Prior TIA	8.6% (26/302)	5.4% (17/312)	0.13
Prior Stroke	12.3% (37/302)	11.2% (35/312)	0.69
Prior Stroke or TIA	18.5% (56/302)	15.7% (49/312)	0.35
Prior Myocardial Infarction	51.7% (156/302)	51.3% (160/312)	0.93
Coronary Artery Disease (CAD)	72.2% (218/302)	73.1% (228/312)	0.80
Hypertension	80.5% (243/302)	80.4% (251/312)	0.99

Demographics and Baseline Characteristics	Summary Statistics*		P-Value†
	Device (N=302)	Control (N=312)	
Hypercholesterolemia	55.0% (166/302)	52.2% (163/312)	0.50
Angina	16.9% (51/302)	23.4% (73/312)	0.04
Chronic Obstructive Pulmonary Disease	23.5% (71/302)	23.1% (72/312)	0.90
Arrhythmia Event History	66.6% (201/302)	64.4% (201/312)	0.58
Ventricular Fibrillation	5.6% (17/302)	8.0% (25/312)	0.24
Ventricular Flutter	0.0% (0/302)	0.0% (0/312)	0.99
Ventricular Tachycardia	24.8% (75/302)	22.4% (70/312)	0.48
Atrial Flutter	10.3% (31/302)	10.9% (34/312)	0.80
Atrial Fibrillation	55.6% (168/302)	51.0% (159/312)	0.25
Atrial Fibrillation or Flutter	57.3% (173/302)	53.2% (166/312)	0.31
Any Hospitalization 12 months prior to enrollment	67.5% (204/302)	65.1% (203/312)	0.51
Heart Failure	58.3% (176/302)	56.1% (175/312)	0.58
Other-Cardiovascular	11.6% (35/302)	9.3% (29/312)	0.35
Non-Cardiovascular	7.9% (24/302)	7.1% (22/312)	0.67
Co-morbidity			
Diabetes	35.1% (106/302)	39.4% (123/312)	0.27
Peripheral Vascular Disease	17.2% (52/302)	18.3% (57/312)	0.73
Renal Disease	57.0% (172/302)	56.7% (177/312)	0.96
History of Anemia	22.5% (68/302)	24.4% (76/312)	0.59
History of Major Bleeds or Bleeding Disorder	7.6% (23/302)	7.1% (22/312)	0.79
STS Replacement Score (%)	7.8 ± 5.5 (302)	8.5 ± 6.2 (312)	0.16
STS Repair Score (%)	5.6 ± 5.6 (302)	6.0 ± 5.4 (312)	0.39
Prior Cardiac Interventions			
Coronary Artery Bypass Craft (CABG)	40.1% (121/302)	40.4% (126/312)	0.94
PTCA/Stents/Atherectomy	43.0% (130/302)	49.0% (153/312)	0.14
Device Implantation			
None	33.1% (100/302)	33.0% (103/312)	0.98
ICD	30.1% (91/302)	32.4% (101/312)	0.55
CRT-P	1.7% (5/302)	1.9% (6/312)	0.80
CRT-D	36.4% (110/302)	33.0% (103/312)	0.37
Pacemaker	6.0% (18/302)	8.0% (25/312)	0.32
Defibrillator (ICD or CRT-D)	62.6% (189/302)	61.5% (192/312)	0.79
Resynchronization (CRT-D or CRT-P)	38.1% (115/302)	34.9% (109/312)	0.42
Pacing (CRT-P or Pacemaker)	7.3% (22/302)	9.9% (31/312)	0.24
Prior Cardiac Valve Interventions			
Aortic Valve Intervention	3.3% (10/302)	4.5% (14/312)	0.45
Pulmonic Valve Intervention	0.0% (0/302)	0.0% (0/312)	0.99

Demographics and Baseline Characteristics	Summary Statistics*		P-Value†
	Device (N=302)	Control (N=312)	
Tricuspid Valve Intervention	0.0% (0/302)	0.0% (0/312)	0.99
Mitral Valve Intervention	0.3% (1/302)	0.0% (0/312)	0.49
Echocardiographic Core Laboratory Measures			
Mitral regurgitation severity			
3+: Moderate-to-Severe	49.0% (148/302)	55.3% (172/311)	0.12
4+: Severe	51.0% (154/302)	44.7% (139/311)	
Effective Regurgitant Orifice Area (EROA, cm <sup>2</sup> )	0.41 ± 0.15 (289)	0.40 ± 0.15 (302)	0.42
Left Ventricular Ejection Fraction (LVEF, %)	31.3 ± 9.1 (281)	31.3 ± 9.6 (294)	0.97
≤ 40 %	82.2% (231/281)	82.0% (241/294)	0.94
Left Ventricular End Systolic Dimension (LVESD, cm)	5.3 ± 0.9 (301)	5.3 ± 0.9 (306)	0.82
Left Ventricular End Diastolic Dimension (LVEDD, cm)	6.2 ± 0.7 (301)	6.2 ± 0.8 (307)	0.80
Left Ventricular End Systolic Volume (LVESV, mL)	135.5 ± 56.1 (281)	134.3 ± 60.3 (294)	0.81
Left Ventricular End Diastolic Volume (LVEDV, mL)	194.4 ± 69.2 (281)	191.0 ± 72.9 (294)	0.57
LVEDV Index (mL/m <sup>2</sup> )	102.3 ± 33.7 (279)	100.6 ± 35.0 (288)	0.56
Right Ventricular Systolic Pressure (RVSP, mmHg)	44.0 ± 13.4 (253)	44.6 ± 14.0 (275)	0.61
Medication Use at Baseline			
Beta-blocker	91.1% (275/302)	89.7% (280/312)	0.58
ACEI, ARB or ARNI	71.5% (216/302)	62.8% (196/312)	0.02
Mineralocorticoid receptor antagonist	50.7% (153/302)	49.7% (155/312)	0.81
Nitrate	6.3% (19/302)	8.0% (25/312)	0.41
Hydralazine	16.6% (50/302)	17.6% (55/312)	0.72
Diuretic	89.4% (270/302)	88.8% (277/312)	0.80
Chronic oral anticoagulant	46.4% (140/302)	40.1% (125/312)	0.12
Aspirin	57.6% (174/302)	64.7% (202/312)	0.07
P2Y12 receptor inhibitor	25.2% (76/302)	22.8% (71/312)	0.48
Statin	62.6% (189/302)	60.6% (189/312)	0.61

\*Continuous measures - Mean ± SD; categorical measures - % (no./total no.).


†P-values are from *t*-test for continuous variables and from Chi-square test or Fisher's exact test when Cochran's rule is not met for categorical variables. All p-values displayed are two-sided and for information only.

## 18.4 Safety and Effectiveness Results

### 18.4.1 Primary Safety Endpoint

The rate of freedom from device-related complications at 12 months was 96.6%, with a lower 95% confidence limit of 94.8%, which was higher than the pre-specified performance goal of 88% ( $p < 0.0001$ ), as shown in Figure 12. As such, the COAPT Trial met its primary safety endpoint. A breakdown of the composite primary safety endpoint events is presented in Table 25.

**Figure 12: Kaplan-Meier Curve of the Primary Safety Endpoint (SA Population)**

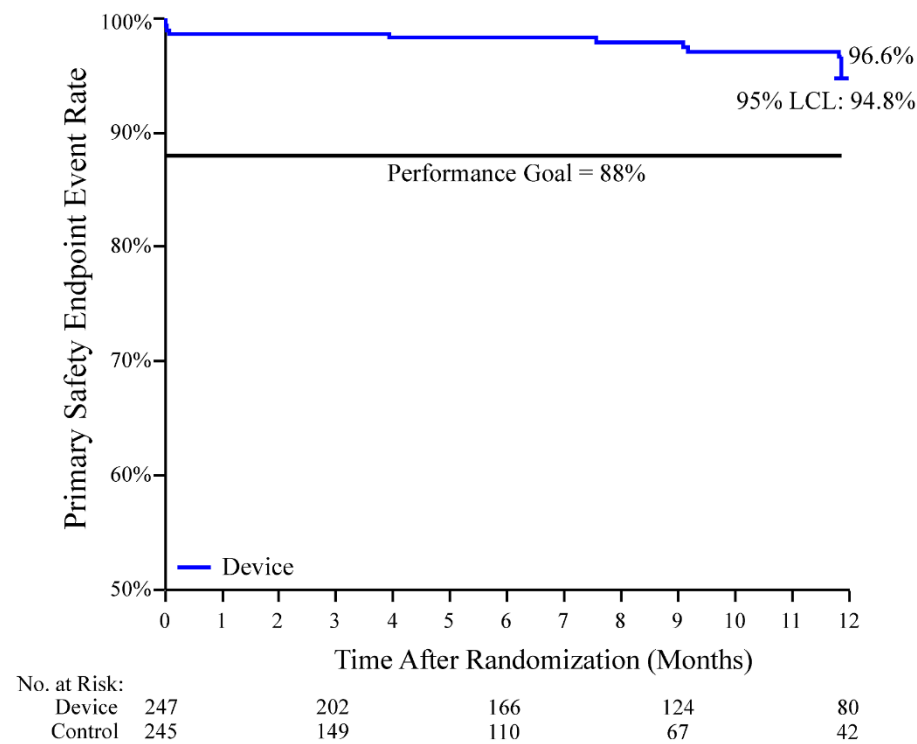

**Table 25: Outcomes of the Primary Safety Endpoint Components (SA Population)**

Event	Summary Statistics* (N = 293)
Device-related complications at 12 months	9 (3.4%)
-- Single leaflet device attachment	2 (0.7%)
-- Device embolization	1 (0.3%)
-- Endocarditis requiring surgery	0 (0.0%)
-- Mitral stenosis requiring surgery	0 (0.0%)
-- LVAD implant	3 (1.2%)
-- Heart transplant	2 (0.8%)
-- Any device-related complication requiring non-elective cardiovascular surgery	1 (0.3%)

\*# events (Kaplan-Meier rate)

#### 18.4.2 Primary Effectiveness Endpoint

A total of 160 and 283 HF hospitalizations occurred within 24 months in the Device and Control groups, respectively. The annualized rates (events per patient-year) of HF hospitalization were 0.358 in the Device group and 0.679 in the Control group, with a hazard ratio (HR) of 0.525 (upper 95% confidence limit: 0.664), representing a 47.5% reduction in the risk of recurrent HF hospitalization by the Joint Frailty Model in favor of the Device ( $p < 0.0001$ ), as summarized in Table 26 and Figure 13. Therefore, the COAPT Trial met its primary effectiveness endpoint. The successes of the primary safety endpoint and the primary effectiveness endpoint were confirmed by the AT analysis, PP analysis, and sensitivity analysis.


**Table 26: Recurrent HF Hospitalization through 24 Months – Primary Effectiveness Endpoint (ITT Population)**

	Device (N=302)	Control (N=312)	Hazard Ratio - Device vs. Control [95% CI]	Relative Risk Reduction - Device vs. Control [95% CI]	P-Value
Number of Subjects*	92 (30.5%)	151 (48.4%)			
Number of Events	160	283			
Total Follow- Up (patient- years)*	446.5	416.8			
Annualized Rate [95% CI]†	0.358 [0.307, 0.418]	0.679 [0.604, 0.763]			
<b>Joint Frailty Model</b>			0.525 [-, 0.664]	0.475 [0.336, -]	< 0.0001

\*The total follow-up in patient-years was calculated as the sum of follow-up patient-years for each subject through 24 months at the time of data cut-off or end of study, whichever was earlier.

†The annualized rate was calculated as total number of HF hospitalization events divided by total follow-up years through 24 months.

Note: (1) Hospitalizations that were adjudicated by the CEC as related to HF using the pre-specified protocol definition were included as events in the analysis; (2) Hospitalizations for MV surgery, LVAD implant or heart transplant during the follow-up period were treated as HF hospitalizations; and (3) For subjects in the Control group who received the MitraClip device due to HF or cardiac symptoms, the hospitalizations for the MitraClip procedure were treated as HF hospitalizations.

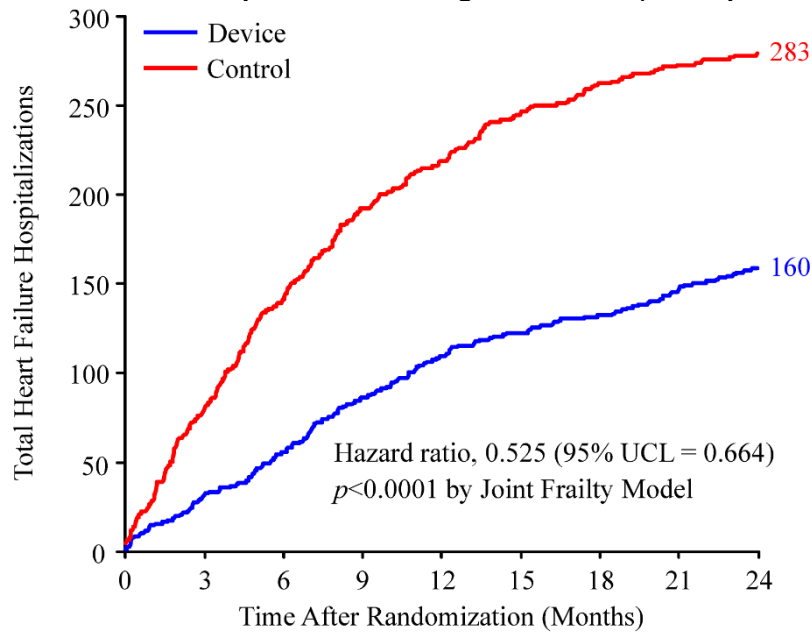**Figure 13: Total HF Hospitalization through 24 Months (ITT Population)**

No. at Risk:

Device	302	286	269	253	236	191	178	161	124
Control	312	294	271	245	219	176	145	121	88

**18.4.3 Powered Secondary Endpoints**

Hypothesis testing was performed on 10 pre-specified secondary endpoints using a hierarchical test procedure, as shown in Table 27.

**Table 27: Summary of Hierarchical Secondary Endpoints (ITT Population)**

Secondary Endpoints (hierarchical order)	Device (N=302)	Control (N=312)	HR or Difference [95% CI]	Lower 95% confidence limit	P-Value*
#1 Proportion of MR Severity $\leq 2+$ at 12 Months; % (no./total no.) [95% CI]	94.8% (199/210) [90.82%, 97.36%]	46.9% (82/175) [39.29%, 54.53%]	-	-	< 0.0001
#2 All-Cause Mortality at 12 Months (Non-inferiority); <sup>†</sup> Kaplan-Meier estimate (SE) of event rate	19.1% (2.3%)	23.2% (2.4%)	0.809 [-, 1.085]	-	0.0003
#3 Finkelstein-Schoenfeld Analysis of a Hierarchical Composite of All-Cause Mortality and Recurrent HF Hospitalization through 24 Months	-	-	-	-	< 0.0001
#4 Change in KCCQ Overall Summary Score at 12 Months over Baseline; least square means (SE) [95% CI]	12.50 (1.82) [8.93, 16.08]	-3.56 (1.85) [-7.21, 0.08]	16.07 [10.97, 21.17]	-	< 0.0001

Secondary Endpoints (hierarchical order)	Device (N=302)	Control (N=312)	HR or Difference [95% CI]	Lower 95% confidence limit	P-Value*
#5 Change in 6MWD at 12 Months over Baseline; least square means (SE) [95% CI]	-2.17 (9.12) [-20.10, 15.76]	-60.03 (8.99) [-77.69, -42.36]	57.86 [32.67, 83.05]	-	< 0.0001
#6 All-Cause Recurrent Hospitalizations through 24 Months;† annualized rate [95% CI]	1.062 [0.970, 1.162]	1.464 [1.352, 1.585]	0.760 [0.602, 0.960]	-	0.0213
#7 Proportion of NYHA Functional Class of I/II at 12-Month; % (no./total no.) [95% CI]	72.2% (171/237) [65.98%, 77.76%]	49.6% (115/232) [42.96%, 56.19%]	-	-	< 0.0001
#8 Change in Left Ventricular End Diastolic Volume at 12 Months over Baseline; least square means (SE) [95% CI]	-3.71 (5.08) [-13.71, 6.28]	17.06 (5.10) [7.03, 27.08]	-20.77 [-34.93, -6.62]	-	0.0041
#9 All-Cause Mortality through 24 Months;† Kaplan-Meier estimate (SE) of event rate	29.1% (2.8%)	46.1% (3.2%)	0.615 [0.463, 0.816]	-	0.0008
#10 Estimate of Freedom from All-Cause Mortality, Stroke, MI or Non-Elective Cardiovascular Surgery for Device-Related Complications at 30 Days; % (no./total no.)	96.9% (284/293)	-	-	94.7%	<0.0001

\*All p-values were tests for superiority, except for the secondary endpoint of mortality at 12 months (#2), which was a test for non-inferiority, and for the secondary endpoint of freedom from composite of all-cause mortality, stroke, MI or non-elective cardiovascular surgery for device-related complications at 30 days (#10), which was compared against a performance goal.


†Analyzed when the last subject completed the 12-month follow-up.

Note: (1) Imputation of worst clinical outcomes for subjects experiencing HF death prior to 12 months for the changes in KCCQ, 6MWD, LVEDV and NYHA class. (2) Continuous endpoints (KCCQ, 6MWD, and LVEDV) were analyzed using Analysis of Covariance (ANCOVA). (3) HR – Hazard Ratio; CI – Confidence Interval; SE – Standard Error.

All powered secondary endpoints were met, as summarized below:


18.4.3.1 There were significantly more subjects with MR severity  $\leq 2+$  in the Device group than in the Control group at 12 months (94.8% vs. 46.9%). The MR severity grades over time in both groups are shown in Figure 14.

**Figure 14: MR Severity Grades over Time (ITT Population)**

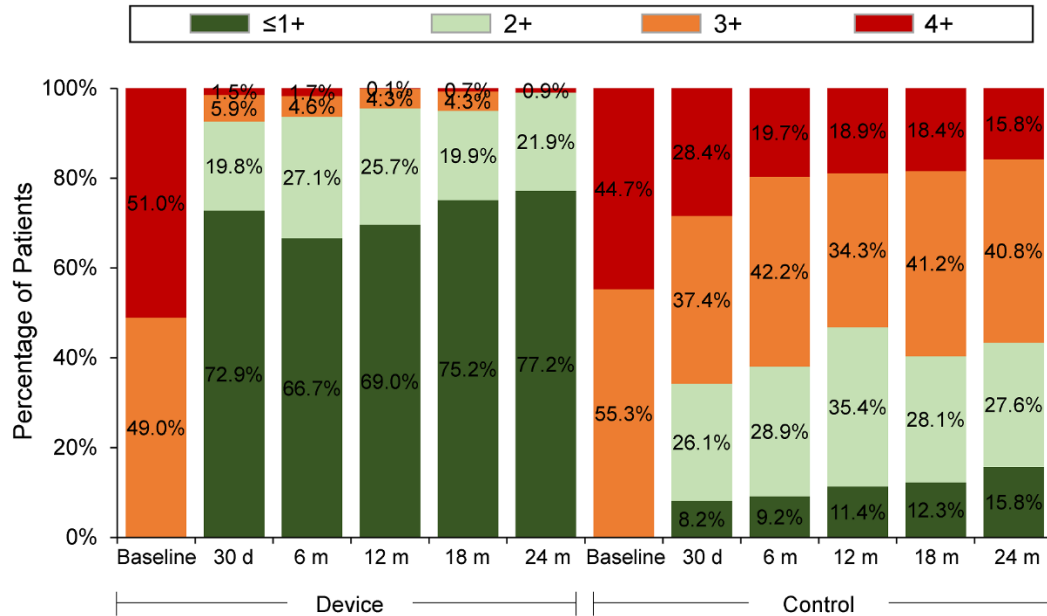

18.4.3.2 The Device group was found to be non-inferior to the Control group in all-cause mortality at 12 months (19.1% vs. 23.2%), as shown in Figure 15.

**Figure 15: Kaplan-Meier Curve of All-Cause Mortality through 12 Months (ITT Population)**

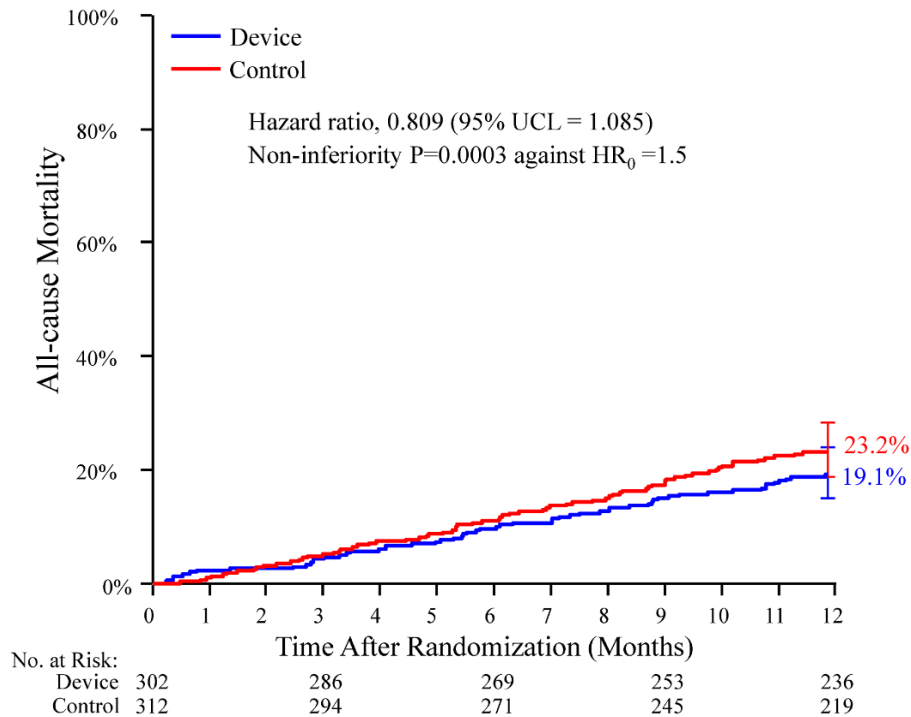




18.4.3.3 Subjects in the Device group experienced a significant reduction in the hierarchical composite of all-cause mortality and recurrent HF hospitalization compared to those in the Control group.


18.4.3.4 Subjects in the Device group experienced a significantly greater improvement in QoL (as assessed by the change in KCCQ Overall Summary Score at 12 months over baseline) compared to those in the Control group (12.50 vs. -3.56), as shown in Figure 16.

**Figure 16: KCCQ Overall Summary Score at Baseline and 12 Months (ITT Population)**

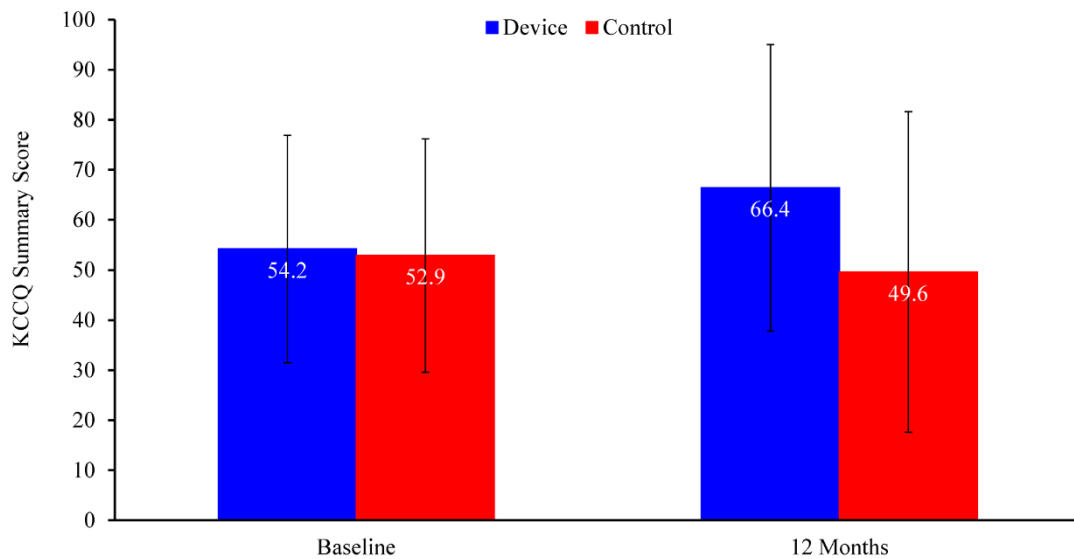


18.4.3.5 Subjects in the Device group experienced significantly greater preservation of functional capacity (as assessed by the change in 6MWD at 12 months over baseline) compared to those in the Control group (-2.17 m vs. -60.03 m), as shown in Figure 17.

**Figure 17: 6MWD at Baseline and 12 Months (ITT Population)**

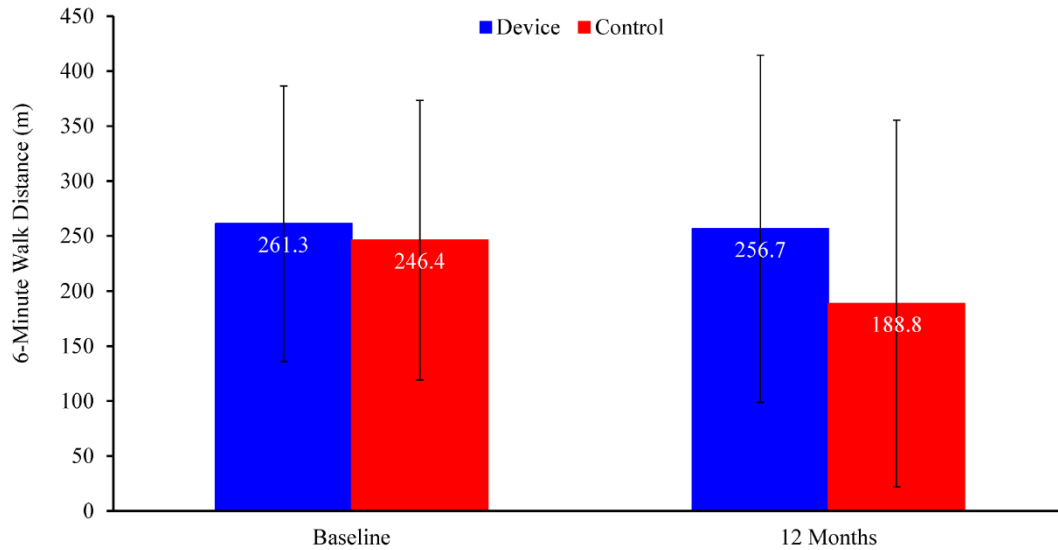


18.4.3.6 Subjects in the Device group experienced a significantly lower annualized rate (events per patient-year) of all-cause hospitalizations compared to those in the Control group (1.062 vs. 1.464). The total all-cause hospitalization through 24 months is shown in Figure 18.

**Figure 18: Total All-Cause Hospitalization through 24 Months (ITT Population)**

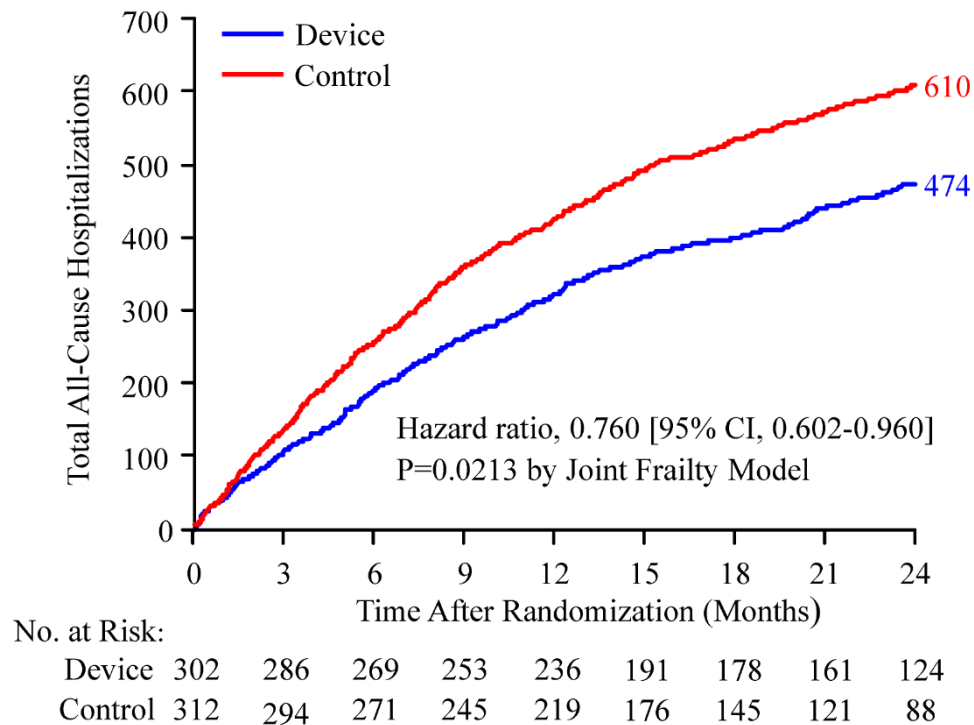


18.4.3.7 Subjects in the Device group experienced a significantly greater improvement in NYHA Functional Class at 12 months compared to those in the Control group (Class I or II: 72.2% vs. 49.6%), as shown in Figure 19A, where subjects who died prior to 12 months were imputed as having NYHA Class IV. The NYHA Functional Class (unimputed) through 24 months is shown in Figure 19B.

**Figure 19A: NYHA Functional Class at 12 Months (ITT Population)**

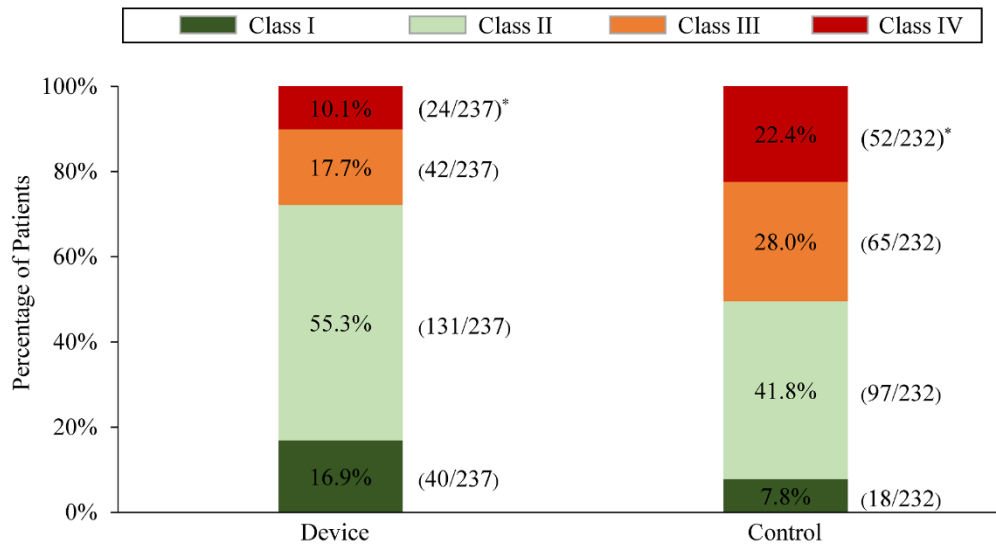

\*Subjects died of HF prior to 12 months were imputed as having NYHA Class IV (Device group: 18; Control group: 41)


**Figure 19B: NYHA Functional Class through 24 Months (ITT Population)**

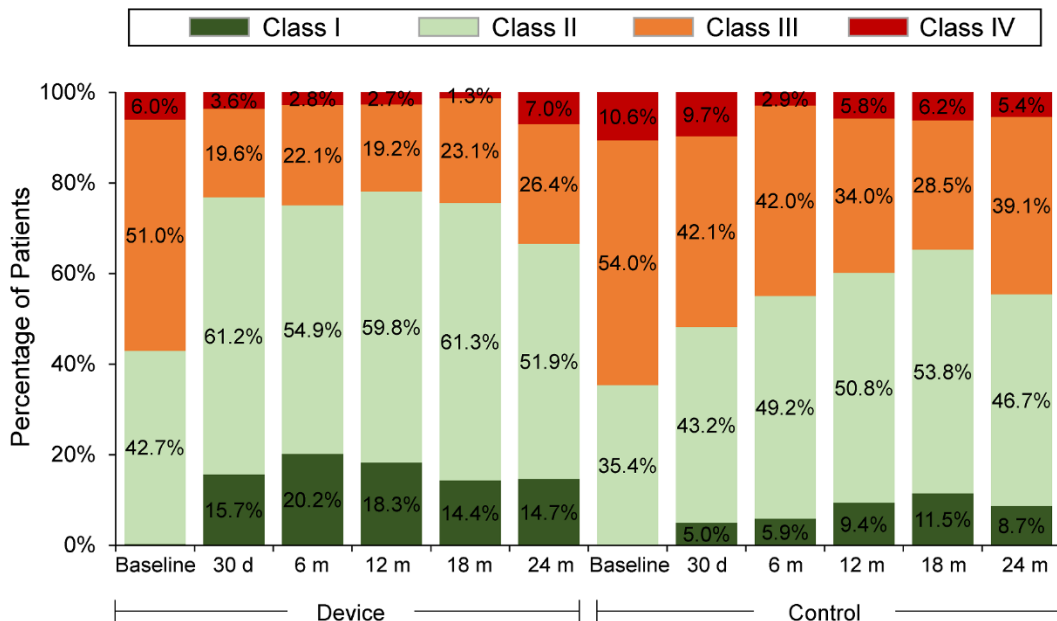

18.4.3.8 Subjects in the Device group experienced significantly greater reduction in LVEDV between baseline and 12 months compared to those in the Control group (-3.71 mL vs. 17.06 mL), as shown in Figure 20. However, while per protocol this endpoint passes, this finding appears to be primarily related to pre-specified imputation of LVEDV values for subjects who died of HF prior to completing the 12-month follow-up where these subjects were assigned the worst LVEDV change between baseline and 12 months observed for any subject in the analysis (126 mL).


Because subjects in the Control group had a numerically higher (41 vs. 18) incidence of HF-related mortality than those in the Device group and the worst change in LVEDV was extreme, calculations for the LVEDV change from baseline in the Control group patients could be skewed mathematically to the larger end. It should be noted that neither clinically nor statistically significant difference in LVEDV change from baseline to 12 months was observed between the Device and Control groups based on un-imputed unpaired and paired analyses or based on a responder analysis.

**Figure 20: LVEDV Change from Baseline to 12 Months (ITT Population)**

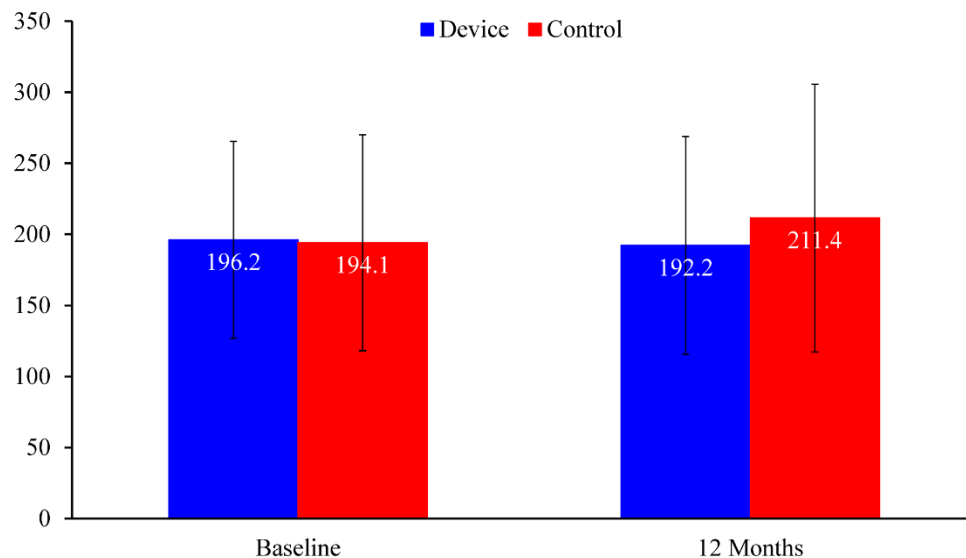


18.4.3.9 Subjects in the Device group experienced significantly lower all-cause mortality at 24 months compared to those in the Control group (Kaplan-Meier estimate: 29.1% vs. 46.1%), as shown in Figure 21. The number needed to treat (NNT) to save one life within 24 months was 5.9 (95% CI: [3.9, 11.7]).

**Figure 21: Kaplan-Meier Curve of All-Cause Mortality through 24 Months (ITT Population)**

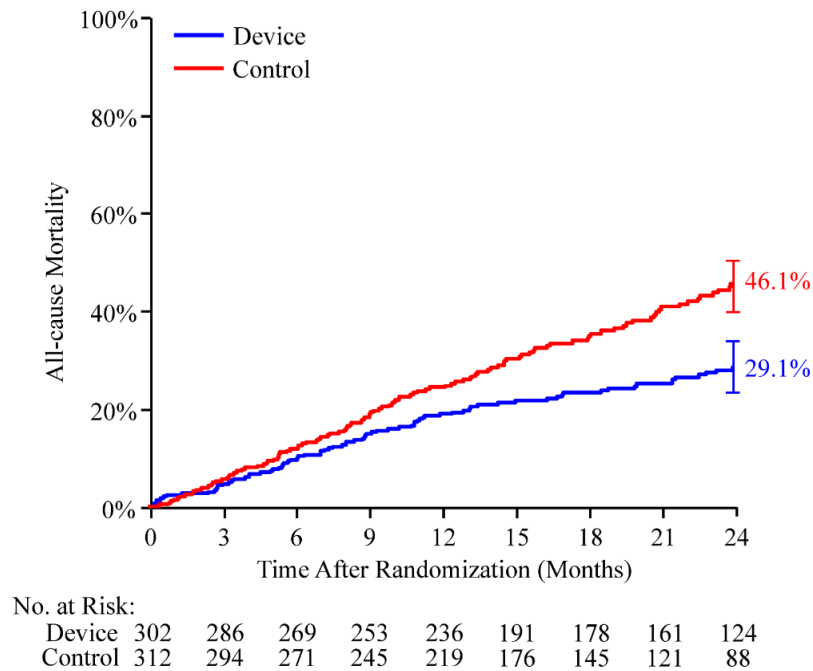

18.4.3.10 The rate of freedom from all-cause mortality, stroke, MI, or non-elective cardiovascular surgery for device-related complications at 30 days was 96.9%, with a lower 95% confidence limit of 94.7%, which met the prespecified performance goal of 80%, as shown in Figure 22.

**Figure 22: Freedom from All-Cause Mortality, Stroke, MI or Non-Elective Cardiovascular Surgery for Device-Related Complications at 30 Days (SA Population)**

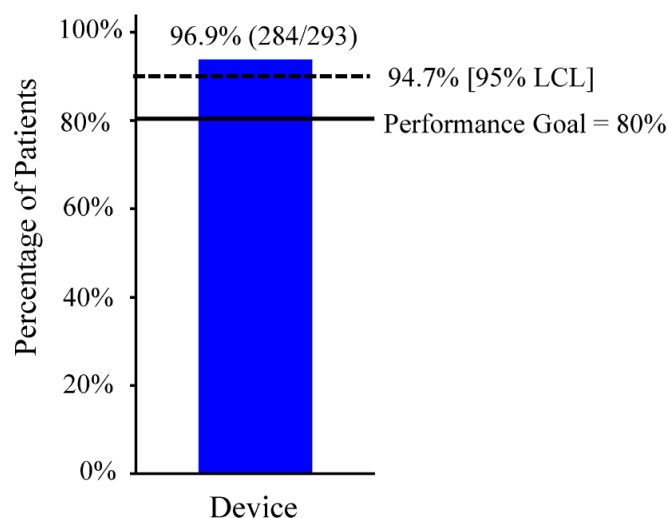

#### 18.4.4 Long-term Endpoints up to 60 Months:

The 60-month COAPT results were evaluated on an As-Treated (AT) population. Subjects in the Control group were allowed to cross over to the Device group after primary endpoint analysis at 24-months. The crossover patients are included in the control group for this analysis.

18.4.4.1 The rate of freedom from device-related complications through 60 months for subjects in the Device group was 89.2%. All occurrences of SLDA, device embolization and device-related complications requiring non-elective CV surgery were within the first 30 days. A breakdown of the composite primary safety endpoint events is presented in Table 28.


**Table 28: Outcomes of Long-Term Safety Endpoint Components (SA Population)**

Event	Summary Statistics* (N = 293)
Device-related complications at 60 months	23 (10.8%)
-- Single leaflet device attachment	2 (0.7%)
-- Device embolization	1 (0.3%)
-- Endocarditis requiring surgery	0 (0.0%)
-- Mitral stenosis requiring surgery	0 (0.0%)
-- LVAD implant	13 (6.5%)
-- Heart transplant	9 (4.7%)
-- Any device-related complication requiring non-elective cardiovascular surgery	1 (0.3%)

\*# events (Kaplan-Meier rate)


18.4.4.2 At 60-months, the Kaplan-Meier estimate of all-cause mortality event rate was 57.3% (freedom from event rate of 42.7%) for the Device group, and 67.2% (freedom from event rate of 32.8%) for the Control group, as shown in Figure 23.

**Figure 23: All-Cause Mortality through 60 Months**

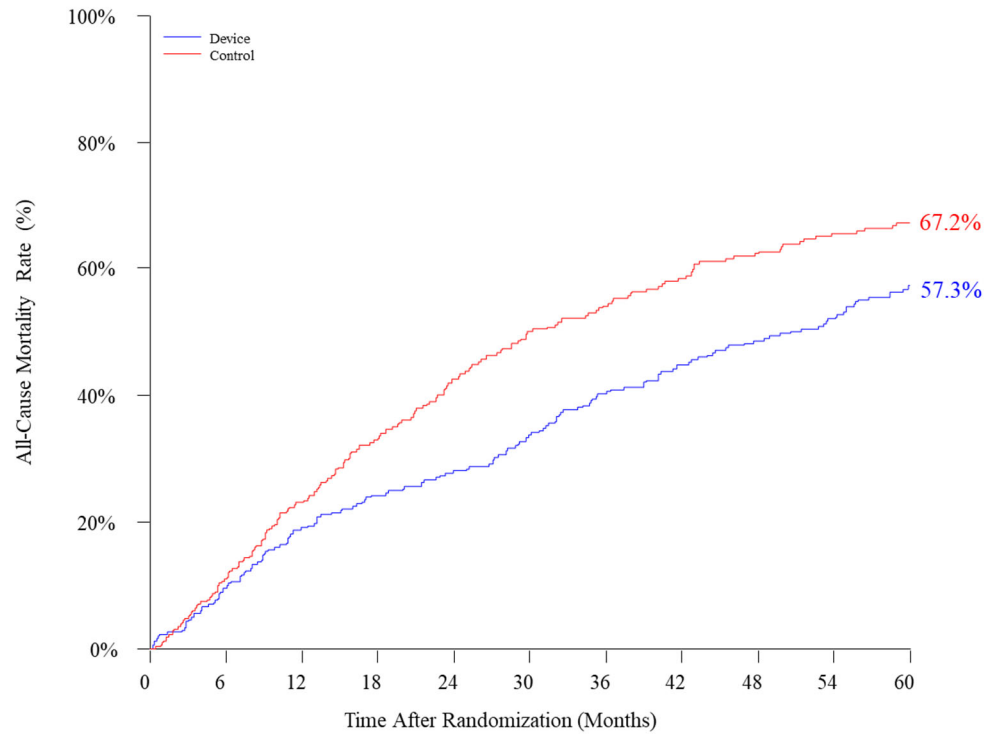


18.4.4.3 At 60-months, the Kaplan-Meier estimate of first heart failure hospitalization event rate was 61.2% (freedom from event rate of 38.8%) for the Device group, compared to 82.4% (freedom from event rate of 17.6%) for the Control group, as shown in Figure 24.

**Figure 24: First Heart Failure Hospitalization through 60 months**

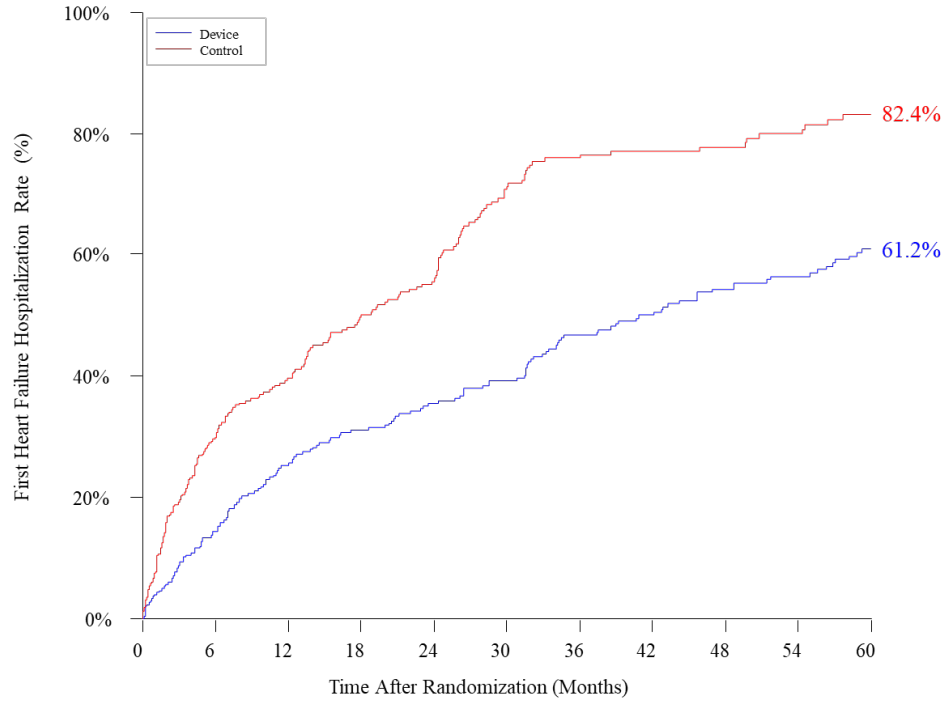

18.4.4.4 At 60 months, 94.7% of subjects in the Device Group achieved  $\leq 2+$  MR compared to 91.3% of subjects in the Control Group, as shown in Figure 25.


**Figure 25: Mitral Regurgitation Severity Through 60 Months**

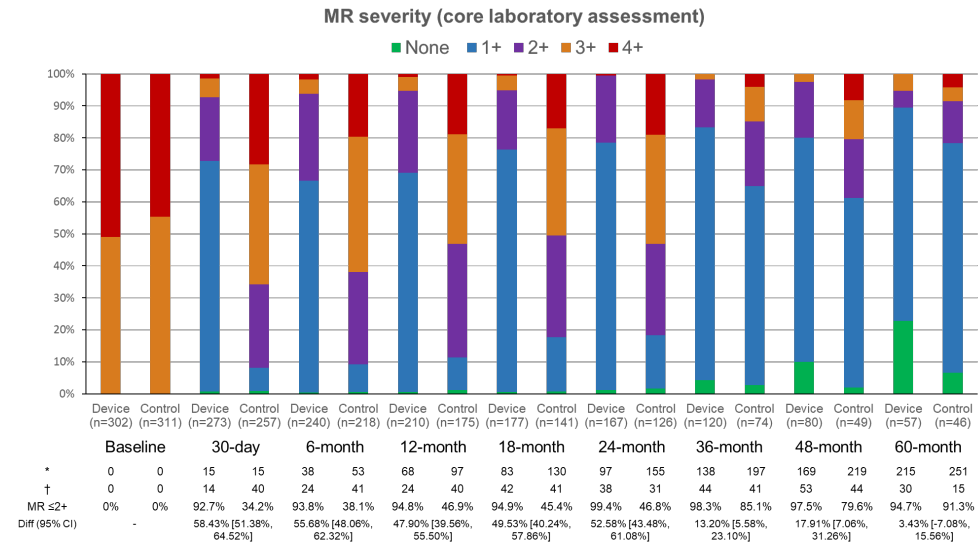



18.4.4.5 At 60 months, subjects in the Device Group were more commonly in NYHA Class I or II compared to those in the Control Group (24.1% vs 15.7%), as shown in Figure 26.

**Figure 26: New York Heart Association Functional Class Through 60 Months**

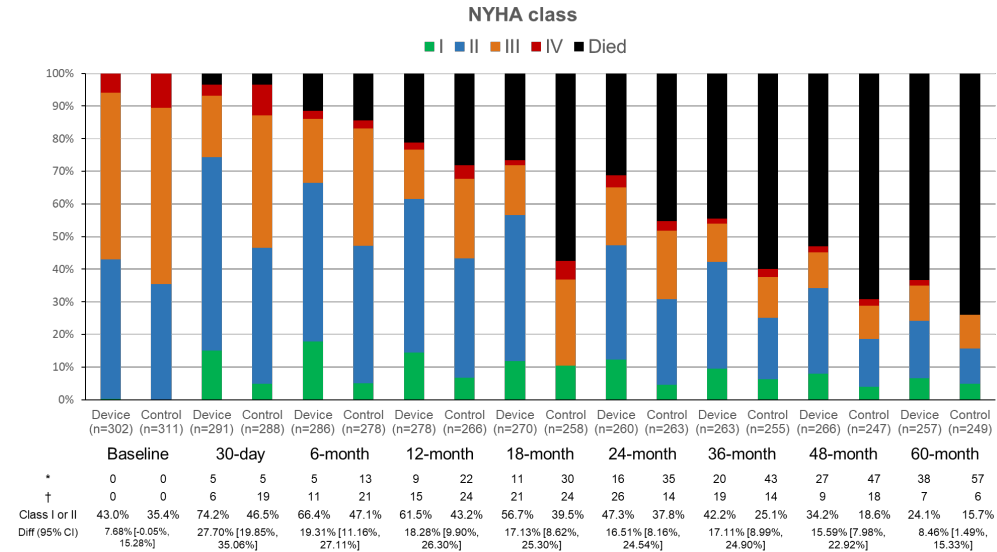

## 18.4.5 Adverse Events

The adverse events that occurred in the trial through 24 months are presented in Table 29.

**Table 29: CEC-Adjudicated Adverse Events through 24 Months (SA Population)**

Events	0-30 Days		0-12 Months		0-24 Months	
	Device	Control	Device	Control	Device	Control
All-cause mortality*	2.3% (7)	1.0% (3)	19.1% (57)	23.2% (70)	29.1% (80)	46.1% (121)
Cardiovascular	2.3% (7)	0.6% (2)	13.8% (40)	19.4% (57)	23.2% (60)	37.0% (93)
Heart failure	0.7% (2)	0.6% (2)	6.2% (17)	13.8% (39)	12.0% (28)	25.9% (61)
Stroke	0.7% (2)	0.0% (0)	2.9% (8)	2.9% (8)	4.4% (11)	5.1% (11)
Transient ischemic attack	0.0% (0)	0.0% (0)	1.1% (3)	1.1% (3)	1.1% (3)	1.1% (3)
Endocarditis requiring surgery	0.0% (0)	0.0% (0)	0.0% (0)	0.0% (0)	0.0% (0)	0.0% (0)
ECL confirmed mitral stenosis requiring surgery	0.0% (0)	0.0% (0)	0.0% (0)	0.0% (0)	0.0% (0)	0.0% (0)
LVAD implant	0.0% (0)	1.0% (3)	0.8% (2)	3.9% (11)	3.0% (6)	7.1% (16)
Heart transplant	0.0% (0)	0.0% (0)	0.8% (2)	2.2% (6)	1.4% (3)	3.6% (8)
Myocardial Infarction†	0.3% (1)	0.0% (0)	NA	NA	NA	NA

Events	0-30 Days		0-12 Months		0-24 Months	
	Device	Control	Device	Control	Device	Control
Major bleeding <sup>†</sup>	5.0% (15)	1.0% (3)	NA	NA	NA	NA
Iatrogenic ASD requiring intervention	0.7% (2)	NA	0.7% (2)	NA	0.7% (2)	NA
Device-related complications requiring non-elective CV surgery	0.3% (1)	NA	0.3% (1)	NA	0.3% (1)	NA

\*Include adjudicated death events and deaths from the national death registry (for subjects who were lost to follow-up or withdrew from the COAPT study).

<sup>†</sup>Events were adjudicated up to 30 days post treatment visit.

Note: (1) Kaplan-Meier rate (# patients with events). Include only each patient's first occurrence of each event. (2) The follow-up duration was calculated from the randomization date. (3) ECL: Echocardiography Core Laboratory; LVAD: Left Ventricular Assists Device; ASD: Atrial Septal Defect; CV: Cardiovascular.

#### 18.4.6 Subgroup Analyses

##### Pre-specified Analyses:

The primary safety and primary effectiveness endpoints were examined across the following 4 subgroups:

- Sex (male vs. female)
- Etiology of cardiomyopathy (ischemic vs. non-ischemic)
- LVEF (> 40% vs. ≤ 40%)
- Extreme surgical risk status (yes vs. no, as determined by the Central Eligibility Committee)

There was no clinically significant difference among the subgroups for the primary safety outcome, and there were no clinically significant interaction effects between treatment and subgroups for the primary effectiveness outcome.

**Post hoc Analyses:**


In light of publication of another study of the MitraClip device in literature, entitled “Mitra-FR Trial” (see reference [4]), a comparison of the baseline characteristics of the subjects enrolled in the COAPT Trial and Mitra-FR Trial was performed to further define the proper patient population for the MitraClip SMR indication. The comparison suggested there were some differences between the two trials as shown in Table 30.

**Table 30: Comparison in EROA and LVEDVi between Mitra-FR and COAPT**

Baseline Characteristics	Mitra-FR	COAPT
EROA (mean $\pm$ SD; cm <sup>2</sup> )	0.31 $\pm$ 0.11	0.41 $\pm$ 0.15
LVEDVi (mean $\pm$ SD; mL/m <sup>2</sup> )	135 $\pm$ 35	101 $\pm$ 34


To explore whether there was any correlation between the clinical outcomes and the baseline EROA and LVEDVi, a *post hoc* subgroup analysis was conducted on the COAPT dataset, by comparing the composite rate of all-cause mortality or HF hospitalization between subjects with an EROA  $\leq$  0.30 cm<sup>2</sup> and an LVEDVi > 96 mL/m<sup>2</sup> and those with an EROA > 0.30 cm<sup>2</sup> or an LVEDVi  $\leq$  96 mL/m<sup>2</sup>, where 0.30 was the lower bound of the EROA defining, along with other parameters, Grade III (or 3+) MR as per the 2017 ASE Recommendation for Noninvasive Evaluation of Native Valvular Regurgitation and 96 was the median LVEDVi value in the COAPT Trial.<sup>5</sup> A total of 22 subjects in the Device group and 34 subjects in the Control group had an EROA  $\leq$  0.30 cm<sup>2</sup> and an LVEDVi > 96 mL/m<sup>2</sup>. The results of the subgroups analysis are shown in Figure 27. COAPT subjects with relatively less severe MR and larger left ventricles (Fig 27a) did not show a clinically meaningful benefit of the Device for all-cause mortality or HF hospitalization at the 12-month timepoint. For the remaining COAPT subjects (those with an EROA > 0.30 cm<sup>2</sup> or an LVEDVi  $\leq$  96 mL/m<sup>2</sup>; Figure 27b), the difference in all-cause mortality or HF hospitalization seen in the overall population was maintained.

**Figure 27: Subgroup Analysis Stratified by EROA and LVEDVi**

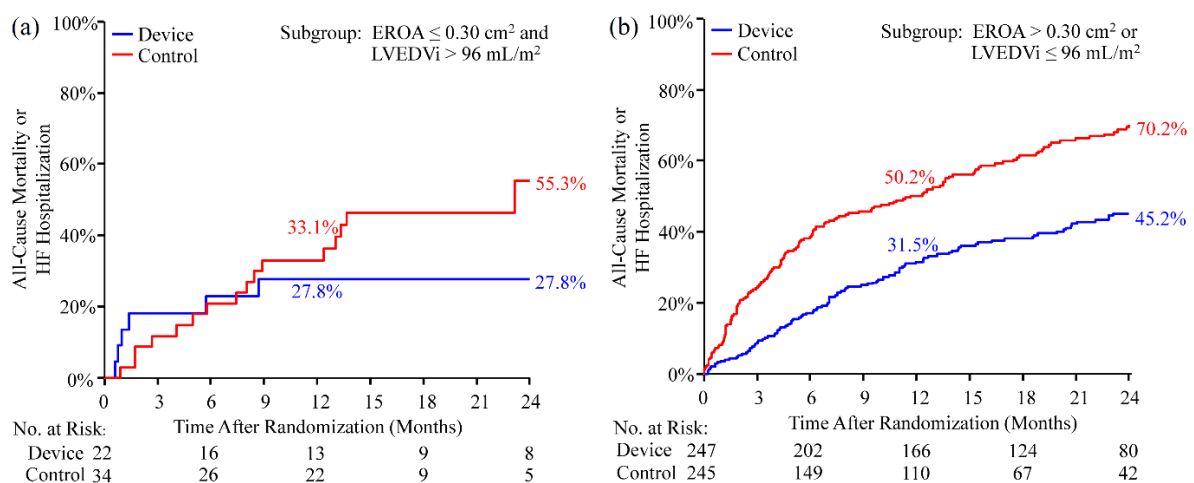


Despite the absence of benefit of reduced all-cause mortality or HF hospitalization in the subgroup with an EROA  $\leq 0.30$  cm<sup>2</sup> and an LVEDVi  $> 96$  mL/m<sup>2</sup>, clinically meaningful improvements in the overall 6MWD (as shown in Figure 28; 11 subjects in the Device group and 26 subjects in the Control group had 6MWD values) and KCCQ (as shown in Figure 29; 15 subjects in the Device group and 27 subjects in the Control group had KCCQ values) compared to baseline were observed in Device group patients, an effect not observed in the same sub-population of the Control group. However, because of the nature of the *post hoc* subgroup analysis and the small sample size in the subgroup with an EROA  $\leq 0.30$  cm<sup>2</sup> and an LVEDVi  $> 96$  mL/m<sup>2</sup>, no statistical or clinical intra-group inferences can be made.

**Figure 28: 6MWD for Subjects with an EROA  $\leq 0.30$  cm<sup>2</sup> and an LVEDVi  $> 96$  mL/m<sup>2</sup>**

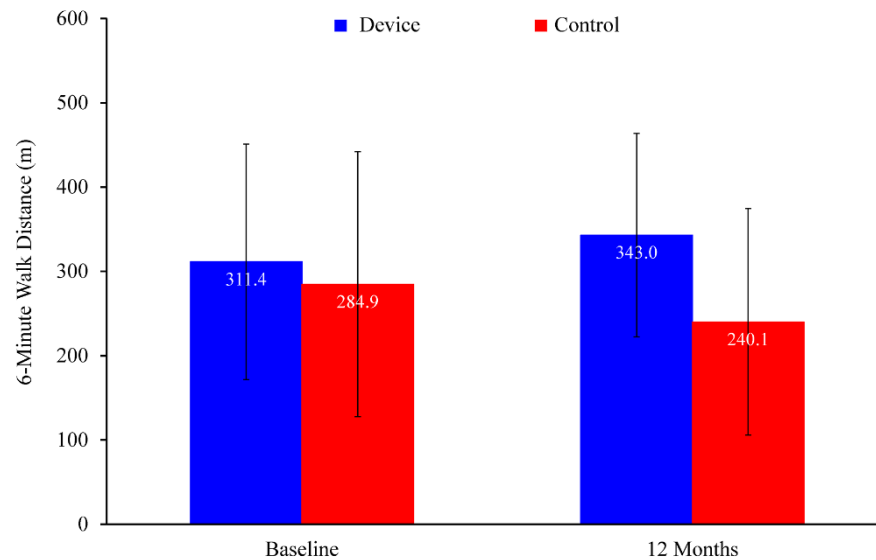

**Figure 29: KCCQ Score for Subjects with an EROA  $\leq 0.30$  cm<sup>2</sup> and an LVEDVi  $> 96$  mL/m<sup>2</sup>**

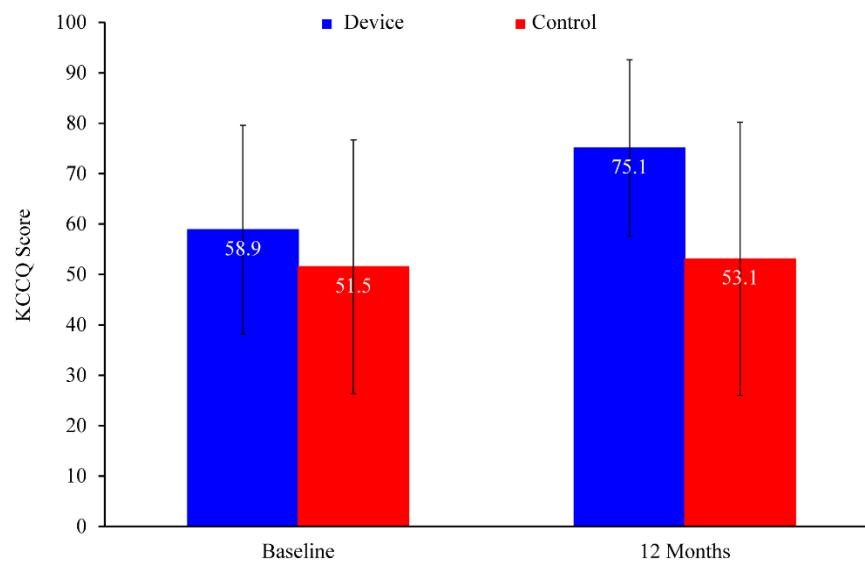

#### 18.4.7 Procedural Data

The procedural data of the Device group are summarized in Table 31.

**Table 31: Procedural Data Summary for Device Subjects – AT Population**

Procedure Data	Device (N=294)
MitraClip Procedure Attempted	100.0%
Implant Rate	98.0%
Number of Clips Implanted	
0 Clip	2.0%
1 Clip	36.4%
2 Clips	53.1%
3 Clips	8.2%
4 Clips	0.3%
Total Number of Clips Implanted	495
Total Procedure Time (min) Mean $\pm$ SD (n) Median (Q1, Q3)	163.0 $\pm$ 117.5 (294) 146.5 (108.0, 199.0)
Device Procedure Time (min) Mean $\pm$ SD (n) Median (Q1, Q3)	118.8 $\pm$ 63.3 (283) 106.0 (73.0, 148.0)
Device Time (min) Mean $\pm$ SD (n) Median (Q1, Q3)	82.6 $\pm$ 80.6 (288) 65.5 (40.0, 100.0)
Fluoroscopy Duration (min) Mean $\pm$ SD (n) Median (Q1, Q3)	33.91 $\pm$ 23.15 (285) 29.50 (18.60, 43.00)

---

#### 18.4.7 References

- [1] Rogers JK, Pocock SJ, McMurray JJV, Granger CB, Michelson EL, Östergren J, et al. Analysing recurrent hospitalizations in heart failure: a review of statistical methodology with application to charm-preserved. *European Journal of Heart Failure* 2014; 16:33–40.
- [2] Rogers JK, Jhund PS, Perez A, Böhm M, Cleland JG, Gullestad L, et al. Effect of rosuvastatin on repeat heart failure hospitalizations: the CORONA trial (controlled rosuvastatin multinational trial in heart failure). *Journal of the American College of Cardiology Heart Failure* 2014; 2:289–297.
- [3] Rogers JK, Yaroshinsky A, Pocock SJ, Stokard D, Pogoda Janice. Analysis of recurrent events with an associated informative dropout time: Application of the joint frailty model. *Statistics in Medicine* 2016; 35:2195-205.
- [4] Obadia JF, Messika-Zeitoun D, Leurent G, Iung B, Bonnet G, Piriou N, et al. Percutaneous Repair or Medical Treatment for Secondary Mitral Regurgitation. *New England Journal of Medicine* 2018; 379:2297-2306.
- [5] Stone GW, Lindenfeld J, Abraham WT, Kar S, Lim DS, Mishell JM, et al. Transcatheter Mitral-Valve Repair in Patients with Heart Failure. *New England Journal of Medicine* 2018;379:2307-2318.
- [6] Zoghbi WA, Adams D, Bonow RO, Enriquez-Sarano M, Foster E, Grayburn PA, et al. Recommendations for Noninvasive Evaluation of Native Valvular Regurgitation: A Report from the American Society of Echocardiography Developed in Collaboration with the Society for Cardiovascular Magnetic Resonance. *Journal of the American Society of Echocardiography* 2017;30:303-371.

---

## 19.0 POST-APPROVAL STUDY – MITRACLIP ANALYSIS COHORT (MAC)

### 19.1 Study Objectives

The primary safety objective of the post approval study, MitraClip Analysis Cohort (MAC) is to demonstrate that freedom from a composite of death and device-related complications at 30 days post-procedure is greater than a performance goal of 80%. The primary effectiveness objective of the post approval study, MitraClip Analysis Cohort (MAC) is to confirm that patients achieve improvement in distance walked on the 6 Minute Walk Test (6MWT) at 12 months over baseline.

### 19.2 Study Design

The Mitral Module of the National Transcatheter Valve Therapy Registry (TVT Registry) housed jointly by the American College of Cardiology Foundation (ACCF) and Society of Thoracic Surgeons (STS) is a registry to track patient safety and real-world outcomes related to transcatheter MV procedures. The Comprehensive, Linked Registry-Based Surveillance MitraClip Registry is a prospective, single-arm, multi-center, observational study of a minimum of 2,000 MitraClip patients consecutively entered into the TVT Registry.

The MitraClip Analysis Cohort (MAC) is a subset of patients entered into the MitraClip Registry who fulfill the following criteria:

- Significant symptomatic MR ( $\geq 3+$ ) due to primary abnormality of the mitral apparatus, otherwise known as Primary (or degenerative) Mitral Regurgitation.
- Prohibitive risk for mitral valve surgery (MV) as determined by a heart team, which includes a cardiologist experienced in MV disease and a cardiac surgeon experienced in MV surgery, and existing comorbidities that would not preclude the expected benefit from reduction of the MR.
- Perform the 6MWT at baseline or unable to perform the 6MWT due to cardiac reasons.

Active follow-up of patients was performed through 12 months with scheduled follow-up visit at discharge 30 days, and 12 months.

### 19.3 Study Population

The data source for MAC consists of a subset of data from TVT registry from sites/centers of excellence and experienced with performing and collecting 6MWT data and KCCQ data (i.e., sites/centers with at least 75% completion rate in terms of collecting both 6MWT and KCCQ data will be included in MAC). The patient demographic in the MAC dataset are as follows: mean age is  $80.6 \pm 9.8$  years, 55.2% male/44.8 % female, 83.5% were NYHA Functional Class III/IV, mean STS valve replacement score of  $11.17 \pm 6.85\%$ , and mean STS valve repair score of  $8.74 \pm 6.74\%$ . All patients had  $\geq 3+$  MR and 90.4% patients had only Primary MR etiology. The most common co-morbidities were heart failure within 2 weeks (79.7%), heart failure hospitalization in the past year (57.4%), hypertension (84.8%), and atrial fibrillation/flutter (66.7%).

---

## 19.4 Data Source

A subset of the first 2,000 MitraClip patients consecutively entered into the STS/ACC TVT Registry.

## 19.5 Key Study Endpoints

The primary safety endpoint is a composite of death and device-related complications at 30 days post-procedure. The device-related complications include single leaflet device attachment (SLDA), device embolization, device thrombosis, and other device-related events (e.g., endocarditis, MV stenosis), which result in MV re-intervention, unplanned other cardiac surgery, or unplanned vascular surgery or intervention through 30 days post-procedure.

The primary effectiveness endpoint is change in 6MWT distance at 12 months over baseline.

## 19.6 Total Number of Enrolled Study Sites and Subjects, Follow-Up Rate

- 50 sites were selected (out of 125 sites) who met the minimum 75% completion rate for 6MWT and KCCQ data entry in the STS/ACC TVT Registry
- A total of 554 patients met eligibility criteria for this study between November 04, 2013, and July 29, 2015
- **Clinical Follow up rate (MAC study):**
  - 30-day clinical follow up rate: 94.8%
  - 12-month clinical follow up rate: 78.2%

## 19.7 Study Visits and Length of Follow-Up

Active follow-up of patients were performed through 12 months with scheduled follow-up visit at discharge, 30 days, and 12 months.

## 19.8 Results from the MAC Cohort

### 19.8.1 Final Safety Findings (Key Endpoints)

The primary safety endpoint of MAC was met. The freedom from composite of death and device-related complications at 30-days post-procedure rate was 95.5% (529/554), with the lower one-sided 95% confidence limit of 93.8%, which was significantly greater than the pre-specified performance goal of 80% ( $p < 0.0001$ , one-sided exact binomial test). The secondary safety endpoints of MAC were met.

- The estimated freedom from all-cause mortality at 12-months was  $80.2\% \pm 1.85\%$ , with the lower limit of one-sided 95% confidence interval of 77.2%, which is greater than the pre-specified performance goal of 66% ( $P < 0.0001$ ).
- Freedom from device-related complications at 12 months was  $98.9\% \pm 0.45\%$ , with the lower limit of one-sided 95% confidence interval of 98.2%, which is greater than the pre-specified performance goal of 90% ( $p < 0.0001$ ).

### 19.8.2 Final Effectiveness Findings (Key Endpoints)

The primary effectiveness endpoint of MAC was also met. The 6MWT distance increased significantly at 12 months over baseline with an average of  $31.6 \pm 115.2\text{m}$  ( $n=184$ , one-sided  $p=0.0001$ , lower 95% CI 17.5m).

The secondary effectiveness endpoints of MAC were met.

- MR Severity at 12 months - Most (86.7%) of the patients went from  $\geq 3+$  MR at baseline to  $\leq 2+$  MR at 12-month follow-up ( $p < 0.0001$ ).



- Change in Left Ventricular End Diastolic Volume (LVEDV) at 12 months Over Baseline – LVEDV showed a mean reduction of  $10.9 \pm 48.9$  ml ( $p=0.0111$ ).
- Change in Left Ventricular Internal Diastolic Dimension (LVIDd) at 12 months Over Baseline – LVIDd showed a mean reduction of  $0.22 \pm 0.72$  cm ( $p<0.0001$ ).
- Change in KCCQ at 12 months over baseline – The KCCQ overall summary score showed a mean improvement of  $30.4 \pm 25.8$  points ( $p<0.0001$ ).
- NYHA Functional Class (FC) Improvement at 12 months – At 12 months post-procedure, 82.6% were NYHA FC I/II, as compared to 20.5% NYHA FC I/II at baseline ( $p<0.0001$ ).

### **19.8.3 Study Strength and Weakness**

#### **Strengths**

- Captured real-world commercial usage and experience with the MitraClip device
- Large sample size allow detection for low frequency safety events

#### **Weaknesses**

- Data entry into the TVT Registry remains voluntary and thus subject to reporting bias
- Data adjudication was limited to specific event categories
- Echocardiographic variables are site reported, not core laboratory adjudicated

## **20.0 REAL-WORLD USE SURVEILLANCE STUDIES – EXPAND AND EXPAND G4**

### **20.1 EXPAND**

#### **20.1.1 Study Objective**

To confirm the safety and performance of the next generation MitraClip™ NTR and MitraClip™ XTR Systems in a contemporary, real-world setting.

#### **20.1.2 Study Design**

The MitraClip EXPAND study was a Prospective, Multi-Center, Single-Arm, International, Post Market, Observational Study designed to collect real-world data on the use of the next generation MitraClip NTR and MitraClip XTR Systems. Approximately 1,000 commercial patients from the EU or US were planned be included in the analysis for *The MitraClip EXPAND Study*.

#### **20.1.3 Study Population**

Consented subjects treated commercially in EU or US with the MitraClip NTR and MitraClip XTR Systems were considered. All subjects were required to meet the below criteria before being included in the MitraClip XTR Study.

##### **Inclusion Criteria**

1. Subjects scheduled to receive the MitraClip per the current approved indications for use
2. Subjects who give consent for their participation
3. Subjects with Symptomatic MR ( $\geq 3+$ )

##### **Exclusion Criteria**

1. Subjects participating in another clinical study that could impact the follow-up or results of this study.

---

#### 20.1.4 Data Source

All data collected was monitored by qualified Abbott personnel. Major adverse events were adjudicated by an independent clinical events committee. Leaflet adverse events were adjudicated by an independent physician committee. Echocardiographic measurements at baseline, discharge, 30-days, and 1-year were assessed by an independent echocardiographic core lab. Results presented here reflect data collected through the final data cutoff of August 5, 2020.

#### 20.1.5 Key Study Endpoints

- **Safety**

Safety was assessed via the evaluation of Major Adverse Events (MAE) at 30 days defined as: *A composite of all-cause Death, Myocardial Infarction, Stroke, or non-elective Cardiovascular (CV) surgery for device related complications (CEC adjudicated)*. The assessment of safety included all occurrences through 30 days post procedure.

- **Performance**

Performance was evaluated as MR Reduction to  $\leq 2+$  at 30-days for this study.

- **Additional Key Endpoints**

- Acute Procedural Success (APS) defined as successful implantation of the MitraClip device with resulting MR severity of 2+ or less on discharge Echocardiogram (30-day echocardiogram will be used if discharge is unavailable or uninterpretable). Subjects who die or undergo mitral valve surgery before discharge are considered to be an APS failure
- MR Reduction to  $\leq 1+$  at Discharge, 1, 12 months
- Functional Improvement Measures (Discharge, 1, 12 months)
  - New York Heart Association (NYHA) functional class improvement
  - Quality of Life (QOL) assessed using the Kansas City Cardiomyopathy Questionnaire (KCCQ)

#### 20.1.6 Total number of Enrolled Study Sites and Subjects, Follow-up Rate

A total of 1045 subjects were enrolled at 58 sites between April 26, 2018 and March 29, 2019. Enrollment in the study was closed on March 29, 2019. Of the 1045 subjects enrolled, 4 subjects did not have any procedure attempted.

#### 20.1.7 Study Visits and Length of Follow-up

Enrolled patients that received a MitraClip implant underwent study follow-up as per standard of care. Follow-up data on clinical assessment was collected for a discharge visit, 30-day visit, 6-month visit (or phone call) and 12-month visit.

#### 20.1.8 Results from the EXPAND Study

##### 20.1.8.1 Thirty-Day Safety Findings

The 30-Day findings in EXPAND demonstrate safe treatment of patients with MR using the MitraClip NTR and XTR System. There was a low incidence of major adverse events and single leaflet device attachment events (Table 32).

**Table 32 Safety Outcomes from EXPAND**

	0-30 Days
Major Adverse Events*	3.9% (40/1036)
All-Cause Death	2.3% (24/1036)
Myocardial Infarction	0.0% (0/1036)
Stroke	0.8% (8/1036)
Non-Elective Cardiovascular Surgery for Device-Related Complications	1.1% (11/1036)
Single Leaflet Device Attachment	1.7% (18/1041)

\*Subjects are only counted once for each type of event. Subjects who terminated from the study prior to the lower limit of the 30-day visit window (Day 16) without any site reported adverse events are excluded from the analysis.

### **20.1.8.2 Thirty-Day Effectiveness Findings**

Treatment with the MitraClip NTR and XTR System in EXPAND was effective (Table 33). There was high acute procedural success rate achieved in 95.8% (985/1028) of subjects with a core lab adjudication of MR  $\leq 2+$  at discharge or 30-days without death or reintervention. There was significant MR reduction 88.8% of subjects having an MR  $\leq 1+$  at 30-days compared to 8.9% with MR  $\leq 1+$  at baseline. Subjects also experienced 30-day symptomatic improvements with 80.1% of subjects in NYHA Class I/II compared to 21.4% at baseline and an increase in KCCQ-OS score of  $19.3 \pm 25.0$  points from baseline to 30 days.

**Table 33 Effectiveness Outcomes from EXPAND**

	30 Days
MR $\leq 2+$	97.8% (845/864)
MR $\leq 1+$	88.8% (767/864)
NYHA Class I/II	80.1% (691/863)
Change in KCCQ-OS Score from Baseline	$19.3 \pm 25.0$ (853)

## **20.2 EXPAND G4**

### **20.2.1 Study Objective**

To evaluate the safety and performance of the MitraClip G4 System in a contemporary, post-market setting.

### **20.2.2 Study Design**

This is a Prospective, Multi-Center, Single-Arm, International, Post Market, Observational Study designed to collect real-world data on the use of the next generation MitraClip G4 System. The study is designed in two phases:

- **Phase I:** Clinical and echocardiographic outcomes data was collected from approximately 100 subjects enrolled at approximately 20 sites in the United States with the goal of evaluating the early safety and procedural outcomes associated with the MitraClip G4. The follow up duration was 30 days. Phase I was completed August 2020.
- **Phase II:** Up to 1,000 post-market, consented patients will be treated with a MitraClip G4 device according to local guidelines and IFU from the Europe, the United States, Canada and Japan will be included in the analysis for The MitraClip EXPAND G4 Study. Follow-up echocardiograms will be collected at Discharge, 30 days and 1 year and 5 year at post-procedure visits. Additional clinical follow-up visits will be at 6 months (phone call), 2, 3, 4 year (office visits). Cardiovascular adverse events will

---

be reported through 5 years to confirm safety of the MitraClip G4 System.

### 20.2.3 Study Population

Consented subjects treated commercially in the US with the MitraClip G4 System were considered. All subjects were required to meet the below criteria before being included in the MitraClip EXPAND G4 Study.

#### Inclusion Criteria

1. Subjects scheduled to receive the MitraClip per the current approved indications for use
2. Subjects who give consent for their participation

#### Exclusion Criteria

1. Subjects participating in another clinical study that could impact the follow-up or results of this study.

### 20.2.4 Data Source

All data collected was monitored by qualified Abbott personnel. Major adverse events were assessed by the site and reviewed by qualified Abbott safety personnel. Leaflet adverse events and echocardiographic measurements at baseline, discharge, 30-days, 1-year, and 5-years were assessed by an independent echocardiographic core lab. Data collected here reflect the 30-day cutoff of May 9, 2022.

### 20.2.5 Key Study Endpoints

- **Safety**

Safety was assessed via the evaluation of Major Adverse Events (MAE) at 30 days defined as: *A composite of all-cause Death, Myocardial Infarction, Stroke, or non-elective Cardiovascular (CV) surgery for device related complications.* The assessment of safety included all occurrences through 30 days post procedure.

- **Performance**

Performance was evaluated as MR Reduction to  $\leq 2+$  at 30-days for this study.

- **Additional Key Endpoints**

- Acute Procedural Success (APS) defined as successful implantation of the MitraClip device with resulting MR severity of 2+ or less on discharge Echocardiogram (30-day echocardiogram will be used if discharge is unavailable or uninterpretable). Subjects who die or undergo mitral valve surgery before discharge are considered to be an APS failure.
- MR Reduction to  $\leq 1+$  (Discharge, 1, 12 months, and 5 years)
- MR Reduction to  $\leq 2+$  (Discharge, 1, 12 months, and 5 years)
- Functional Improvement Measures (Baseline, Discharge, 12 months, and 2, 3, 4, 5 years)
  - New York Heart Association (NYHA) functional class improvement
  - Quality of Life (QOL) assessed using the Kansas City Cardiomyopathy Questionnaire (KCCQ)

### 20.2.6 Total number of Enrolled Study Sites and Subjects, Follow-up Rate

From January 21, 2020 to February 4, 2022, a total of 1164 subjects at 60 clinical sites were enrolled into the MitraClip EXPAND G4 Study: 101 from Phase I and 1063 from Phase II.

### 20.2.7 Study Visits and Length of Follow-up

Follow-up data on clinical assessment are submitted at discharge, 30-day visit, 6-month phone call, 1 year, 2 year, 3 year, 4 year and 5 year visits. Echocardiograms are required to be submitted for the baseline, discharge, 30-day and 12-month and 5 year time points.

### 20.2.8 Results from the EXPAND G4 Study

#### 20.2.8.1 Thirty-Day Safety Findings

The 30-Day findings in EXPAND G4 demonstrate safe treatment of patients with MR using the MitraClip G4 System. There was a low incidence of major adverse events and single leaflet device attachment events (Table 34).

**Table 34 Safety Outcomes from EXPAND G4**

	<b>0-30 Days</b>
Major Adverse Events*	2.7% (31/1158)
All-Cause Death	1.3% (15/1158)
Myocardial Infarction	0.2% (2/1158)
Stroke	0.5% (6/1158)
Mitral Valve Surgical Reintervention	0.9% (10/1158)
Single Leaflet Device Attachment	1.0% (12/1158)

\*Subjects are only counted once for each type of event. Subjects who terminated from the study prior to the lower limit of the 30-day visit window (Day 16) without any site reported adverse events are excluded from the analysis.

#### 20.2.8.2 Thirty-Day Effectiveness Findings

Treatment with the MitraClip G4 System in EXPAND G4 was effective (Table 35). There was a high acute procedural success rate achieved in 96.2% (1099/1143) of subjects with a core lab adjudication of MR severity  $\leq 2+$  at discharge or 30-days without death or reintervention. There was significant MR reduction to  $\leq 1+$  at 30-days compared to 5.8% of subjects with MR  $\leq 1+$  at baseline. Subjects also experienced 30-day symptomatic improvements with a greater proportion in NYHA Class I/II compared to 30.8% at baseline and an increase in KCCQ-OS score of  $18.5 \pm 24.7$  points from baseline to 30-days.

**Table 35 Effectiveness Outcomes from EXPAND G4**

	<b>30 Days</b>
MR $\leq 2+$	97.9% (808/825)
MR $\leq 1+$	90.8% (749/825)
NYHA Class I/II	83.3% (862/1035)
Change in KCCQ-OS Score from Baseline	$18.5 \pm 24.7$ (982)

## 21.0 MITRACLIP™ G4 PROCEDURE STEP-BY-STEP INSTRUCTIONS

### 21.1 DEFINITION OF TERMS

*Defined Terms are in italics throughout document.*

TERM	DEFINITION AND RELATED TECHNIQUE
<i>Lock the Clip</i>	<ol style="list-style-type: none"> <li>1. Rotate the Lock Lever outward.</li> <li>2. Fully advance the Lock Lever.</li> <li>3. Rotate the Lock Lever inward to engage the lever.</li> </ol>
<i>Unlock the Clip</i>	<ol style="list-style-type: none"> <li>1. Rotate the Lock Lever outward and then retract the lever until the mark on the lever is fully exposed.</li> <li>2. Rotate the Lock Lever inward to engage the lever.</li> </ol>
<i>Open the Clip Arms</i>	<ol style="list-style-type: none"> <li>1. Confirm the Clip is unlocked.</li> <li>2. Turn the Arm Positioner at least 1/2 turn in the “Close” (clockwise) direction.</li> <li>3. Turn the Arm Positioner in the “Open” (counter-clockwise) direction until the desired <i>Clip Arm Angle</i> is achieved.</li> </ol> <p><b>NOTE 1:</b> If Clip does not open smoothly, retract the Lock Lever farther, then repeat steps 2 – 3.</p> <p><b>NOTE 2:</b> If the Clip Arms fail to open visibly (as observed under fluoroscopic guidance), use the following techniques in the order provided, as needed:</p> <ol style="list-style-type: none"> <li>A. Stop and return <i>Arm Positioner to Neutral</i>. Retract Lock Lever farther, then turn the Arm Positioner farther in the “Close” direction before turning in the “Open” direction. Advance the lock lever just enough so that the mark on the lever is still fully exposed.</li> <li>B. Turn the <i>Arm Positioner to Neutral</i>, then incrementally iterate the amount of Arm Positioner rotation in the “Close” direction followed by rotation in the “Open” direction. Iterate until Clip opens or until it is no longer possible to rotate the Arm Positioner in the “Close” direction. Advance the lock lever just enough so that the mark on the lever is still fully exposed.</li> <li>C. Turn the <i>Arm Positioner to Neutral</i>, iterate the amount of Lock Lever retraction past the mark in 5 mm increments, and rotate the Arm Positioner fully in the “Close” direction, before rotating in the “Open” direction, until Clip opens. Advance the lock lever just enough so that the mark on the lever is still fully exposed.</li> <li>D. Advance the Gripper Lever and repeat NOTE 2, Step C. Retract the Gripper Lever after Clip opens.</li> <li>E. If in the LA and free of tissue, release the DC Fastener, then release the Sleeve curves and repeat NOTE 2, Step C.</li> </ol> <p><b>WARNING:</b> Do release the DC Fastener before releasing Sleeve curves, otherwise it may result in device damage and / or device or component embolization.</p> <ol style="list-style-type: none"> <li>F. If the Clip does not open after performing all steps in NOTE 2, DO NOT use the device.</li> </ol>
<i>Arm Positioner to Neutral</i>	Turn the Arm Positioner in the “Close” or “Open” direction until no resistance to turning is noted.
<i>Invert the Clip Arms</i>	<ol style="list-style-type: none"> <li>1. Confirm the Clip is unlocked.</li> <li>2. Turn the Arm Positioner at least 1/2 turn in the “Close” direction.</li> </ol>

TERM	DEFINITION AND RELATED TECHNIQUE
	<p>3. Turn the Arm Positioner in the “Open” direction until a <i>Clip Arm Angle</i> of 180° is observed under fluoroscopic guidance. Note the orientation of the blue line on the Arm Positioner.</p> <p>4. Continue turning the Arm Positioner in the “Open” direction until the Clip Arms invert, no more than 1 full turn from 180°. DO NOT over-invert the Clip Arms. DO NOT turn Arm Positioner more than 1 full turn past a <i>Clip Arm Angle</i> of 180° or past when resistance is first noted.</p> <p><b>WARNING: Turning the Arm Positioner in the “Open” direction more than 1 full turn past a <i>Clip Arm Angle</i> of 180° or turning past when resistance is first noted may result in device damage which could cause the Clip to become non-functional and lead to embolization, and / or conversion to surgical intervention.</b></p>
<i>Raise the Gripper(s)</i>	<p>(a) Simultaneously- Confirm Gripper Levers are latched. Slowly retract the Gripper Levers until a hard stop is reached (under fluoroscopic and echocardiographic observation).</p> <p style="text-align: center;"><b>or</b></p> <p>(b) Independently- Unlatch the Gripper Levers and slowly retract the desired Gripper Lever until a hard stop is reached (under fluoroscopic and echocardiographic observation). Re-latch the Gripper Levers.</p>
<i>Lower the Gripper(s)</i>	<p>(a) Simultaneously- Confirm the Gripper Levers are latched and fully advance the Gripper Levers.</p> <p style="text-align: center;"><b>or</b></p> <p>(b) Independently- Unlatch the Gripper Levers and fully advance the desired Gripper Lever(s). Re-latch the Gripper Levers.</p>
<i>Identify Gripper Orientation</i>	<p>1. Unlatch the Gripper Levers.</p> <p>2. Advance and retract the Gripper Lever with the tactile marker under imaging (echocardiography) to identify Gripper Lever to the corresponding leaflet.</p> <p>3. Once Gripper(s) are identified, <i>Raise the Gripper(s)</i> until ready for leaflet capture.</p>
<i>Clip Arm Angle</i>	<ul style="list-style-type: none"> <li>Angle between the inner edges of both Clip Arms.</li> <li>All <i>Clip Arm Angles</i> are measured using fluoroscopy with optimal view allowing clear observation of the tip of the Clip and both arms in the same plane so they appear as a “V” (see Figure 6).</li> </ul>
<i>Grasping Arm Angle</i>	A <i>Clip Arm Angle</i> of approximately 120 degrees.
<i>Fully Close the Clip Arms</i>	<p>Turn the Arm Positioner in the “Close” direction until the Clip Arms contact the DC.</p> <ul style="list-style-type: none"> <li>Under direct visualization, the Clip is fully closed when the Clip Covering contacts the DC.</li> <li>Under fluoroscopic observation, the Clip is fully closed when the inner edges of the Clip Arms are parallel.</li> </ul>

TERM	DEFINITION AND RELATED TECHNIQUE
<i>Establish Final Arm Angle</i>	<p>Verification step to confirm that the pre-deployment <i>Clip Arm Angle</i> will reflect the <i>Clip Arm Angle</i> post-deployment.</p> <ol style="list-style-type: none"> <li>1. With the Lock Lever fully advanced, and the <i>Arm Positioner to Neutral</i> (note the orientation of the blue line on the Arm Positioner), turn the Arm Positioner 1 turn in the “Open” direction (confirm blue line has returned to the original orientation). The Clip Arms may open slightly (~5°) and then remain in a stable position.</li> </ol> <p><b>NOTE:</b> If continued opening of the Clip Arms is noted, reconfirm that the Lock Lever is completely advanced. Close the Clip Arms and <i>Establish Final Arm Angle</i>.</p> <p><b>WARNING:</b> <b>DO NOT</b> turn the Arm Positioner more than 1 turn in the “Open” direction from neutral. Failure to stop turning the Arm Positioner at 1 turn in the “Open” direction past neutral may result in Clip opening or device damage which could cause the Clip to become non-functional and lead to embolization and / or conversion to surgical intervention.</p>
<i>Left Atrial Pressure Capability</i>	<p>If using left atrial pressure capability, connect a fluid filled pressure monitoring system (not supplied) to the luer port of the Guide valve housing. Calibrate the pressure monitoring system and obtain measurements via the manufacturer’s instructions for use.</p> <p><b>NOTE:</b> Ensure no bubbles are present after inserting the CDS into the Guide. When measuring left atrial pressure, ensure the Clip is not obstructing the Guide Lumen, the Guide RO Tip Ring is between the RO Alignment Markers of the Sleeve, and the Sleeve is not deflected more than 90 degrees as this may impair Guide lumen patency.</p>

## 22.0 PATIENT PREPARATION

- 22.1 Prepare the patient per institution’s standard practice for transseptal catheterization.
- 22.2 Place support plate under patient’s leg in the region between the area of the upper leg and the knee and place the Lift over the ipsilateral lower extremity prior to draping the patient.
- 22.3 Place the Lift on the Support Plate such that the front edge (i.e., the edge that corresponds with the shorter legs of the Lift) is approximately 80 cm from the patient’s mid sternum.
- 22.4 Adjust the height of the Lift so that the front edge of the Lift is close to the patient’s leg but is not impinging on it. Adjust the back legs to be 2 or 3 notches above the front legs (i.e., the back legs of the Lift are taller than the front legs).

**CAUTION:** Ensure the Lift and Support Plate are covered completely by sterile drape during the procedure. Use towels as necessary to minimize direct contact between the patient and all surfaces of both the Lift and Support Plate.

- 22.5 Prepare the patient for invasive hemodynamic monitoring.

## 23.0 MITRACLIP™ G4 SYSTEM PREPARATION BEFORE USE

**WARNING:** **DO NOT** use the MitraClip™ G4 System after the “Use By” date stated on the package label, and never reuse or re-sterilize the system. Use of expired, reused, or re-sterilized devices may result in infection.



---

**WARNING:** Always inspect the MitraClip™ G4 System and its packaging to verify no damage has occurred as a result of shipping and handling and that the sterile barrier has not been compromised. **DO NOT** use the device if damage is detected. Use of product with a compromised sterile barrier may result in infection. Use of damaged product may result in patient injury.

- DO NOT remove the protective cover placed over the Clip.

**WARNING:** **DO NOT** handle the Clip directly; leave it in the protective cover to avoid potential contamination. Removal of the protective cover may result in infection. Removal of the protective cover may result in damaged product which may result in patient injury.

- The preparation is most easily accomplished with the aid of an assistant.

### **23.1 Steerable Guide Catheter Preparation**

**WARNING:** All lumens contain air when shipped. Use proper de-airing techniques before and during use to minimize the risk of air embolism.

23.1.1 Carefully remove the white Guide tip shape retainer and transparent protective tubing from the Guide tip.

23.1.2 Inspect Steerable Guide Catheter and Dilator to verify they are undamaged.

**WARNING:** **DO NOT** use if damage is detected. Use of damaged product may result in air embolism, vascular and / or cardiac injury.

23.1.3 Remove the sterile package containing Fasteners and Silicone Pad from the Steerable Guide Catheter tray.

23.1.4 Fill a basin with heparinized saline.

23.1.5 Flush and de-air the Guide and Dilator with heparinized saline:

23.1.5.1 Connect 3-way stopcocks to the Guide and Dilator flush ports, de-air the Dilator, and close the stopcock and the Rotating Hemostatic Valve.

23.1.5.2 Hydrate 5-10 cm of the distal end of the Dilator with heparinized saline and insert the Dilator approximately 10 cm into Guide then remove.

23.1.5.3 Connect high pressure tubing and a 50–60 cc syringe filled with heparinized saline to the Guide flush port.

23.1.5.4 De-air the Guide.

23.1.5.4.1 With the tip raised, displace all air from the Guide while tapping along the length of the catheter shaft.

23.1.5.4.2 Cover the Guide tip with finger once heparinized saline exits the Guide.

---

23.1.5.4.3 Close the Guide stopcock.

23.1.6 Submerge the Guide tip in the basin of heparinized saline.

23.1.7 While the Guide tip is submerged in the basin of heparinized saline, remove finger from Guide tip and check the Guide valve for leaks by raising the handle to a vertical position for a minimum of 30 seconds.

23.1.8 Hydrate 5-10 cm of the distal end of the Dilator with heparinized saline.

23.1.9 Cover the Guide tip with finger and insert the Dilator into the Guide while Guide tip remains submerged in the basin of heparinized saline.

23.1.9.1 While advancing the Dilator, continually watch for air in the Guide Hemostasis Valve housing. If needed, remove finger from Guide tip and aspirate while assuring the Guide tip is submerged.

23.1.10 Remove finger from tip of Guide when the Dilator tip approaches the Guide tip.

23.1.11 Position the Dilator to create a smooth transition.

## **23.2 Stabilizer Preparation**

Assemble the sterilized Stabilizer. Set the Stabilizer aside in a protected sterile environment for later use.

## **23.3 Clip Delivery System Preparation**

23.3.1 Release the DC Fastener and rotate the DC Handle 90 degrees clockwise when the CDS is removed from packaging. Secure the DC Fastener.

23.3.2 Inspect all Clip Delivery System parts, including the Clip, DC shaft, and Sleeve to verify they are undamaged.

**WARNING: DO NOT use the device if damage is detected. Use of damaged product may result in air embolism, Implant or device component embolization, vascular and / or cardiac injury.**

### **Sleeve Preparation**

23.3.3 Connect 3-way stopcocks to the Sleeve flush port and bottom DC flush port.

23.3.4 Remove the cap from the Clip Introducer and place the cap on the top flush port of the DC Handle.

23.3.5 Connect a 3-way stopcock to the Clip Introducer flush port.

23.3.6 Attach a 50-60 cc syringe filled with heparinized saline to the 3-way stopcock on the Clip Introducer, de-air the Clip Introducer, and close the stopcock.

23.3.7 Connect one high pressure tube to each drip line from the pressurized bags with sterile heparinized saline; flush and de-air the lines.

---

23.3.8 Connect one high pressure tube to the 3-way stopcock on the bottom flush port of the DC Handle and one high pressure tube to the 3-way stopcock on the flush port of the Sleeve Handle.

23.3.9 Flush and de-air the Sleeve with heparinized saline.

23.3.9.1 While flushing, release the DC Fastener, retract and advance the DC Handle to remove residual air from the lumen.

**WARNING: DO NOT use excessive force when pulling the DC (Delivery Catheter) Radiopaque Ring against the Sleeve tip, while translating the DC shaft. It may result in device damage including distal tip embolization.**

23.3.10 Secure the DC Fastener with DC Handle fully advanced.

#### **Delivery Catheter Preparation**

**CAUTION: DO NOT handle the Clip directly; leave in the protective cover to avoid potential contamination. Removal of the protective cover may result in damaged product which may result in patient injury.**

23.3.11 Temporarily remove the cap from top flush port of the DC Handle to flush and de-air DC Handle and all lumens of the DC with heparinized saline.

23.3.12 After de-airing the DC Handle chamber, replace the cap to close off top flush port of the DC Handle.

23.3.13 Retract and advance the Lock Lever several times to remove residual air from the lumens.

**WARNING: DO NOT retract the Lock Lever forcefully. It may result in the inability to lock or unlock the Clip. Damage could occur causing the Clip to not unlock or open. Inability to open the Clip may result in valve injury or lead to deployment of the Clip in an unintended location.**

23.3.14 Loosen the Lock Lever Cap to de-air. DO NOT turn Lock Lever Cap more than 1/2 turn in the "Open" direction. After de-airing, tighten the Lock Lever Cap.

23.3.15 With the tip raised and the shaft held taut, displace all air from the DC and sleeve while tapping along the length of the catheter shaft.

23.3.16 Confirm continuous flow from the distal end of the DC and Sleeve.

#### **Delivery Catheter and Clip Inspection**

**WARNING: DO NOT handle the Clip directly, leave in the protective cover to avoid potential contamination. Removal of the protective cover may result in infection. Removal of the protective cover may result in damaged product which may result in patient injury.**

23.3.17 *Raise the Gripper(s), Unlock the Clip, and confirm the Clip is at Grasping Arm Angle.*

---

**NOTE:** If *Clip Arm Angle* is greater than *Grasping Arm Angle*, close the Clip to *Grasping Arm Angle*; if *Clip Arm Angle* is less than *Grasping Arm Angle*, *Open the Clip Arms to Grasping Arm Angle*.

23.3.18 *Lower the Gripper(s)* once to de-air the lumens.

23.3.19 Close the Clip Arms to approximately 60 degrees and *Lock the Clip*.

23.3.20 Close the Clip Arms to approximately 20 degrees and *Establish Final Arm Angle*.

**WARNING:** **DO NOT** turn the Arm Positioner more than 1 turn in the “Open” direction from neutral. Failure to stop turning the Arm Positioner at 1 turn in the “Open” direction past neutral may result in Clip opening or device damage which could cause the Clip to become non-functional and lead to embolization and / or conversion to surgical intervention.

23.3.21 Return the *Arm Positioner to Neutral*.

23.3.22 *Raise the Gripper(s), Unlock the Clip, and Invert the Clip Arms*.

**WARNING:** **DO NOT** retract the Lock Lever forcefully. It may result in the inability to lock or unlock the Clip. Damage could occur causing the Clip to not unlock or open. Inability to open the Clip may result in valve injury or lead to deployment of the Clip in an unintended location.

**WARNING:** Turning the Arm Positioner in the “Open” direction more than 1 full turn past a *Clip Arm Angle* of 180° or turning past when resistance is first noted may result in device damage which could cause the Clip to become non-functional and lead to embolization, and / or conversion to surgical intervention.

23.3.23 *Lock the Clip, Fully Close the Clip Arms, and Lower the Gripper(s)*.

23.3.24 Release the DC Fastener and torque the DC Handle clockwise and counter-clockwise ¼ turn while translating the shaft.

23.3.25 Retract the DC fully against the sleeve and secure the DC Fastener.

**WARNING:** Do not use excessive force when pulling the DC (Delivery Catheter) Radiopaque Ring against the Sleeve tip while translating the DC shaft. It may result in device damage including distal tip embolization.

23.3.26 Without removing the protective cover, carefully slide the Clip Introducer over the Clip and stop when the tip of the Clip is just proximal to the tip of the Clip Introducer.

**WARNING:** **DO NOT** compress the Clip Arms. Compressing the Clip Arms may result in inability to open the Clip. Inability to open the Clip may result in valve injury or lead to deployment of the Clip in an unintended location.

---

**WARNING:** Fully Close the Clip Arms before insertion or retraction into the Clip Introducer. Failure to do so may result in difficulty or inability to advance or retract the Clip, which may result in vascular and / or cardiac injury, air embolism, and / or the need for surgical intervention.

23.3.27 Return the *Arm Positioner* to *Neutral*.

23.3.28 Temporarily discontinue heparinized saline flushes.

23.3.29 Re-start heparinized saline flushes just before the use of CDS.

**WARNING:**

- Heparinized saline flush should be continuous throughout the procedure.
- Discontinuing flush may result in air embolism and / or thrombus formation.
- Ensure flow is visible through the drip chamber and that the tubing is free from kinks and / or obstruction.
- Ensure pressure of 300 mmHg is maintained.

## 24.0 ACCESS TO THE MITRAL VALVE

**NOTE:** This is a suggested sequence for the procedure. Variations may be used based upon patient anatomy.

24.1 Access the LA to accommodate the Guide tip using transvenous, transseptal techniques and equipment.

24.2 Heparinize the patient.

**WARNING:** Failure to administer heparin once transseptal access has been achieved may result in thrombus formation.

24.3 Carefully place an exchange length supportive guidewire in the left upper pulmonary vein or LA. Dilate the subcutaneous tissue and femoral vein to accommodate the Guide shaft using standard dilation technique.

## 25.0 STEERABLE GUIDE CATHETER INSERTION

**WARNING:** Confirm a smooth transition between the Dilator and the tip of the Guide to minimize the risk of vascular and / or cardiac injury.

**CAUTION:** Always use pressure monitoring, echocardiography and fluoroscopy for guidance and observation during use of the MitraClip™ G4 System.

**WARNING:** Always use a careful, deliberate, and iterative approach to positioning the MitraClip™ G4 System. It is recommended to make multiple small adjustments rather than single large adjustments. Large adjustments may result in vascular and / or cardiac injury.

25.1 Rotate the +/- Knob in the “-” direction until the Guide curve is substantially straightened.

25.2 Wet the surface of the Guide shaft with sterile saline.

- 
- 25.3 Insert the Guide-Dilator assembly over the stationary guidewire into the femoral vein.
- WARNING:** DO NOT use excessive force to advance or manipulate the Guide-Dilator assembly. If resistance is encountered, use echocardiography and / or fluoroscopy to assess before proceeding. Use of excessive force may result in arrhythmias, vascular and / or cardiac injury, including creation of a clinically significant atrial septal defect.
- 25.4 Advance the Guide-Dilator assembly to the RA. Rotate the +/- Knob to Neutral, then place tip of the Dilator partially across the atrial septum.
- 25.5 Slowly dilate the atrial septum by gradually advancing the tip of the Guide-Dilator assembly.
- WARNING:** DO NOT rapidly advance the Guide-Dilator assembly across the atrial septum. Rapid advancement may result in vascular and / or cardiac injury.
- 25.6 Advance the Guide-Dilator assembly until the tip of the Guide extends approximately 3 cm in the LA.
- 25.7 Adjust Guide deflection and torque to position the tip away from adjacent tissues.
- 25.8 Place the Silicone Pad on the sterile drape over the Lift. Place the Stabilizer onto the Silicone Pad.
- 25.9 Secure the Guide in the Stabilizer slot using the Fastener. Ensure the Fastener engages the metallic tube on the Guide shaft. The Guide handle should be immediately adjacent to the Stabilizer, such that they are in contact with each other.
- 25.10 Retract the Dilator approximately 5 cm into the Guide, leaving the guide wire in the left upper pulmonary vein or LA.
- CAUTION:** Always loosen the Fastener before torquing the Guide to prevent stripping the screw.
- 25.11 Retract the guidewire into the tip of the Dilator. Remove the Dilator and guidewire while gently aspirating the Guide (starting when the Dilator is approximately halfway retracted into the Guide, approximately 40 cm) using a 50–60 cc syringe. Cover Guide Hemostasis Valve with finger upon Dilator removal.
- NOTE:** Avoid contacting tissue or creating a vacuum in the Guide lumen. If necessary, position the Guide handle below the level of the LA to allow blood to fill the Guide lumen.
- WARNING:** DO NOT create a vacuum while removing the dilator from the Guide; air may enter the lumen of the Guide which may result in air embolism.
- WARNING:** Failure to fully retract guidewire into the Dilator may result in air embolism.

---

## 26.0 CLIP DELIVERY SYSTEM INSERTION

- 26.1 Confirm the Guide lumen is completely de-aired.

**WARNING:** To minimize the potential of air embolism, **DO NOT** introduce the CDS into the Guide until the Guide lumen has been completely de-aired.

- 26.2 Confirm there is a slow, continuous heparinized saline flush through both the Sleeve and the DC and that the tip of the Clip is just proximal to the tip of the Clip Introducer.

**CAUTION:** Failure to continuously flush the CDS with heparinized saline may reduce device performance.

**WARNING:**

- Heparinized saline flush should be continuous throughout the procedure.
- Discontinuing flush may result in air embolism and / or thrombus formation.
- Ensure flow is visible through the drip chamber and that tubing is free from kinks and / or obstruction.
- Ensure pressure of 300 mmHg is maintained.

- 26.3 Carefully remove the protective cover surrounding the Clip and the Clip Introducer.

- 26.4 Confirm that the stopcock on the Clip Introducer flush port is closed and that the Clip Introducer is de-aired.

- 26.5 While flushing heparinized saline on the Guide Hemostasis Valve, place the tip of the Clip Introducer against the Guide Hemostasis Valve and advance the Clip Introducer straight into the valve in a continuous motion while rotating the Clip Introducer in small clockwise and counterclockwise motions until the Clip can be observed distal to the valve.

**WARNING:** **DO NOT** continue to advance the Clip Introducer if resistance is felt; the Guide Hemostasis Valve, Clip Introducer or the Clip may be damaged. Damage to these components may result in air embolism, vascular or cardiac injury.

**WARNING:** To minimize the potential of air embolism, ensure proper de-airing when inserting the Clip Introducer into the Guide Hemostasis Valve.

- 26.6 Leave the Clip Introducer fully inserted in the Guide Hemostasis Valve throughout the procedure.

- 26.7 Align the Longitudinal Alignment Marker on the Sleeve shaft with the Alignment Marker on the Guide Hemostasis Valve.

- 26.8 Turn the +/- Knob to neutral then carefully advance the CDS through the Guide under fluoroscopic guidance. Stop when tip of the Clip is even with the tip of the Guide.

---

**NOTE:** If resistance to CDS advancement is felt, reduce Guide deflection.

26.9 Under echocardiographic guidance, advance the CDS and retract the Guide iteratively as needed while maintaining the Guide in the LA. Stop when the Guide RO Tip Ring is between the RO Alignment Markers of the Sleeve, as confirmed under fluoroscopic guidance.

26.10 Position the Sleeve Handle in the Stabilizer slot.

26.11 Confirm that the Clip is free from the left atrial wall and valve tissue.

**WARNING:** Failure to confirm that the Clip is free from the left atrial wall and valve tissue may result in cardiac injury.

## **27.0 INITIAL MITRACLIP™ G4 SYSTEM POSITIONING IN THE LEFT ATRIUM**

**NOTE:** Positioning is achieved with iterative adjustments of the Guide and CDS using torque, translation, and knob adjustments. The goals of positioning are:

A. Positioning the Clip centrally over the valve with respect to anterior-posterior and medial-lateral directions.

B. Aligning the Clip so the DC Shaft is perpendicular to the plane of the mitral valve.

C. Positioning the distal tip of the Clip at least 1 cm above the leaflets.

**WARNING:** Excessive torque on the Guide and translation of the MitraClip™ G4 System may inadvertently displace the tip of the Guide from the LA, which may result in arrhythmias or cardiac injury.

**WARNING:** DO NOT continue to rotate or manipulate any of the knobs if significant resistance is noted; device damage may occur and result in cardiac injury.

27.1 Adjust the Guide position as necessary to maintain that the Clip is free from adjacent tissue.

27.2 Adjust Sleeve deflection using the M/L Knob and / or the A/P Knob to deflect the Clip towards the apex. Retract the DC Radiopaque Ring against the Sleeve tip as necessary.

27.3 During Sleeve deflections confirm that the Guide RO Tip Ring is between the RO Alignment Markers of the Sleeve prior to making maximum Sleeve deflections.

**WARNING:** DO NOT deflect the Sleeve tip more than 90 degrees as device damage may occur. Use of damaged product may result in cardiac injury.

27.4 Secure the Sleeve handle in the Stabilizer using the Fastener.

27.5 To reposition the MitraClip™ G4 System, move the Stabilizer with the system until positioning is adequate.



- 
- 27.6 Adjust the MitraClip™ G4 System position to maintain adequate height above the mitral valve in the LA.

**WARNING:** Maintain the Clip above the leaflets until ready to grasp to minimize the risk of Clip entanglement in the chordal apparatus. Clip entanglement may result in cardiac injury, worsening mitral regurgitation, difficulty or inability to remove the Clip and conversion to surgical intervention.

## 28.0 FINAL MITRACLIP™ G4 SYSTEM POSITIONING

- 28.1 *Raise the Gripper(s).*

**CAUTION:** Raising the Grippers more often than needed or retracting the Gripper Lever forcefully may damage the Gripper cover or Gripper line and impair CDS performance.

- 28.2 *Unlock the Clip and Open the Clip Arms* to approximately 180 degrees.

**WARNING:** DO NOT RETRACT THE LOCK LEVER FORCEFULLY. Retracting the Lock Lever forcefully may result in the inability to unlock Clip. Inability to open the Clip may result in valve injury or lead to deployment of the Clip in an unintended location.

- 28.3 Adjust the MitraClip™ G4 System to reposition the Clip as necessary. Confirm that the distal tip of the Clip is at least 1 cm above the leaflets.

- 28.4 Rotate the DC handle to align the Clip Arms perpendicular to the line of coaptation. DO NOT rotate the Clip more than 90 degrees in each direction.

- 28.5 Carefully translate the DC shaft multiple times to release stored torque. Fully retract the DC.

**WARNING:** Fully release stored torque. If not done, it may result in unwanted Clip Arm orientation changes during grasping. Torque of the DC Handle more than 180 degrees may result in DC damage and cardiac injury.

- 28.6 *Identify Gripper Orientation.*

- 28.7 Close the Clip to a *Clip Arm Angle* of approximately 60 degrees.

- 28.8 Complete final MitraClip™ G4 System positioning in the LA using multiple imaging planes. Re-secure the Guide and Sleeve Fasteners as appropriate.

## 29.0 GRASPING THE LEAFLETS AND VERIFYING THE GRASP

- 29.1 Advance the DC distally to position the Clip approximately 2 cm below the valve. Ensure that the Clip Arms are oriented perpendicular to the line of coaptation.

**WARNING:** Do confirm that the Clip Arms are perpendicular to the line of coaptation. Failure to do so may result in loss of leaflet capture and insertion. Loss of leaflet capture and insertion may result in inadequate reduction of mitral regurgitation and may result in a single leaflet device attachment (SLDA).

- 
- WARNING:** DO NOT make substantial Clip Arm orientation adjustment in the LV. Clip entanglement in chordae may result in cardiac injury and worsening mitral regurgitation; and may result in difficulty or inability to remove the Clip, and conversion to surgical intervention.
- WARNING:** Always ensure that either the Grippers are raised or that the Clip is closed while in the LV to avoid potential cardiac injury.
- 29.2 Open the Clip Arms to the Grasping Arm Angle.
- 29.3 Without using excessive force, retract the DC to grasp the anterior and posterior leaflets.
- WARNING:** An improper grasp will allow one or both leaflets to move freely. Closing and deploying the Clip in this situation may result in loss of leaflet capture and insertion. Loss of leaflet capture and insertion may result in inadequate reduction of mitral regurgitation and may result in a single leaflet device attachment (SLDA).
- 29.4 If the grasp appears satisfactory, Lower the Gripper(s) onto the leaflets.
- NOTE:** Simultaneous leaflet capture with both Grippers should be attempted first. If unsuccessful, Raise the Gripper(s) to release leaflet capture and Lower the Gripper(s) to capture leaflets. Raising and lowering Gripper(s) can be done simultaneously or independently.
- WARNING:** DO NOT advance the DC handle or adjust the position of the MitraClip™ G4 System in a way that increases tension on the leaflets after grasping the leaflets, as valve injury may occur.
- 29.4.1 If both Grippers have not lowered:
- 29.4.1.1 Lock the Clip.
- 29.4.1.2 Confirm both Grippers have lowered.
- 29.4.1.3 Unlock the Clip.
- WARNING:** Failure to confirm that both Grippers have been lowered onto the leaflets prior to closing the Clip may result in loss of leaflet capture and insertion. Loss of leaflet capture and insertion may result in inadequate reduction of mitral regurgitation and may result in a single leaflet device attachment (SLDA).
- 29.5 Close the Clip until the Clip Arm Angle is approximately 60 degrees. Release tension on the DC and secure the DC Fastener.
- 29.6 Use echocardiographic imaging to verify insertion of both leaflets and satisfactory grasp by observation of:
- Leaflet immobilization
  - Single or multiple valve orifice(s)
  - Limited leaflet mobility relative to the tips of both Clip Arms
  - Adequate MR reduction.

---

29.6.1 If grasping fails to hold both leaflets and the Clip retracts to the LA, reposition the MitraClip™ G4 System.

29.6.1.1 *Open the Clip Arms* and reorient the Clip Arms in the LA, as needed, then repeat grasping steps (refer to Section 28.0 and 29.0).

29.6.1.1.1 If significant repositioning is necessary, *Fully Close the Clip Arms* and *Lower the Gripper(s)* then repeat positioning and grasping steps.

29.6.2 If the Sleeve limits DC travel during grasping, an inadequate grasp may require repositioning of the MitraClip™ G4 System.

28.6.2.1 *Raise the Gripper(s)* and *Open the Clip Arms* to approximately 180 degrees and advance the DC handle. Repeat positioning and grasping steps as necessary (refer to Section 28.0 and 29.0).

### 30.0 CLOSING THE CLIP AND EVALUATING CLIP POSITION

30.1 *Lock the Clip.*

**WARNING:** Failure to Lock the Clip may result in loss of leaflet capture and insertion. Loss of leaflet capture and insertion may result in inadequate reduction of mitral regurgitation and may result in a single leaflet device attachment (SLDA).

30.2 Slowly close the Clip just until the leaflets are coapted and MR is sufficiently reduced. The Clip should maintain a distinct “V” shape.

**WARNING:** DO NOT use excessive force to close the Clip further than is necessary to adequately reduce MR. Leaflet injury may occur. DO NOT close the Clip too tightly as it may result in leaflet injury or an inability to deploy the Clip. Inability to deploy the Clip may result in worsening mitral regurgitation, cardiac injury, a single leaflet device attachment (SLDA), and / or conversion to surgical intervention.

**WARNING:** Failure to turn the Arm Positioner at least ½ turn in the “Close” direction after locking the Clip may result in loss of leaflet capture and insertion. Loss of leaflet capture and insertion may result in inadequate reduction of mitral regurgitation and may result in a single leaflet device attachment (SLDA).

30.3 Use echocardiographic imaging to verify valve function, satisfactory coaptation, and insertion of both leaflets by observation of:

- Leaflet immobilization
- Single or multiple valve orifice(s)
- Limited leaflet mobility relative to the tips of both Clip Arms
- Adequate MR reduction.

---

30.3.1 If the Clip position is not satisfactory, *Raise the Gripper(s), Unlock the Clip and Invert the Clip Arms.*

**WARNING:** Turning the Arm Positioner in the “Open” direction more than 1 full turn past a *Clip Arm Angle* of 180° or turning past when resistance is first noted may result in device damage which could cause the Clip to become non-functional and lead to embolization, and / or conversion to surgical intervention.

30.3.2 Retract the inverted Clip into the LA.

30.3.3 Confirm both leaflets move freely.

30.3.4 Repeat positioning steps, as necessary, then repeat grasping steps.

### 31.0 MITRACLIP™ G4 IMPLANT PRE-DEPLOYMENT CLIP ASSESSMENT

31.1 Confirm DC Handle is secure.

**WARNING:** Do secure the DC Handle. If not done, it may result in leaflet injury or loss of leaflet insertion with resultant worsening mitral regurgitation, single leaflet device attachment (SLDA), and / or conversion to surgical intervention.

31.2 *Establish Final Arm Angle.*

**WARNING:** DO NOT turn the Arm Positioner more than 1 turn in the “Open” direction from neutral. Failure to stop turning the Arm Positioner at 1 turn in the “Open” direction past neutral may result in Clip opening or device damage which could cause the Clip to become non-functional and lead to embolization and / or conversion to surgical intervention.

31.3 Turn the Arm Positioner to the “closed” side of the neutral position.

31.4 Perform mean pressure gradient assessment prior to proceeding to deployment.

### 32.0 CLIP DEPLOYMENT

#### 32.1 Deployment Step 1: Lock Line Removal

32.1.1 While holding the ends of the Lock Line remove the Lock Lever Cap and “O” ring. Unwrap the two ends of the Lock Line in a counterclockwise direction. Separate the ends of the Lock Line and remove the plastic cover from the lines so that no twists or knots are present.

**WARNING:** Do not let Line unravel freely. Do not remove Lock Line or plastic covers if line is bunched. Letting Line unravel freely may result in knots in the line. Removing Line if it is bunched may result in difficulty or inability to remove line due to knots or twists.

32.1.2 Grasp one of the free ends of the Lock Line, confirm the line moves freely, and slowly remove the Lock Line. Pull the Lock Line coaxial to the Lock Lever. If resistance is noted, stop and pull on the other free end to remove the Lock Line.

---

### 32.1.3 *Establish Final Arm Angle.*

**NOTE:** The Clip Arms may open slightly (~5°) and then remain in a stable position. If Arms open more than slightly, close the Clip to the desired Arm position and re-*Establish Final Arm Angle*.

**WARNING:** **DO NOT** turn the Arm Positioner more than 1 turn in the “Open” direction from neutral. Failure to stop turning the Arm Positioner at 1 turn in the “Open” direction past neutral may result in Implant opening or device damage which could cause the Implant to become non-functional and lead to embolization and / or conversion to surgical intervention.

### 32.1.4 *Turn the Arm Positioner to Neutral.*

## 32.2 **Deployment Step 2: Delivery Catheter Shaft Detachment**

32.2.1 Confirm that the Arm Positioner is Neutral. Remove the Release Pin from the DC Handle.

32.2.2 Turn the Arm Positioner in the “Open” direction until the Release Pin groove is fully exposed.

**NOTE:** After the Release Pin is removed, turning the Arm Positioner in the “Open” direction will not open the Clip Arms.

32.2.3 Turn the Actuator Knob of the DC approximately 8 turns in the direction of the arrow printed on the Actuator Knob.

If it is difficult to turn the Actuator Knob, STOP and confirm that the Arm Positioner has been turned in the “Open” direction, such that the Release Pin groove is fully exposed.

**WARNING:** **Stop turning the Actuator Knob when resistance is felt, otherwise it may result in inability to deploy the Clip. Inability to deploy the Clip may result in worsening mitral regurgitation, cardiac injury, a single leaflet device attachment (SLDA), and / or conversion to surgical intervention.**

32.2.4 Retract the Actuator Knob approximately 0.5 cm after it is fully unthreaded.

32.2.5 Fully retract the Gripper Levers.

**WARNING:** **Do retract the Gripper Levers fully prior to retracting the DC Handle. Failure to retract the Gripper Levers may result in higher forces to deploy the clip. This may result in worsening mitral regurgitation, cardiac injury or a single leaflet device attachment (SLDA).**

32.2.6 Release the DC Fastener, and slowly retract the DC Handle until the DC Radiopaque Ring is against the tip of the Sleeve.

---

32.2.6.1 If resistance is felt during DC Detachment and Clip separation from the DC can be confirmed by fluoroscopy, confirm the Gripper Levers are fully retracted. If resistance is still felt, access the Gripper Lines.

32.2.6.2 If resistance is felt during DC Detachment and Clip separation from the DC cannot be confirmed by fluoroscopy:

32.2.6.2.1 Confirm that the Actuator Knob is retracted approximately 0.5 cm beyond the fully exposed Release Pin Groove and the Gripper Levers are fully retracted.

32.2.6.2.2 If resistance is still felt during DC Detachment, confirm by fluoroscopy that the DC Shaft and the Clip are coaxially aligned. If they are not aligned, secure the DC fastener and slowly translate the Stabilizer and / or rotate the Guide until the angulation between the DC Shaft and the Clip is reduced or eliminated (as observed under fluoroscopic imaging).

32.2.6.2.3 If resistance is still felt, release the DC Fastener, retract the Actuator Knob an additional 0.5 cm.

32.2.6.2.4 If resistance is still felt, secure the DC Fastener and fully advance both Gripper Levers.

32.2.6.2.5 If resistance is still felt, access the Gripper Lines.

32.2.6.2.6 If resistance is still felt, secure the DC Fastener and repeat alignment steps to reduce / eliminate the angulation between the DC shaft and the Clip (as observed under fluoroscopic imaging).

32.2.7 Confirm that the DC Fastener is secure.

**32.2.8 Confirm that the Clip position is stable.**

32.2.9 Use echocardiographic imaging to verify valve function, satisfactory coaptation, and insertion of both leaflets by observation of:

- Leaflet immobilization
- Single or multiple valve orifice(s)
- Limited leaflet mobility relative to the tips of both Clip Arms
- Adequate MR reduction.

**32.2.10 If placing an additional Clip proceed to Section 33.0. If not placing an additional Clip proceed to Section 34.0.**

### **33.0 ADDITIONAL MITRACLIP™ G4 IMPLANT PLACEMENT**

**WARNING:** Use caution not to displace or dislodge an implanted Clip when placing an additional Clip; Clip detachment from leaflet(s) may occur which may result in a single leaflet device attachment (SLDA) or device embolization.

- 
- 33.1 When placing an additional Clip, the following are recommended:
- 33.1.1 In the LA, ensure Clip Arms are oriented perpendicular to the line of coaptation and Grippers are raised.
  - 33.1.2 Use both fluoroscopy and echocardiography when crossing into the LV and during grasping.
  - 33.1.3 Cross into the LV with a *Clip Arm Angle* of < 60 degrees.

**WARNING:** DO NOT use excessive force or retraction distance during grasping. This may compromise leaflet capture and insertion. Loss of leaflet capture and insertion may result in inadequate reduction of mitral regurgitation and may result in a single leaflet device attachment (SLDA).

## 34.0 MITRACLIP™ G4 SYSTEM REMOVAL

**WARNING:** During MitraClip™ G4 System removal always retract the CDS by pulling only on the Sleeve Handle. Retracting the CDS by pulling on the DC Handle may result in device damage and / or device or component embolization and may result in vascular and / or cardiac injury.

**WARNING:** Do release the DC Fastener before releasing Sleeve curves otherwise it may result in device damage and / or device or component embolization.

**WARNING:** Use echocardiographic guidance while releasing Sleeve deflection. Failure to do so may result in cardiac injury.

### 34.1 MitraClip™ G4 System Removal After Clip Deployment

#### 34.1.1 Removal of the CDS While Leaving the Guide in Place.

- 34.1.1.1 Release the DC Fastener.
- 34.1.1.2 Slowly release Sleeve deflection by rotating the M/L Knob and the A/P Knob to neutral.
- 34.1.1.3 Secure DC Fastener once Sleeve curves are released.
- 34.1.1.4 Straighten the Guide with the +/- Knob when the Delivery Catheter tip is free from the left atrial wall and the mitral valve.
- 34.1.1.5 Release the Sleeve Fastener and retract the CDS approximately 10 cm into the Guide by pulling only on the Sleeve Handle.
- 34.1.1.6 Confirm that the Clip Introducer is still fully advanced in the Guide Hemostasis Valve.
- 34.1.1.7 Retract the CDS by pulling only on the Sleeve Handle and position the Delivery Catheter tip inside the Clip Introducer. Begin gently aspirating the Guide (starting when the CDS is approximately halfway into the Guide, approximately 40 cm retracted) using a 50–60 cc syringe.



- 
- 34.1.1.8 Remove the CDS and the Clip Introducer simultaneously from the Guide by pulling on the Sleeve shaft and Clip Introducer. Ensure the Delivery Catheter tip is inside the Clip Introducer by visualizing the Proximal Sleeve alignment marker just outside the Clip Introducer. Aspirate the Guide during removal of the CDS and Clip Introducer. Cover Guide Hemostasis Valve with finger upon CDS removal. If necessary, position the Guide Handle below the level of the LA to allow blood to fill the Guide Lumen.

**WARNING: DO NOT remove the tip of the CDS from the Guide without removing the Clip Introducer simultaneously. Failure to remove the Clip Introducer simultaneously may result in air embolism.**

**WARNING: DO NOT create a vacuum while removing the CDS from the Guide; air may enter the lumen of the Guide which may result in air embolism.**

- 34.1.1.9 Aspirate using a 50–60 cc syringe to remove any remaining air from the Guide.

34.1.2 Removal of the CDS and Guide simultaneously.

- 34.1.2.1 Release the DC Fastener.
- 34.1.2.2 Slowly release Sleeve curves by rotating the M/L Knob and the A/P Knob to neutral.
- 34.1.2.3 Secure the DC Fastener once Sleeve curves are released.
- 34.1.2.4 Straighten the Guide with the +/- Knob when the Delivery Catheter tip is free from the left atrial wall and the mitral valve.
- 34.1.2.5 Release the Sleeve Fastener and retract the CDS approximately 10 cm into the Guide by pulling only on the Sleeve Handle.
- 34.1.2.6 Carefully retract the Guide tip into the RA. The Guide may be straightened further with the +/- Knob if desired.
- 34.1.2.7 Remove the MitraClip™ G4 System from the femoral vein, while providing hemostasis.

**34.2 MitraClip™ G4 System Removal With Clip Attached**

34.2.1 Removal of the CDS while leaving the Guide in place.

- 34.2.1.1 Confirm Clip is locked.
- 34.2.1.2 *Fully Close the Clip Arms* and continue to turn the Arm Positioner until it is no longer possible to rotate the Arm Positioner. Return the *Arm Positioner to Neutral*.



---

**WARNING: Failure to follow Step 34.2.1.2 prior to retraction into the Guide may result in device damage, inability to remove the CDS and / or vascular and cardiac injury.**

34.2.1.3 *Lower the Gripper(s).*

34.2.1.4 Release the DC Fastener and retract the DC Handle until the DC Radiopaque Ring is fully against the tip of the Sleeve.

34.2.1.5 Slowly release Sleeve deflection by rotating the M/L Knob and the A/P Knob to neutral.

34.2.1.6 Rotate DC handle such that the clip arms are perpendicular to the guide curve plane.

34.2.1.7 Secure the DC Fastener once Sleeve curves are released.

34.2.1.8 Straighten the Guide with the +/- Knob when the tip of the MitraClip™ G4 Implant is free from the left atrial wall and the mitral valve.

**WARNING: Straighten the Guide prior to retracting the Clip into the Guide. If not done, it may result in device damage, inability to remove the CDS and / or vascular and cardiac injury.**

34.2.1.9 Release the Sleeve Fastener and retract the CDS into the Guide by pulling only on the Sleeve Handle.

**NOTE:** If resistance is noted, advance and rotate the Clip by rotating the DC Handle then retract the CDS into the Guide. The Guide and / or Sleeve position may also be adjusted to facilitate Clip entry into the Guide. If necessary, retract the Sleeve or advance the Clip to create a 2-3cm separation to facilitate Clip entry into the Guide.

**WARNING: Use fluoroscopic guidance while retracting the CDS into the Guide. Failure to do so may result in device damage, inability to remove the CDS and / or vascular and cardiac injury.**

34.2.1.10 Confirm that the Clip Introducer is still fully advanced in the Guide Hemostasis Valve.

34.2.1.11 Retract the CDS by pulling only on the Sleeve Handle and position the Clip inside the Clip Introducer. Begin gently aspirating the Guide (starting when the CDS is approximately halfway into the Guide, approximately 40 cm retracted) using a 50–60 cc syringe.

---

34.2.1.12 Remove CDS and Clip Introducer simultaneously from the Guide by pulling on the Sleeve shaft and Clip Introducer. Ensure the Clip is inside the Clip Introducer by visualizing the Proximal Sleeve alignment marker just outside the Clip Introducer. Aspirate the Guide during removal of the CDS and Clip Introducer. If necessary, position the Guide Handle below the level of the LA to allow blood to fill the Guide lumen.

**WARNING:** DO NOT remove the tip of the CDS from the Guide without removing the Clip Introducer simultaneously and with the Clip inside the Clip Introducer. Failure to remove the Clip Introducer simultaneously may result in air embolism.

**WARNING:** DO NOT create a vacuum while removing the CDS from the Guide; air may enter the lumen of the Guide which may result in air embolism.

**WARNING:** DO NOT re-use the CDS after removal. Replace the CDS with a new device. Reinserting the CDS after removal may result in inability to open the Clip. Inability to open the Clip may result in valve injury or lead to deployment of the Clip in an unintended location.

34.2.1.13 Aspirate using a 50–60 cc syringe to remove any remaining air from the Guide.

34.2.2 Simultaneous removal of CDS and Guide.

34.2.2.1 Confirm Clip is locked.

34.2.2.2 *Fully Close the Clip Arms* and continue to turn the Arm Positioner until it is no longer possible to rotate the Arm Positioner. Return the *Arm Positioner to Neutral*.

**WARNING:** Failure to follow Step 34.2.2.2 prior to retraction into the Guide may result in device damage, inability to remove the CDS and / or vascular and cardiac injury.

34.2.2.3 *Lower the Gripper(s).*

34.2.2.4 Release the DC Fastener and retract the DC Handle until the DC Radiopaque Ring is fully against the tip of the Sleeve.

34.2.2.5 Slowly release Sleeve deflection by rotating the M/L Knob and the A/P Knob to neutral.

34.2.2.6 Rotate DC handle such that the Clip Arms are perpendicular to the guide curve plane.

34.2.2.7 Secure the DC Fastener once Sleeve curves are released.

- 
- 34.2.2.8 Straighten the Guide with the +/- Knob when the tip of the MitraClip™ G4 Implant is free from the left atrial wall and the mitral valve.

**WARNING:** Straighten the Guide prior to retracting the Clip into the Guide. Failure to do so may result in device damage, inability to remove the CDS and / or vascular and cardiac injury.

- 34.2.2.9 Release the Sleeve Fastener and retract the CDS approximately 10 cm into the Guide by pulling only on the Sleeve Handle.

**NOTE:** If resistance is noted, advance and rotate the Clip by rotating the DC Handle then retract the CDS into the Guide. The Guide and / or Sleeve position may also be adjusted to facilitate Clip entry into the Guide. If necessary, retract the Sleeve or advance the Clip to create a 2-3 cm separation to facilitate Clip entry into the Guide.

**WARNING:** Use fluoroscopic guidance while retracting the CDS into the Guide. Failure to do so may result in device damage, inability to remove the CDS and / or vascular and cardiac injury.

- 34.2.2.10 Carefully retract the Guide tip into the RA. The Guide may be straightened further with the +/- Knob if desired.

- 34.2.2.11 Remove the MitraClip™ G4 System from the femoral vein, while providing hemostasis.

™ Indicates a trademark of the Abbott group of companies.

Pat. [www.abbott.com/patents](http://www.abbott.com/patents)

© 2023 Abbott. All Rights Reserved.



**Abbott Medical**

177 County Road B East  
St. Paul, MN 55117 USA

TEL: +1 855 478 5833





















FAX: +1 651 756 5833

Outside USA TEL: +1 855 478 5833

Outside USA FAX: +1 651 756 5833

## GRAPHICAL SYMBOLS FOR MEDICAL DEVICE LABELING

Source reference for recognized symbols retained at the manufacturer in the Medical Device File.

	Batch code		Do not re-sterilize
	Catalogue number		Do not re-use
	Use-by date		Non-pyrogenic
	Sterilized using ethylene oxide		Do not use if package is damaged and consult instructions for use
	Single sterile barrier system		Manufacturer
	Consult instructions for use or consult electronic instructions for use		Date of manufacture
	Refer to instruction manual/booklet (Symbol Color: Blue)		Packaging unit
	MR Conditional		Keep away from sunlight
	Unique device identifier		Keep dry
	Medical device		CAUTION: Federal law restricts this device to sale by or on the order of a physician.